



Aug 14,  
2019

## FLASH v2.0 V.2


Amy Lyden<sup>1</sup>, Gloria Castañeda<sup>1</sup>, Emily Crawford<sup>1</sup>, Saharai Caldera<sup>2</sup>, Jenai Quan<sup>1</sup>, Katherine Travisano<sup>1</sup>

<sup>1</sup>CZ Biohub, <sup>2</sup>CZ Biohub, UCSF

1 Works for me dx.doi.org/10.17504/protocols.io.6iwhcfe

Amy Lyden ⚡ ⚡

### ABSTRACT

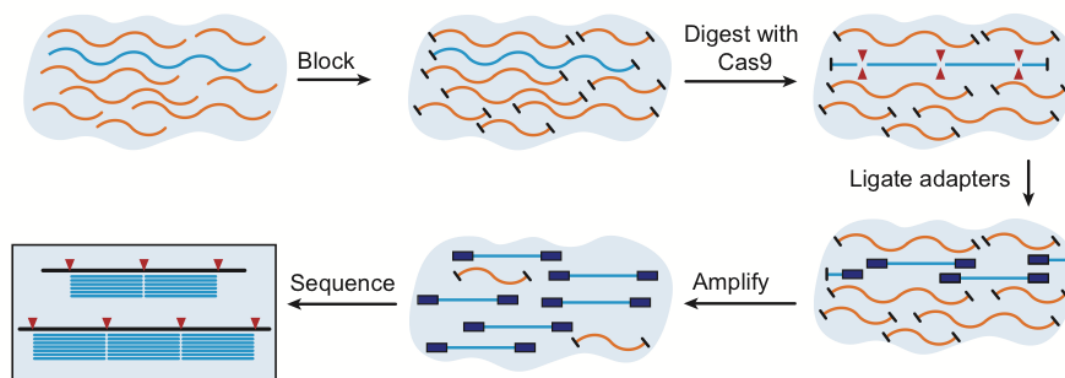

FLASH is a crispr-cas9

### FLASH workflow

technology that enriches for targeted sequences in sequencing libraries. The initial DNA sample undergoes a blocking step that removes of the 5' phosphoryl groups of the DNA fragments, resulting in a product that is not amendable for downstream adaptor ligation or amplification via standard Illumina-based library preparation. The subsequent incorporation of targeted CRISPR-cas9 library exposes the desired regions of interest, allowing them to be processed into a library. For more information on methods and results, please see the [FLASH paper](#). For FLASH guide RNA design help, please see our [github](#).

### GUIDELINES

1. This protocol has been used most successfully with starting inputs of DNA ranging from 10 pg – 100 ng. Limited success may be achieved with as little as 100 fg DNA.
2. Keep all enzymes on a chilled enzyme block. Immediately before use, allow them to come to room temperature for 5 minutes and then vortex briefly.
3. The NEBNext® dA-Tailing Module includes the Klenow fragment and the dA-Tailing buffer. The Klenow fragment allows for the 5' ends to remain dephosphorylated after the initial blocking step. Do NOT use the NEBNext Ultra II End Prep kit for this step.
4. AmpureXP beads or other magnetic beads may be used instead of SPRI beads. Adjust the ratios of beads:sample accordingly to ensure proper removal of unwanted products. Different beads yield varying size selection cut-offs. Refer to this table to adjust ratios according to desired effect.

SPRI Beads : sample ratio	Step(s)	Desired cut-off/effect
1.7 : 1	Post-FLASH cleanup	No size selection; cleanup of unwanted buffers and deactivated cas9 protein from reaction
1 : 1	Post-adaptor ligation clean up	Removal of buffers AND Stringent removal of all adaptor dimers

1: 0.9	post-indexing Q5 PCR cleanup	Removal of buffers AND Stringent removal of all primer dimers
1: 1	post-KAPA amplification cleanup	Removal of buffers AND Stringent removal of all primer dimers

SPRI bead cut-offs

## MATERIALS

NAME	CATALOG #	VENDOR
NEBNext dA-Tailing Module - 100 rxns	E6053L	New England Biolabs
Sodium Orthovanadate (Vanadate) - 1 ml	P0758S	New England Biolabs
USER Enzyme - 250 units	M5505L	New England Biolabs
NEBNext Adaptor for Illumina	E7337 in Kits E7335, E7500, E771	New England Biolabs
NEBNext Ultra II Ligation Module - 96 rxns	E7595L	New England Biolabs
NEBNext Ultra II Q5 Master Mix - 250 rxns	M0544L	New England Biolabs
Qubit dsDNA HS Assay kit	Q32854	Thermo Fisher Scientific
Nuclease-free water	AM9932	Ambion
Bioanalyzer chips and reagents (DNA High Sensitivity kit)		Agilent Technologies
PCR Thermocycler		
TruSeq i7/i5 Indexing Primers - Custom (or NEBNext® Multiplex Oligos for Illumina)	E7500L	New England Biolabs
Proteinase K	P8107S	New England Biolabs
SPRI beads (homemade) or Ampure XP beads	View	
Kapa HiFi Real-Time Amplification Kit	KK2702	Kapa Biosystems
Magnetic rack for PCR strips	View	
rAPID alkaline phosphatase enzyme and buffer	4898133001	Sigma Aldrich
cas9 4µM or higher concentration		Homemade
Dual guide RNAs (4µM - targeted to genes or regions to be depleted - crisprRNA and tracr RNA - quantified by RNA Qubit)		Homemade
Illumina P5 and P7 primers 5uM combined; P5: 5' AATGATACGGCGACCACCGAGATCT P7: 5' CAAGCAGAGACGGCATACGAGAT	View	
10X Cas9 Activity Buffer (500mM Tris pH 8.0 - 1M NaCl - 100mM MgCL2 - 10mM TCEP)		

## STEPS MATERIALS

NAME	CATALOG #	VENDOR
rAPID alkaline phosphatase enzyme and buffer	4898133001	Sigma Aldrich
Dual guide RNAs (4µM - targeted to genes or regions to be depleted - crisprRNA and tracr RNA - quantified by RNA Qubit)		Homemade
cas9 4µM or higher concentration		Homemade
10X Cas9 Activity Buffer (500mM Tris pH 8.0 - 1M NaCl - 100mM MgCL2 - 10mM TCEP)		

NAME ▾	CATALOG # ▾	VENDOR ▾
Proteinase K	P8107S	New England Biolabs
SPRI beads (homemade) or Ampure XP beads	<a href="#">View</a>	
Klenow Fragment (3'-5' exo-) - 1,000 units	M0212L	New England Biolabs
NEBNext Adaptor for Illumina	<a href="#">View</a>	New England Biolabs
NEBNext Ultra II Ligation Module - 96 rxns	E7595L	New England Biolabs
80% Ethanol		
TruSeq i7/i5 Indexing Primers - Custom (or NEBNext® Multiplex Oligos for Illumina)	E7500L	New England Biolabs
SPRI beads (homemade) or Ampure XP beads	<a href="#">View</a>	
80% Ethanol		
NEBNext Ultra II Q5 Master Mix - 250 rxns	M0544L	New England Biolabs
USER Enzyme - 250 units	M5505L	New England Biolabs
SPRI beads (homemade) or Ampure XP beads	<a href="#">View</a>	
80% Ethanol		
Kapa HiFi Real-Time Amplification Kit	KK2702	Kapa Biosystems
Illumina P5 and P7 primers 5uM combined; P5: 5' AATGATACGGCGACCACCGAGATCT P7: 5' CAAGCAGAAGACGGCATACGAGAT	<a href="#">View</a>	
SPRI beads (homemade) or Ampure XP beads	<a href="#">View</a>	
80% Ethanol		
Qubit 1X dsDNA High Sensitivity Assay Kit	Q33230	Thermo Fisher Scientific
Bioanalyzer chips and reagents (DNA High Sensitivity kit)		Agilent Technologies
Sodium Orthovanadate (Vanadate) - 1 ml	P0758S	New England Biolabs

#### BEFORE STARTING


Ensure that you are working in a PCR hood in a pre-PCR space if you are working with metagenomic samples or in a PCR hood if you are working with isolate samples. FLASH is very sensitive to environmental contamination.


Please refer to the guidelines section in this protocol if you are using Ampure beads or other SPRI beads to ensure you use the correct cut-offs.

For mixing of sample prior to PCR, avoid vortexing to keep DNA intact, and instead mix by pipetting and tapping sides of tube.

#### Dephosphorylation

- 1 Normalize your cDNA or gDNA to anywhere between 10pg-100ng. For most samples, we recommend 5-10ng input.




- 2 Prepare a reaction for each cDNA or gDNA sample and mix well. Add the components in the order specified below. You can make a master mix (MM) of the rapid alkaline phosphatase (RAP) buffer and enzyme. If you are using the RAP MM, mix  **3 µl** of MM with each sample. Mix thoroughly with a pipette or by tapping to avoid shearing.




**rAPID alkaline phosphatase enzyme and buffer**  
by [Sigma Aldrich](#)  
Catalog #: [4898133001](#)

Component	1X
DNA, 10pg - 100ng	x µL
rAPid Alkaline Phosphatase Buffer	2 µL
rAPid Alkaline Phosphatase	1 µL
H2O	up to 20 µL



RAP MM.


- 3 Incubate at  **37 °C** for  **00:30:00** with heated lid OFF.
- 4 Add  **1 µl** sodium orthovanadate (competitive inhibitor of phosphatases) to the quench the reaction and mix well with a pipette or by tapping.



**Sodium Orthovanadate (Vanadate) - 1 ml**  
by [New England Biolabs](#)  
Catalog #: [P0758S](#)


#### Cas9 Treatment

- 5 Prepare dual-guide RNAs by annealing crRNA and tracrRNA at an equimolar amount at  **95 °C** for  **00:00:30** then allowing the mixture to cool to room temperature on the bench. (Note. If crRNA and tracrRNA have been previously annealed and stored, it is recommended to re-anneal the dual-guide RNAs under the above conditions)




**Dual guide RNAs (4µM - targeted to genes or regions to be depleted - crisprRNA and tracr RNA - quantified by RNA Qubit)**  
by [Homemade](#)

- 6 If your starting stock of Cas9 is more than 4μM, you must dilute your stock of Cas9 to 4μM by using 1X Cas9 activity buffer.



cas9 4μM or higher concentration  
by Homemade



10X Cas9 Activity Buffer (500mM Tris pH 8.0 - 1M NaCl - 100mM MgCL2 - 10mM TCEP)

- 7 Make a Cas9 master mix as described below. **Add the components in the order specified to prevent precipitation.**





Component	1X	___ X
10x Cas9 Activity buffer	3 μL	___ μL
Cas9 4μM*	2.5 μL	___ μL
dgRNAs 4μM**	3 μL	___ μL
H2O	0.5 μL	___ μL
Total	9 μL	___ μL




Cas9 MM


*\*Remember to dilute your Cas9 stock to 4uM if not already at 4uM. Use 1x Cas9 activity buffer to dilute your Cas9 enzyme if you do need to dilute your stock.*

*\*\*Remember to dilute your dgRNAs to 4uM using water if needed.*

Note: Most experiments for the 2018 FLASH paper were performed at this Cas9 concentration. However, we have demonstrated that lower concentrations work equally well on bacterial isolate DNA. Consult the manuscript for more details.

- 8 Add  9 μl of the master mix to each of your  21 μl blocked DNA samples.
- 9 Mix well by pipetting or tapping the PCR tubes and incubate the reaction at  37 °C for  02:00:00 .


- 10 Deactivate the Cas9 by adding  **1  $\mu$ l** of Proteinase K to each of your sample tubes, mixing again by pipetting or tapping. Incubate at  **37 °C** for  **00:15:00**.









**Proteinase K**  
by [New England Biolabs](#)  
Catalog #: [P8107S](#)


#### SPRI Clean-up at 1.7X

- 11 Equilibrate clean SPRI beads to room temperature and vortex well to mix.







**SPRI beads (homemade) or Ampure XP beads**  
[View](#)

- 12 Add beads equivalent to 1.7X the sample volume to each sample tube (for  **31  $\mu$ l** of sample, add  **53  $\mu$ l** beads).
- 13 Mix well by pipetting or tapping the tubes. Pulse-spin in a picofuge for no more than 2 seconds.
- 14 Incubate for  **00:05:00** at room temperature, then put the tubes on the magnetic rack. Allow beads to separate on the magnet for 3-5 minutes, or until the supernatant is clear.
- 15 Keeping the tubes on the magnet, carefully remove and discard the supernatant.
- 16 Add  **200  $\mu$ l**  **80 % volume** ethanol (prepared fresh). Incubate beads for  **00:01:00** and then remove the ethanol



**80% Ethanol**

- 17 Repeat the above ethanol wash step.

- 18 Allow the beads to air dry for  **00:05:00**. Do not overdry. Dry beads should appear matte (rather than glossy), but should not have a cracked appearance. Overdried beads, as indicated by a cracked appearance, may not resuspend or elute well.
- 19 Remove tubes from magnet and resuspend in  **53 µl** nuclease-free H<sub>2</sub>O by pipetting up and down.
- 20 Continue to resuspend by tapping the tubes and then spin down briefly in a picofuge.
- 21 Allow  **00:02:00** for DNA to elute from beads, then transfer tubes back to magnet
- 22 Allow the beads to separate for at least 2 minutes.
- 23 Collect  **50.4 µl** of supernatant to clean PCR tubes.

#### dA tailing

- 24 Prepare the following mixture for each sample. You can prepare the dA-tailing buffer and Klenow fragment ahead of time as a MM. If you prepared dA-tailing MM, add 9.6µL of MM to each sample.




Component	1X
FLASHed sample	50.4 µL
dA-Tailing buffer	6 µL
Klenow fragment	3.6 µL
<i>Total</i>	<i>60</i> <i>µL</i>



**Klenow Fragment (3'-5' exo-) - 1,000 units**

by New England Biolabs

Catalog #: [M0212L](#)

- 25 Mix well by pipetting up and down several times with a P200 set to 40µL.
- 26 Incubate all tubes at  **37 °C** for  **00:30:00** with heated lid OFF
- 27 Cool all tubes to  **4 °C** and **proceed with the next part as soon as possible.**

28



! The reagents used in this step are very viscous and must be mixed well before using.

! Do NOT make a master mix for this step, although the Ligation Master Mix and Ligation Enhancer may be mixed up to 4 hours before and kept at 4°C.

- 29 Prepare the following mixture for each sample. Alternatively, if you prepared an adaptor MM with Ligation MM and Ligation Enhancer, add 31uL of MM to each sample and then add 2.5uL of adaptor to each sample.

Component	1X
dA-tailed sample from part IV	60 µL
NEB Ultra II Ligation Master Mix	30 µL
NEBNext Ligation Enhancer	1 µL
NEBNext Adaptor 1:100 or 1:300 dilution*	2.5 µL
<i>Total</i>	<i>93.5 µL</i>

*Note: Adaptor dilution of 1:100 works for an initial DNA input up to 100 ng. A 1:300 dilution is recommended for an initial DNA input of under 10 ng. The subsequent cleanup step should remove all extra adaptor and adaptor dimers.*



NEBNext Adaptor for Illumina

by New England Biolabs

[View](#)



NEBNext Ultra II Ligation Module - 96 rxns

by New England Biolabs

Catalog #: [E7595L](#)

- 30 Prepare the above mixture and mix well by pipetting up and down several times with a P200 set to 50 µl.

- 31 Incubate at 20 °C for 00:15:00 in a thermocycler with the heated lid OFF.



## SPRI Clean Up 1X + Addition of TruSeq Indexing Primers

- 32 Thaw a TruSeq i5/i7 barcode plate or other TruSeq primers, and choose barcodes for each sample. Take note of plate color/barcodes to be used. **DO NOT** use the same barcode for more than one sample in a sequencing run. Dual unique TruSeq barcodes are preferable.



TruSeq i7/i5 Indexing Primers - Custom (or NEBNext® Multiplex Oligos for Illumina)  
by New England Biolabs  
Catalog #: E7500L

- 33 Equilibrate clean SPRI beads to room temperature and vortex well to mix.



SPRI beads (homemade) or Ampure XP beads  
[View](#)

- 34 Add beads equivalent to 1X the sample volume to each sample tube (for 93.5 µl of sample, add 93.5 µl beads).



- 35 Mix well by pipetting or tapping the tubes. Pulse-spin in a picofuge for no more than 2 seconds.

- 36 Incubate for 00:15:00 at room temperature, then put the tubes on the magnetic rack. Allow beads to separate on the magnet for 3-5 minutes, or until the supernatant is clear.




Longer incubation time during SPRI addition is because it may be helpful in removing unwanted fragments/adaptor dimers


- 37 Keeping the tubes on the magnet, carefully remove and discard the supernatant.

38 Add  200  $\mu$ l [M]80 % volume ethanol (prepared fresh). Incubate beads for  00:01:00 and then remove the ethanol




39 Repeat the above ethanol wash step.

40 Allow the beads to air dry for  00:05:00 . Do not overdry. Dry beads should appear matte (rather than glossy), but should not have a cracked appearance. Overdried beads may not resuspend or elute well. Ensure the beads are fully dry before eluting.

41 Remove tubes from magnet and resuspend in  17  $\mu$ l nuclease-free H<sub>2</sub>O.


42 Resuspend well by tapping the tubes and spin down briefly in a picofuge.

43 Allow  00:02:00 for DNA to elute from beads, then transfer tubes back to magnet

44 Allow the beads to separate for at least 2 minutes.

45 Collect  15  $\mu$ l of elution.


46 Mix  10  $\mu$ l of the elution with the appropriate TruSeq indexing primer barcodes as planned.

 STOPPING POINT: If necessary, samples may be stored at -20C.


#### USER Enzyme + Indexing PCR

47 Spin down plate or tube strips briefly to collect liquid to bottom of the well.

- 48 Prepare the following mixture for each sample. You can make a master mix of USER enzyme and Q5. If doing this, add 28µL of MM to your 25µL sample.



**NEBNext Ultra II Q5 Master Mix - 250 rxns**  
by [New England Biolabs](#)  
Catalog #: [M0544L](#)



**USER Enzyme - 250 units**  
by [New England Biolabs](#)  
Catalog #: [M5505L](#)

Component	1X
Ligated sample + TruSeq indexing primers from part VI	25 µL
NEBNext Ultra II 2X Q5 PCR Master Mix	25 µL
NEB USER Enzyme	3 µL
<i>Total</i>	<i>53</i> <i>µL</i>

User/Q5


- 49 Set up the following cycling conditions for USER enzyme cutting and indexing PCR in a post-PCR room:

Temperature	Time	Cycles
37°C	15mins	1
98°C	30 sec	1
98°C	10 sec	12 cycles
65°C	75 sec	
65°C	5 min	1
4°C	-	-












Set the lid heat to **ON** at 105°C.

SPRI clean-up 0.9X

- 50 Equilibrate clean SPRI beads to room temperature and vortex well to mix.



**SPRI beads (homemade) or Ampure XP beads**  
[View](#)

- 51 Add beads equivalent to 1:0.9x the sample volume to each sample tube (for  53 µl of sample, add  47.7 µl beads).
- 52 Mix well by pipetting or tapping the tubes. Pulse-spin in a picofuge for no more than 2 seconds.
- 53 Incubate for  00:05:00 at room temperature, then put the tubes on the magnetic rack. Allow beads to separate on the magnet for 3-5 minutes, or until the supernatant is clear.
- 54 Keeping the tubes on the magnet, carefully remove and discard the supernatant.
- 55 Add  200 µl  80 % volume ethanol (prepared fresh). Incubate beads for  00:01:00 and then remove the ethanol
-  80% Ethanol
- 56 Repeat the above ethanol wash step.
- 57 Allow the beads to air dry for  00:05:00 . Do not overdry. Dry beads should appear matte (rather than glossy), but should not have a cracked appearance. Overdried beads may not resuspend or elute well.
- 58 Remove tubes from magnet and resuspend in  53 µl nuclease-free H<sub>2</sub>O.
- 59 Resuspend well by tapping the tubes and spin down briefly in a picofuge.
- 60 Allow  00:02:00 for DNA to elute from beads, then transfer tubes back to magnet
- 61 Allow the beads to separate for at least 2 minutes.
- 62 Collect  23 µl of supernatant to clean PCR tubes.

- 63 Using optical PCR strip tubes which are separated from each other so that they can be removed from the thermocycler one at a time, add KAPA amplification MM and Illumina P7 and P5 (5sol-20 and 5sol-21) primers at 5µM to your samples as below. You can make a master mix of KAPA master mix and primers.



#### Kapa HiFi Real-Time Amplification Kit

by Kapa Biosystems

Catalog #: [KK2702](#)



#### Illumina P5 and P7 primers 5uM combined;


P5: 5' AATGATACGGCGACCACCGAGATCT

P7: 5' CAAGCAGAAGACGGCATACGAGAT

[View](#)

Component	1X
Amplified and indexed DNA	23 µL
Kapa amplification master mix	25 µL
Illumina P5 and P7 5uM primers	2 µL
<i>Total</i>	<i>50</i> <i>µL</i>

KAPA amplification

- 64 Add  50 µl of STD 2 to a clean optical PCR tube.
- 65 Cap all tubes with optical caps. Do not write on the caps.
- 66 Place your labeled samples in the RT-PCR thermocycler.

## 67 Set up the following PCR conditions:

Temperature	Time	Cycles
98°C	45 sec	1
98°C	15 sec	20
60°C	30 sec	
72°C	1 min 30 sec	
Plate read		
72°C	30 sec	

Thermocycling Conditions for KAPA Amplification

Run the program, and watch until your sample either:

1. Crosses the standard (STD) 2 threshold
2. Starts to plateau

Then pull your sample out **DURING THE 72C 30sec INCUBATION** after the plate read. It is critical to pull it out during this step, and not when it is denaturing or annealing.

If your sample still has not reached the STD 2 by 20 cycles, let the program finish.

Alternatively, if you have many samples, you can simply apply 10 cycles to all, and evaluate if you need further amplification later.

You may want to take note of the unique number of cycles each sample needed to reach the STD 2 threshold.

### SPRI Clean-up 1X

## 68 Equilibrate clean SPRI beads to room temperature and vortex well to mix.






SPRI beads (homemade) or Ampure XP beads  
[View](#)

## 69 Add beads equivalent to 1X the sample volume to each sample tube (for 50 µl of sample, add 50 µl beads).

## 70 Mix well by pipetting or tapping the tubes. Pulse-spin in a microfuge for no more than 2 seconds.


## 71 Incubate for 00:05:00 at room temperature, then put the tubes on the magnetic rack. Allow beads to separate on the magnet for 3-5 minutes, or until the supernatant is clear.


## 72 Keeping the tubes on the magnet, carefully remove and discard the supernatant.

73 Add  200 µl  80 % volume ethanol (prepared fresh). Incubate beads for  00:01:00 and then remove the ethanol




74 Repeat the above ethanol wash step.

75 Allow the beads to air dry for  00:05:00 . Do not overdry. Dry beads should appear matte (rather than glossy) but should not have a cracked appearance. Overdried beads may not resuspend or elute well.

76 Remove tubes from magnet and resuspend in  27 µl nuclease-free H<sub>2</sub>O.

77 Resuspend well by tapping the tubes and spin down briefly in a picofuge.


78 Allow  00:02:00 for DNA to elute from beads, then transfer tubes back to magnet

79 Allow the beads to separate for at least 2 minutes.


80 Collect  25 µl of supernatant to clean PCR tubes.

#### Library Analysis

81 Quantify by HS DNA Qubit

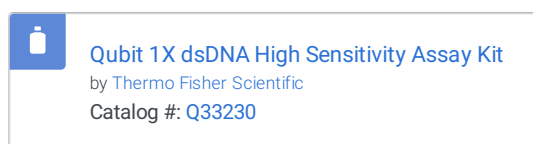

- 82 Run a fragment analysis, such as with the HS DNA Bioanalyzer. You are expecting a fragment trace that is characterized by sharp peaks in the 250-650 range. A successful FLASHed sample trace looks approximately like this: (lower marker in green, upper marker in purple)



**Bioanalyzer chips and reagents (DNA High Sensitivity kit)**  
by Agilent Technologies

