

In vitro transcription of guide RNAs version 8

Mark Dewitt, Julia Wong

Abstract

sgRNA template assembly, in vitro T7 transcription, and SPRI bead cleanup

Citation: Mark Dewitt, Julia Wong In vitro transcription of guide RNAs. **protocols.io**

dx.doi.org/10.17504/protocols.io.hdrb256

Published: 17 Mar 2017

Guidelines

The primers used are: one long, variable oligo that carries the T7 promoter and desired guide sequence; an 82-nt constant oligo that carries the 3' end of the sgRNA; two short external primers for amplification.

Assembly Oligos:

T7FwdVar:

GGATCCTAATACGACTCACTATAG---guide-sequence---GTTTGTAGAGCTAGAA


T7RevLong:

AAAAAAGCACCGACTCGGTGCCACTTTTTCAAGTTGATAACGGACTAGCCTTATTTAACTTGCTATTTCTAGCTCTAAAC

Amplification Primers:

T7FwdAmp: GGATCCTAATACGACTCACTATAG

T7RevAmp: AAAAAAGCACCGACTCGG

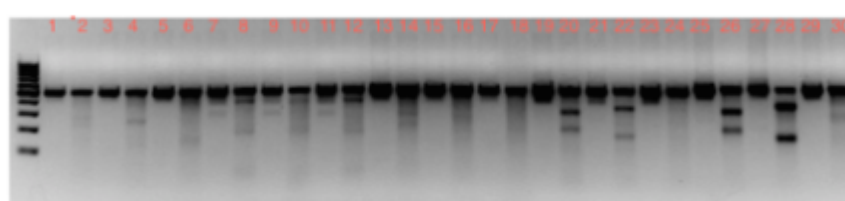

Odd numbers: just PCR product

Even numbers: with T7E1 treatment

Here's a gel that shows some good guides and some that didn't work| examples of more promising looking guides would be #'s 20, 22, 26, and 28.

Designing the guide sequence: The T7 transcript starts with the G before the dashes in the Fwd assembly oligo. If your target calls for a G at the 5' end of the guide, use the one already included in T7FwdVar design. If your target calls for a different base at the 5' end, put the whole guide sequence in T7FwdVar; there will be an extra G on the end that should have minimal effect, unless your guide is quite short.

Making the template: we have used the Phusion-HF DNA polymerase. Although the amplification primers have quite different Tm's, they work together fine.

Protocol

Design sgRNA and order PCR oligos.

Step 1.

Add the desired protospacer sequence to the T7FwdVar oligo and order the oligo from your favorite oligonucleotide supplier. There are many programs available for protospacer design that attempt to optimize on- and/or off-target activity. Which program is most useful depends upon many factors including type of editing, organism being edited, etc. Choice of protospacer design program is beyond the scope of this protocol.

The transcript will start with the bolded G just 5' of the dashes in the T7FwdVar oligo. T7 RNA polymerase requires a 5' G for proper transcript initiation. If your protospacer has a G at the 5' end, you can omit it from the T7FwdVar design to avoid duplication of the G. If your protospacer has a C, T, or A at the 5' end, add the whole protospacer sequence to T7FwdVar. In this case, there will be an extra G added to the 5' end of the protospacer, but literature indicates this will have minimal effect unless your guide is very short.

Primers:

T7FwdVar oligo (5'-GGATCCTAATACGACTCACTATAG^{**G**}--protospacer sequence—GTTT TAGAGCTAGAA-3')
T7RevLong oligo (5'-
AAAAAAGCACCGACTCGGTGCCACTTTTTCAAGTTGATAACGGACTAGCCTTATTTAACTTGCTATTTCTAGCTC
TAAAAC-3')
T7FwdAmp primer (5'-GGATCCTAATACGACTCACTATAG-3')
T7RevAmp primer (5'-AAAAAAGCACCGACTCGG-3')

Making in vitro transcription DNA template

Step 2.

For each T7FwdVar oligo you designed, set up the following PCR (total volume should be 20.0 µL). Make sure everything is RNase free and filter tips are used. Furthermore, wipe down everything (in every step of protocol) with RNase Away to ensure no contamination.

13.4 µl DEPC-treated H₂O
4.0 µL 5x Phusion HF Buffer
0.8 µl 10 mM dNTPs
0.4 µl T7FwdVar (1 µM)
0.4 µl T7RevLong (1 µM)
0.4 µl T7FwdAmp (100 µM)
0.4 µl T7RevAmp (100 µM)
0.2 µl Phusion HF DNA polymerase (2u/µl)

If making multiple sgRNA templates, prepare a master mix with all components except T7FwdVar. Include a no template control (omit T7FwdVar).

Primers:

T7FwdVar oligo (5'-GGATCCTAATACGACTCACTATAG--protospacer sequence—GTTT TAGAGCTAGAA-3')

T7RevLong oligo (5'-

AAAAAAGCACCGACTCGGTGCCACTTTTTCAAGTTGATAACGGACTAGCCTATTTTAACTTGCTATTTCTAGCTC
TAAAC-3')

T7FwdAmp primer (5'-GGATCCTAATACGACTCACTATAG-3')

T7RevAmp primer (5'-AAAAAAGCACCGACTCGG-3')

Making in vitro transcription DNA template

Step 3.

Run PCR:

95° 30 sec

95° 10 sec

57° 10 sec

72° 10 sec

34x steps 2-4

72° 2 min

4° hold

No PCR cleanup necessary at this point

In vitro T7 transcription

Step 4.

Mix the following to make **20 µl** total T7 transcription mix

volume	reagent
2 µl	10x Buffer 1x
2 µl	ATP (100 mM), 10 mM
2 µl	GTP (100 mM), 10 mM
2 µl	CTP (100 mM), 10 mM
2 µl	UTP (100 mM), 10 mM
8 µl	DNA template
2 µl	T7 RNA polymerase mix

In vitro T7 transcription

Step 5.

Incubate transcription mix for 18 hours at 37° in a thermocycler

 **DURATION**

18:00:00

Step 6.

Add 1 µl of RNase-free DNase; incubated 20 min, room Temp

 DURATION

00:20:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 7.

Bring volume to 150 µL with 100% EtOH (this helps binding of small fragments)

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 8.

Add 5X SPRI ([we use homemade SeraPure beads for RNA binding](#))

5*10 (IVT sgRNA)= 50 µL of SPRI Beads

5*20 (IVT sgRNA)= 100 µL SPRI Beads

 REAGENTS

Agencourt AMPure XP [A63880](#) by [Beckman Coulter](#)

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 9.

Pipette to mix 10 times

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 10.

Incubate 5 minutes at room temperature

 DURATION

00:05:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 11.

Place on magnetic stand, 5 min

 DURATION

00:05:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 12.

Discard supernatant

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 13.

Wash#1 Add 200 µL, 80% EtOH. Wait 2 min. Remove EtOH.

 DURATION

00:02:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 14.

Wash #2: Add 200 µL, 80% EtOH. Wait 2 min. Remove EtOH.

 DURATION

00:02:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 15.

Air dry 5-10 min (pellet will change from a glossy/wet to matte/dry.)

 DURATION

00:10:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 16.

Elute 20 µL of water or TE. Pipette mix 10 times.

Step 17.

Incubate 2 minutes at room temperature

 DURATION

00:02:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 18.

Place on magnetic stand, 5 min

 DURATION

00:05:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 19.

Keep Supernatant. Transfer to a new plate / tubes.