



May 10,  
2019

Working

## U Mass - Non-esterified fatty acids [↗](#)

Jason Kim<sup>1</sup>

<sup>1</sup>University of Massachusetts

[dx.doi.org/10.17504/protocols.io.x3hfqj6](https://doi.org/10.17504/protocols.io.x3hfqj6)

Mouse Metabolic Phenotyping Centers

Tech. support email: [info@mmpc.org](mailto:info@mmpc.org)

Lili Liang

### ABSTRACT

#### Summary:

This experiment measures serum and plasma concentrations of non-esterified fatty acids using a 96-well kit. The experiment involves a coupled reaction to measure non-esterified fatty acids (NEFA) which ultimately forms a purple product that absorbs light at 550nm. This allows the concentration of NEFA to be determined from the optical density measured at 540~550nm. Serum fatty acids levels reflect systemic lipid metabolism, lipid digestion/absorption, and lipid clearance. Serum fatty acids levels are altered in obesity, insulin resistance, and type 2 diabetes.

### EXTERNAL LINK

<http://mmpc.org/shared/document.aspx?id=168&docType=Protocol>

### MATERIALS

NAME	CATALOG #	VENDOR	CAS NUMBER	RRID
96-well assay plate blank	SFA-1	Zen-Bio		
Dilution Buffer	SFA-1	Zen-Bio		
FFA Standard	SFA-1	Zen-Bio		
FFA Diluent A	SFA-1	Zen-Bio		
FFA Diluent B	SFA-1	Zen-Bio		
FFA Reagent A	SFA-1	Zen-Bio		
FFA Reagent B	SFA-1	Zen-Bio		
Multichannel Pipette Tray	SFA-1	Zen-Bio		

### MATERIALS TEXT

#### Additional Items

- Multi-channel Pipet , single channel pipet and pipet tips
- Plate reader with a filter of 540 nm
- Incubator at 37°C
- Tubes for dilution of standards


#### Reagent Preparation:

##### Reagent 1:

Preparation of standard curve using the Standard Solution:

1. Briefly spin down the contents of the FFA standard tube before reconstitution.
2. Standard FFA concentrations are 0, 1.4, 4.1, 12.3, 37, 111, and 333  $\mu$ M.

3. The kit standard solution is the 1.0 mM standard concentration.
4. Pipette 120  $\mu\text{l}$  of Dilution Buffer into 6 tubes.
5. Pipette 60  $\mu\text{l}$  of the FFA Standard Stock solution into a tube labeled 333  $\mu\text{M}$ .
6. Prepare a dilution series as depicted below.
7. Mix each new dilution thoroughly before proceeding to the next solution.
8. The Dilution Buffer alone serves as the zero standard solution.

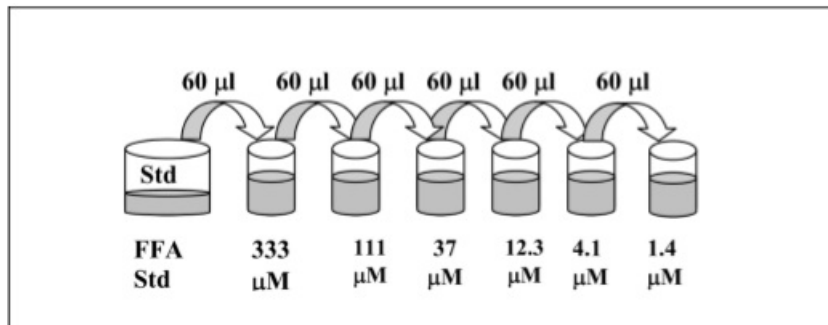

#### Reagent 2:

Preparation of FFA Reagent A:

1. Add 10.5 ml FFA Diluent A per bottle, and gently invert. Do not vortex.
2. Store any remaining solution at 2~8°C. The reagent solution is stable for 10 days after reconstitution when refrigerated at 2~8°C.

#### Reagent 3:

Preparation of FFA Reagent A:

1. Add 5.5 ml FFA Diluent B per bottle, and gently invert. Do not vortex.
2. Store any remaining solution at 2~8°C. The reagent solution is stable for 10 days after reconstitution when refrigerated at 2~8°C.

#### BEFORE STARTING

##### Notes:

- ✓ Freshly prepared blood or plasma samples are recommended. If storing samples, keep blood and plasma samples at -20°C or at -70°C for long-term storage. Avoid freeze/thaw cycles.
- ✓ Avoid using samples with gross hemolysis or lipemia.
- ✓ Allow all reagents to come to room temperature before measurement.

- 1 Add 5  $\mu\text{l}$  (or 1~10  $\mu\text{l}$ ) of serum or plasma sample to a well of Plate A.
- 2 Add dilution buffer to each well to reach a total sample volume of 50  $\mu\text{l}$ .

Addition of 5  $\mu\text{l}$  results in a 10x dilution of sample (5  $\mu\text{l}$  of serum/plasma sample in 50  $\mu\text{l}$  total sample volume).

3

4 Add 50 µl of each standard to empty wells. Use Plate B if necessary.

5 Add 10.5 ml of the reconstituted FFA Reagent A to one of the disposable trays provided with the kit.

6 Add 100 µl of FFA Reagent A to each well.

7 Gently shake the plate to ensure thorough mixing.

8 Place in a 37°C incubator for 10 minutes.

9 Add 5.5 ml of the reconstituted FFA Reagent B to the other disposable tray provided with the kit.

10 Add 50 µl of FFA Reagent B to each well.

11 Gently shake the plate to ensure thorough mixing.

12 Place in a 37°C incubator for 10 minutes.

13 Allow the plate to equilibrate to room temperature for 5 minutes. During this time, ensure that there are no bubbles in the solution mixture.

14 Use a large gauge needle or clean pipet tip to pop any bubbles as this will result in inaccurate absorbance readings.

15 Measure the optical density of each well at 540 nm.



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited