



Oct 20, 2019

cDNA synthesis

[iGEM Dusseldorf](#)¹¹Heinrich-Heine Universität Düsseldorf

1

Works for me

[dx.doi.org/10.17504/protocols.io.8hhht36](https://doi.org/10.17504/protocols.io.8hhht36)


iGEM Dusseldorf ⚡

MATERIALS TEXT

200 ng of RNA

preferred cDNA synthesis Kit

- 1 for cDNA synthesis use synthesis Kit for e.g. the BioRad iScript cDNA synthesis Kit (Cat.# 1708891)
- 2 use 200 ng of RNA for reaction
- 3 the reaction protocol reverse transcription was carried out at 42° C for 30 min
- 4 The reverse transcriptase inactivation was performed at 85° C for 5 min
- 5 used directly or stored at -20° C

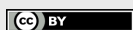

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited