



Mar 21,
2019

In devel.

Second-derivative UV Spectroscopy

Moriah Beck¹

¹Wichita State University

[dx.doi.org/10.17504/protocols.io.q7qdzmw](https://doi.org/10.17504/protocols.io.q7qdzmw)

Beck Lab



Moriah Beck

Wichita State University



ABSTRACT

Second-derivative spectroscopy has proven to be an effective analytical tool because of its ability to resolve overlapping bands in the normal spectrum.

This protocol is used to determine the degree of tyrosine exposure in proteins using second-derivative UV spectroscopy.

THIS PROTOCOL ACCOMPANIES THE FOLLOWING PUBLICATION

Ragone, R, Colonna, G, Balestrieri, C, Servillo, L, and Irace, G
Biochemistry (1984) 23: 1871-75

TyrExposure_2ndDerUV.p
df

2ndDer_struct_Dynamics.
pdf

Second-
derivative_spectroscopy_of
_protei.pdf

PROTOCOL STATUS

In development

We are still developing and optimizing this protocol

SAFETY WARNINGS

Protein Concentrations

- 1 Determine protein concentrations either using absorption coefficients (Abs @ 280 nm) or using BCA, Bradford, etc. colorimetric assays.


Proteins should be used at concentration where $A_{280} \sim 0.2$

UV Scan

- 2 Scan from 270-300 nm

Obtain Second-derivative

- 3 Copy UV scan data from instrument into a Prism file with columns for wavelength and absorbance values.

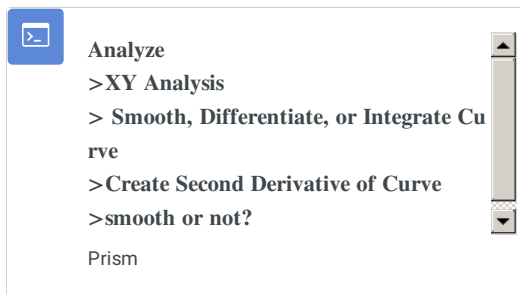

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited