OPEN ACCESS



Nuclear DNA purification from recalcitrant plant species for long-read sequencing

Ashley Jones, Justin Borevitz

Abstract

Evolution has driven genetic diversity of life on Earth, but also created highly complex genomes that are difficult to sequence. Current draft genomes can have thousands to hundreds of thousands of contigs rather than chromosomes, containing incorrect assemblies, gaps and errors. With rapid advances in long-read technologies, it is becoming possible to resolve complex genomes, including repetitive, polyploid plant genomes. Despite the technology being available, a challenge persists: the extraction of pure high molecular weight DNA suitable for long-read sequencing. This is particularly true of recalcitrant native Australian trees such as Eucalypts and Acacias. To resolve this, firstly we optimised a density gradient and detergent based nuclei extraction to limit reads from high copy count plastid genomes. Secondly, we optimised a gentle high molecular weight DNA extraction free of columns and high centrifugation, to limit DNA fragmentation. Finally, the DNA was purified and size selected by gel electrophoresis. For sequencing, we adopted the portable MinION sequencer from Oxford Nanopore Technologies. Using these approaches, we have been obtaining >9 gigabases of sequencing from a single MinION flowcell, including reads over 100 kb in length. Such ultra-long reads assist the assembly of high quality genomes, from telomere to telomere.

Citation: Ashley Jones, Justin Borevitz Nuclear DNA purification from recalcitrant plant species for long-read sequencing. **protocols.io**

dx.doi.org/10.17504/protocols.io.ssveee6

Published: 30 Aug 2018

Guidelines

This protocol is a modified and optimised combination of the following two publications. When citing, please also note the original publications below.

Bolger *et al.* (2014). The genome of the stress-tolerant wild tomato species Solanum pennellii. Nature Genetics **46,** 1034.

Mayjonade *et al.* (2016). Extraction of high-molecular-weight genomic DNA for long-read sequencing of single molecules. BioTechniques **61**, 203-205.

Eucalyptus tissue used in this protocol was kindly provided by the Australian National Botanic

Before start

This protocol is designed for two sample nuclei isolations and DNA extractions, although upscaling is possible. Prepare the following for two samples:

- Autoclave x2 500 mL Nalgene bottles, x10 1 L Schott bottles, x2 beakers and x2 funnels.
- Collect a branch of Eucalyptus leaves, harvest and wash leaves with tap water.
- Weigh 50 g of leaves into 50 g takeaway containers.
- Pre-chill the blender in 4°C cold room.
- For two sample extractions, on the day prepare 1 L nuclei isolation buffer and 1 L of nuclei wash buffer, and cool both down to 4°C.
- Prepare lysis buffer fresh on the day of use.
- Prepare all other solutions as listed below, storing at room temperature until use.

Nuclei isolation buffer (1 L for two samples)

- Adjust to pH 6 (add HCl; estimate listed below).
- Cool down to 4°C

Component MW	Stock	Quantity (1 L)
1 M Hexylene gylcol 118.17 (2-Methyl-2,4-pentanediol)	7.51 M (shipped liquid)	133 mL
10 mM PIPES 302.4	0.5 M	20 mL
10 mM MgCl ₂ 203.3	1 M	10 mL
10 mM Sodium metabisulfite 190.11	10%= 0.53 M	20 mL
6 mM ethylene glycol tetraacetic acid (EGTA) pH 7	0.5 M	12 mL
pH 7 equilibration with HCl 36.46	5 M	4 mL
0.5% (w/v) sodium diethyldithiocarbamate 171.26	10%	50 mL
4% (w/v) PVP-10 10,000	powder	20 g
200 mM L-lysine 146.19	powder	36.53 g
1 mM Dithiothreitol (DTT) 154.25	powder	0.154 g

Nuclei wash buffer (1 L for two samples)

- Adjust to pH 7 (add NaOH; estimate listed below).
- Cool down to 4°C

Component	MW	Stock	Quantity (1 L)
0.5 M Hexylene gylcol (2-Methyl-2,4-pentanediol)	118.17	7.51 M (shipped liquid)	66.66 mL
10 mM PIPES	302.4	0.5 M	20 mL
10 mM MgCl ₂	203.3	1 M	10 mL
10 mM Sodium metabisulfite	190.11	10%= 0.53 M	20 mL
6 mM ethylene glycol tetraacetic acid (EGTA) pH 7	380.4	0.5 M	12 mL
pH 7 equilibration with NaOH	40.00	5 M	3 mL
0.5% (v/v) Triton X-100	647	10%	50 mL
200 mM L-lysine	182.65	powder	36.53
1 mM Dithiothreitol (DTT)	154.25	powder	0.154 g

Lysis buffer (10 mL per sample)

• Prepare a fresh solution for optimal results. The solution should be clear before use.

Component	[Stock]	Quantity (10 mL)
1% Polyvinylpyrrolidone 40 (PVP-40)	10%	1 mL
1% Sodium metabisulfite	10%	1 mL
0.5 M NaCl	5 M	1 mL
100 mM TRIS-HCl pH 8.0	1 M	1 mL
5 mM EDTA pH 8.0	0.5 M	1 mL
3% Sodium dodecyl sulfate (SDS)	20%	1.5 mL
Water	-	3.5 mL

2% Sera-Mag beads solution (10 mL stock)

- Sera-Mag SpeedBead magnetic carboxylate modified particles (Thermo Scientific 65152105050250).
- Store at 4°C for up to 6 months.
- First prepare buffer without the beads. Let the Sera-Mag beads come to room temperature.

Component	[Stock]	Quantity (10 mL)
18% Polyethylene glycol (PEG) 8,000	25%	7,200 μL
1 M NaCl	5 M	2,000 μL
10 mM TRIS-HCl pH 8.0	1 M	100 μL
1 mM EDTA pH 8.0	0.5 M	20 μL
0.05% Tween 20	10%	50 μL
Water	-	430 μL

- Homogenise Sera-Mag beads thoroughly by shaking and swirling.
- Prepare 2% Sera-Mag beads (200 μ L for 10 mL buffer) by washing 4 times with water to remove sodium azide. Magnetise, remove supernatant, add ~500 μ L H₂O, flick tube, magnetise etc.
- Resuspend the clean Sera-Mag beads in the buffer prepared.
- Note: Potentially, AMPURE XP (Beckman Coulter) can be used but they must be washed 4 times with water and resuspended in their initial buffer.

Binding buffer solution (10 mL per sample)

- 2 g of PEG 8,000 (or 8 mL of 25% stock) [final = 20%]
- 75 g of sodium chloride [final = 3 M]
- Adjust to 10 mL with molecular biology grade water (approx 2 mL).

Mix until the solution becomes clear. If PEG8000 is not dissolved, it can lead to a poor yield as PEG 8000 makes gDNA to bind to the beads.

Other solutions

- Triton X-100 (ideally 10% solution)
- Ethanol (100% and 70%)
- RNAse A (e.g. 20 mg/mL Invitrogen PureLink)
- 5 M Potassium acetate

Special equipment (for two samples)

- 500 mL Nalgene bottle (x2)
- 1 L Schott bottles (x10)
- Beakers (x2)
- Funnel (preferably glass) (x2)
- Large Fine hair paint brushes (x2) (e.g. J.Burrows Mop Goat Hair 24 #9341694325125)
- Magnetic Rack (for Falcon and Eppendorf tubes)
- 15 mL Falcon tubes
- 50 mL Falcon tubes
- DNA LoBind Eppendorf tubes (e.g. 022431021)
- Sieve or mesh, 300 μm (e.g. Kartell disk membrane; ART 844)
- Miracloth or stretcher sheets (e.g Drager 334201)
- 200 μL wide-bore pipette tips (e.g. Vertex 4290-00)

Protocol

BLENDER LYSIS (4°C cold room)

Step 1.

Place 50 g of leaves into a blender with enough nuclei isolation buffer to cover the blades (approximately 150 mL).

NOTES

Optimised for Eucalyptus, a diploid ~500 Mb genome. Larger genomes or polyploid plants can have less input accordingly.

Step 2.

Homogenise the tissue on max speed. Transfer to a beaker, repeatedly washing out blender with nuclei isolation buffer (500 mL buffer is allocated per extraction).

Step 3.

Clean the blender with tap water, 70% ethanol and then MQW; process the other sample (repeating steps 1 and 2).

FILTRATION (ideally 4°C cold room)

Step 4.

Filter homogenate using a sieve (or disk membrane mesh placed in a funnel) into a 1 L Schott bottle. Forcibly squeeze out as much residual homogenate from the leaf debris as possible, maximising nuclei capture.

Step 5.

Filter homogenate through 1 layer of Miracloth using a funnel and 1 L Schott bottle. Gently squeeze residual homogenate through the Miracloth.

NOTES

May need to change the layer half way through (very dirty and clogged).

Step 6.

Repeat Filtration through 2 layers of Miracloth then 4 layers (using new Schott bottles).

Step 7.

Finally, gravity filter through 8 layers of Miracloth (no squeezing) into a 500 mL Nalgene bottle. Repeat steps 4-6 for the second sample.

NUCLEI ISOLATION

Step 8.

Add 2.5 mL of 100% Triton X-100, or 25 mL of 10% (final concentration 0.5%).

Step 9.

Incubate mixture on an ice bath with gentle rocking for 30 min.

NOTES

Use this time to clean the mess made with the blender and filtering

Step 10.

Weigh the Nalgene bottles with contents and lids to ensure they are equal for balancing purposes.

Also ensure the Nalgene bottle has an intact o-ring to avoid leakage during centrifugation.

Step 11.

Centrifuge at 600 rcf and 4 °C for 20 minutes.

P NOTES

4,600 rpm on a Sorvall RC5C using rotor SLA-3000 (00).

Step 12.

Discard the supernatant and add 200 mL of nuclei wash buffer.

Step 13.

Gently resuspend the pellet using a large fine hair paint brush, soaked in freshly made pre-chilled nuclei wash buffer.

P NOTES

Ensure the bottles are balanced and have an o-ring as previously described.

Step 14.

Centrifuge at 600 rcf and 4°C for 20 minutes.

Step 15.

Discard supernatant and repeat with another 200 mL of nuclei wash buffer. The pellet should become grey-white with no traces of green. If not, make more nuclei wash buffer and repeat the resuspending and washing steps.

NOTES

1 L of nuclei wash buffer is suitable for x2 samples to be washed x2 times with 200 mL each, then a final resuspension with 100 mL each.

Step 16.

After the final spin down, discard the supernatant and resuspend the pellet with 50 mL of nuclei wash buffer. Transfer evenly across two 50 mL Falcon tubes. Repeat the resuspension with another 50 mL of buffer and add to the two Falcon tubes.

Step 17.

Centrifuge at 600 rcf and 4°C for 30 minutes.

Step 18.

Discard the supernatant. The nuclei pellets can now be stored -80°C.

DNA EXTRACTION FROM NUCLEI

Step 19.

Prepare 10 mL of fresh lysis buffer per sample. Ideally, pre-heat at 50-65°C.

NOTES

The solution should be clear before use.

Step 20.

To a 3-5 mL nuclei pellet, add 10 mL of lysis buffer.

NOTES

Therefore, the SDS will be approximately 2% final concentration. Usually just one of the two Falcon tubes of a nuclei prep per sample are used.

Step 21.

Add 200 µL of RNAse A to each sample (20 mg/mL Invitrogen PureLink).

Step 22.

Incubate the samples at 50-60°C for 1 h, shaking at 400 - 900 rpm if possible.

NOTES

Higher temperatures and longer incubations lead to DNA damage.

Step 23.

Add 1/3 volume of 5 M Potassium Acetate and mix by inverting to precipitate the proteins and the polysaccharides that will complex with SDS.

Step 24.

Incubate on ice (4°C) for 10 min (don't rotate, DNA vulnerable).

Step 25.

Centrifuge at 5,000 rcf for 5 min at 4°C.

Step 26.

Transfer supernatant to a new tube, centrifuge again at 5,000 rcf for 10 min at 4°C.

Step 27.

Transfer to a new tube and add 1.2x binding buffer. Ideally, split across x2 15 mL Falcon tubes.

Step 28.

Add 1 mL of 2% Sera-Mag beads (500 µL per Falcon tube if split).

NOTES

Beads are in excess and could be reduced.

Step 29.

Mix by inverting the tube 20 times. Incubate with gentle agitation using a rotator or a rocker platform for 15 min at room temperature.

Step 30.

Place the tube in a magnetic rack for 30 min or more (until the solution becomes clear, can be over an hour).

Step 31.

Remove the supernatant without disturbing the beads.

Step 32.

Wash the beads by filling the tube with 70% ethanol, let beads settle if disturbed, and pour out ethanol.

Step 33.

Repeat the ethanol wash another 3-5 times, until satisfied the beads are clean.

NOTES

If the beads are very dirty, remove the tube from the magnetic rack, resuspend beads by flicking the tube, magnetise until clear, remove supernatant, repeat washing.

Step 34.

Transfer the beads to an Eppendorf tube. This can be done in a series of small volume ethanol washes that resuspend the beads, transfer to Eppendorf tube, then place on the magnetic rack, remove supernatant and repeat into the same tube.

Step 35.

Remove all traces of ethanol, and let the beads air dry for 1-4 min. Do not let the beads dry

completely, they will crack and significantly reduce DNA recovery.

Step 36.

Using a wide-bore pipette tip, add 220 μ L of ultra-pure H_2O to the beads, gently resuspending. Gently tapping the tube is also suitable.

P NOTES

Larger volume based on PippinHT input across x1 whole cassette.

Step 37.

Place the tubes in the magnetic rack for 30 min. Let the solution become clear.

P NOTES

Highly concentrated DNA will take a long time. The tube can be left on the magnetic rack overnight, or add more water.

Step 38.

Transfer 200 µL of supernatant (contains DNA) to a new tube, avoiding any carry-over of beads.

GEL PURIFICATION: PIPPIN PREP

Step 39.

For some species, the DNA may be pure enough for sequencing at this point. However, for most recalcitrant plants, the DNA will still have impurities and DNA fragmentation is inevitable. Gel purification is an ideal solution to both problems.

Step 40.

Using a PippinHT (Sage Science) or similar automated electrophoresis product, gel purify approximately 30 μ g, following the manufacturer's instructions. A 15 kb high pass separation is recommended, however if DNA is plentiful, 20 kb high pass is more suitable. The PippinHT has 12 lanes, however a lane pair needs be dedicated to an external ladder, leaving 10 lanes for samples. 20 μ L of DNA goes into each lane (therefore 200 μ L elution in previous section). The manufacturer recommends a maximum 1.5 μ g per lane (15 μ g total per cassette), however, can be successfully overloaded to 3 μ g per lane (perhaps more). The size selection will not be as precise, but is of no concern here.

NOTES

Alternatively, a chloroform: isoamyl alcohol 24:1 clean-up can be performed to remove residual phenols and proteins. Other possibilities (untested), are digestion with proteinase K and another bead clean up. Also would be interesting to experiment with adding PVP to the binding buffer to further reduce polyphenols.

Step 41.

After separation, wait at least 45 min (hours or overnight is suitable), to aid elution and recovery.

Step 42.

Collect the contents of all elution wells into a DNA LoBind Eppendorf tube (approx. 300 µL).

Step 43.

Add 30 μ L of 0.1% tween solution to each elution well (provided in kit). Wait at least 10 min and then transfer the contents to the same DNA LoBind Eppendorf tube (another 300 μ L, tube total is approx. 600 μ L).

Step 44.

Add 1.2x binding buffer (approx. 720 μ L), and 100 μ L of 2% Sera-Mag beads to the DNA LoBind Eppendorf tube. Incubate, magnetise and ethanol wash as previously described in 'DNA extraction from nuclei'.

Step 45.

After drying the beads elute with 65 μ L of ultra-pure water. Incubate, magnetise and transfer to a new DNA LoBind Eppendorf tube using a wide-bore pipette tip. Expect 20-30% recovery relative to total input (6-9 μ g).

SEQUENCING LIBRARY PREPARATION

Step 46.


For sequencing, we adopted the portable MinION sequencer from Oxford Nanopore Technologies. There are two native genomic DNA library preparations available; a rapid transposase based method (SQK-RAD004) and ligation based method (SQK-LSK108 /109). Following the manufacturer's instructions, prepare a library and load the MinION for sequencing. Note that Oxford Nanopore recommends a mass of DNA optimised for 0.2 pmol (193 fmol). However, this is based on an average length of 8 kb. Therefore, the amount of input DNA needs to be adjusted. The following is recommended based on PippinHT size selection:

Fragment sizes	Rapid SQK-RAD004	Ligation SQK-LSK108 /109
No size selection, 8 kb	400 ng	1,000 ng
15 kb high pass	800 ng	3,000 ng
20 kb high pass	1,000 ng	4,000 ng

Using the protocol described, we have been obtaining clean high molecular weight DNA (Figure 1, Table 1). DNA fragment size has been predominantly 20-140 kb in length (Figure 2). During sequencing, we have been obtaining >9 gigabases from a single MinION flowcell for recalcitrant Eucalyptus species, including reads over 100 kb in length (Table 2, Figure 3).

(A) Eucalyptus melliodora

(B) Eucalyptus marginata

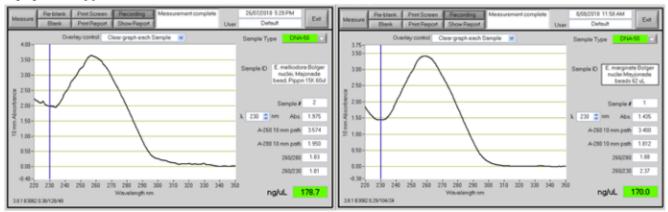

(C) Eucalyptus albens

(D) Eucalyptus sideroxylon

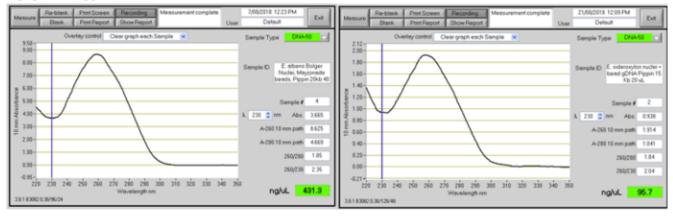


Figure 1: Spectrophotometer results of four recalcitrant *Eucalyptus* species. Readings taken using $1 \mu L$ on a Thermo Scientific Nanodrop 1000.

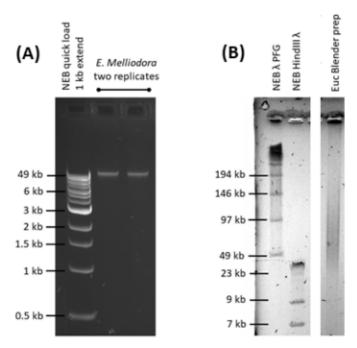


Figure 2: Gel electrophoresis analysis of DNA quality. *Eucalyptus melliodora* is shown as a representative example. (A) 50 ng separated on a 1% agarose gel. (B) 300 ng separated by pulsed field gel electrophoresis.

Table 1: DNA quantification and pippin prep. One of two nuclei pellets processed.

Sample	Approx DNA μg	Pippin input	Pippin setting	Qubit ng/μL	Nano ng/μL	260/ 280	260/ 230	Recovery µg
E. melliodora	36	18	15 kb	134	178.7	1.83	1.81	8.0 (44%)
E. marginata	29	29	20 kb	139	170.0	1.88	2.37	8.6 (30%)
E. albens	32	32	20 kb	259	431.3	1.85	2.35	8.8 (28%)
E. sideroxylon	7.6	7.6	15 kb	95.6	95.7	1.84	2.04	1.9 (25%)

Table 2: Sequencing results with a MinION flow cell per sample (FLO-MIN 106 R9).

Sample	DNA length	Library input ng	Library preparation	Pores	Output Gb	Average length	N50
E. melliodora	15 kb	800	Rapid SQK-RAD004	1,313	10.20	~10 kb	TBD
E. marginata	20 kb	1,000	Rapid SQK-RAD004	1,384	10.10	~10 kb	TBD
E. albens	20 kb	1,000	Rapid SQK-RAD004	1,208	9.26	~10 kb	TBD
E. sideroxylon	15 kb	1,000	Rapid SQK-RAD004	1,584	13.15	~8 kb	TBD

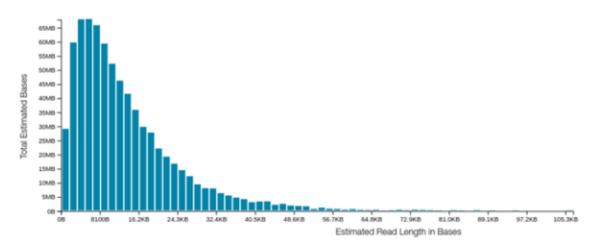


Figure 3: Expected read length histogram on MinKNOW 2.0. Here, *Eucalyptus melliodora* is shown, having an average read length of ~10 kb and reads exceeding 100 kb in length (low numbers but present).