



Sep 20, 2019

Multiplex Immunofluorescence on Fresh Frozen Tissue

Maya Brewer¹, Yuantee Zhu¹, Elizabeth Neumann², Danielle Gutierrez², Jeff Spraggins², Mark de Caestecker³

¹Vanderbilt University Medical Center, ²Vanderbilt University, ³Division of Nephrology, Vanderbilt University Medical Center

1 Works for me dx.doi.org/10.17504/protocols.io.665hhg6

VU Biomolecular Multimodal Imaging Center

Human BioMolecular Atlas Program (HuBMAP) Method Development Community

Maya Brewer

ABSTRACT

Scope:

To describe the procedure for multiple cycles of immunofluorescence on human kidney tissue embedded in carboxymethylcellulose.

Expected Outcome:

Kidney tissue sections that have been tagged with antibodies for imaging microscopy.


MATERIALS

NAME	CATALOG #	VENDOR
1X PBS (Phosphate-buffered saline)		
Phosphate Buffered Saline	28374	Thermo Fisher Scientific
1N NaOH		
Hydrophobic Barrier Pen	H-4000	Vector Laboratories
Glycine	410225	Sigma Aldrich
10X Power Block Universal Blocking Solution	HK085-5K	BioGenex
Antibody Diluent Reagent Solution	003218	Thermo Fisher Scientific
Hoechst 33342	62249	Thermo Fisher Scientific
100% Glycerol	G33	Fisher Scientific
30% Hydrogen Peroxide	216763	Sigma – Aldrich

MATERIALS TEXT

Equipment:

- Moisture/Humidified Chamber
 - 100-Slide Storage Box (Fisher Scientific, 03-448-1)
 - Kimwipes, 8.4 in x 4.4 in (Fisher Scientific, 06-666)
 - ddH₂O in a wash bottle
- LED Cabinet Light (Sears, SPM11582738325)
- Square Cell Culture Plate OR any container with clear bottom

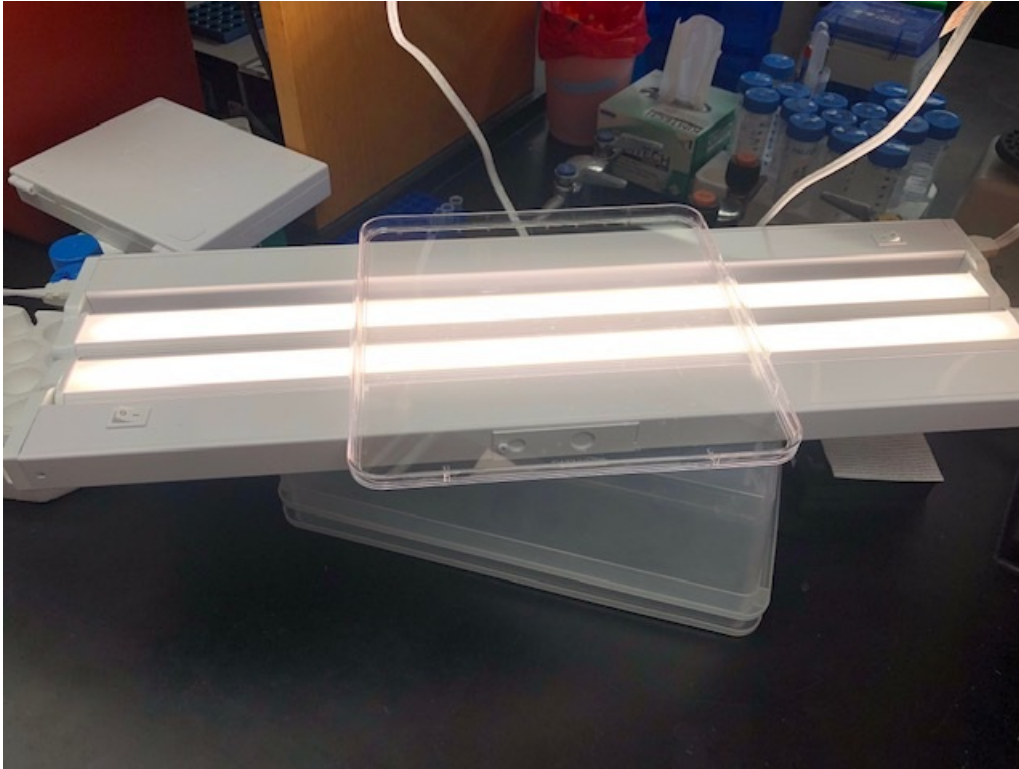


LED light system for fluorophore inactivation.

Immunofluorescence

- 1 If sections are frozen, allow them to equilibrate to room temperature (~15 minutes). Place slides in a humidified chamber. They will remain in the chamber throughout this protocol.
- 2 While slide is dry, using a hydrophobic pen, draw a large barrier around the section. **Do not allow pen to touch the section.**
- 3 Post-fix sections in 10% formalin for 5 minutes.

- 4 Remove fixative by tilting the slide, allowing the solution to flow from the section into the humidified chamber.

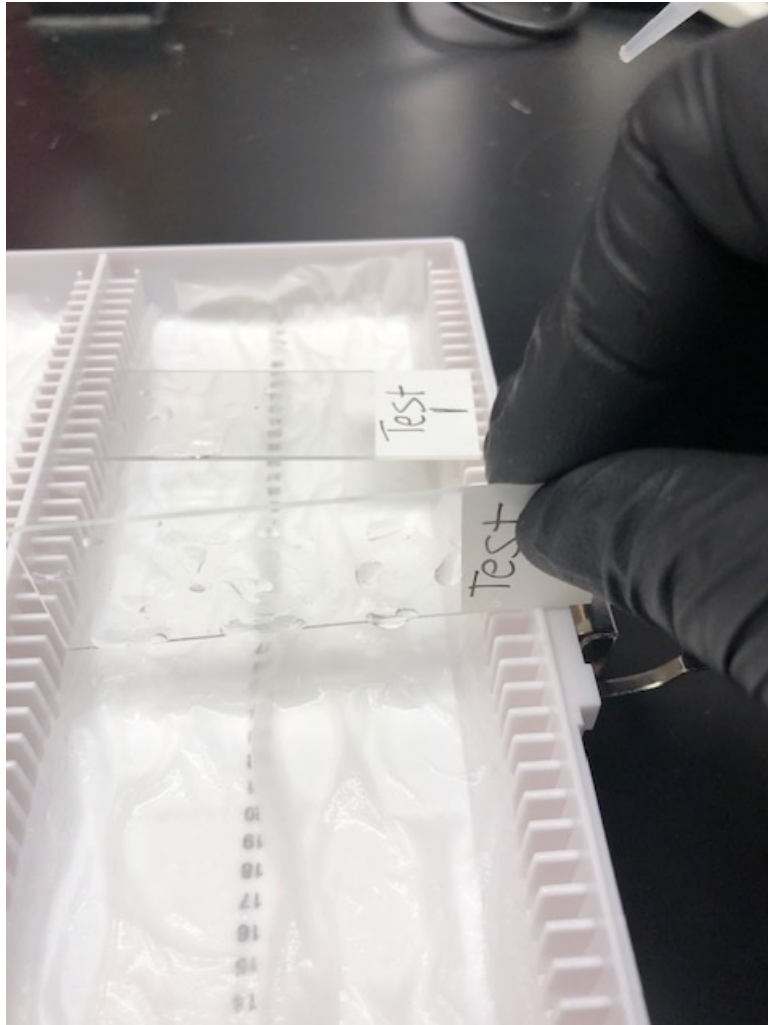

Removing solutions from slides in humidified chamber.

- 5 Wash sections in 1X PBS for 5 minutes three times. For this, tip solution off the slide into the chamber, add PBS to the slide using a pipettor, tip again and repeat.
- 6 Incubate sections for 5 minutes with 50 mM glycine (dilute stock in 1X PBS). This reduces autofluorescence by reducing free aldehyde groups.
- 7 Remove glycine, and wash sections in PBS for 5 minutes twice.
- 8 Incubate sections with 3% hydrogen peroxide (dilute 30% H_2O_2 in 1X PBS) at room temperature for 10 minutes to further reduce autofluorescence. This reduces non-specific fluorescence signals
- 9 Remove hydrogen peroxide, and wash sections in PBS for 5 minutes four times.
- 10 Block sections for 60 minutes with 1X Universal Blocking Reagent (UBR) at room temperature.

- 10.1 Dilute 10X blocking reagent to 1X using 9-parts ddH₂O to 1-part UBR.
- 11 Dilute primary antibody to desired working concentration in Antibody Diluent Reagent during blocking step.
- 12 Add diluted antibody to section and incubate overnight at 4°C.
- 13 Remove solution, and wash sections with PBS for 5 minutes twice.
- 14 If primary antibody is directly conjugated with a fluorophore, skip to #16.
- 15 If using indirect immunofluorescence, dilute fluorophore-conjugated secondary antibody using Antibody Diluent Reagent.
- 16 Add antibody solution to sections and incubate for 60 minutes at room temperature.
- 17 Remove solution, and wash sections with PBS for 5 minutes twice.
- 18 Incubate sections in Hoechst 33342 (1:10,000 dilution of 20mM solution in 1X PBS) for 10 minutes.
- 19 Remove Hoechst 33342, and wash sections with PBS for 5 minutes twice.
- 20 Mount slides in 70% glycerol in PBS. Do not seal the coverslips as they will need to be removed later.
- 20.1 Because slides are not sealed, they must be kept horizontal to prevent the coverslip from falling off, and in the humidified chamber to keep them from drying out.
- 21 Image
- 22 Store slides at 4°C in a moisture chamber.


De-Coverslipping

- 23 Remove coverslip from sections by incubating slides in a vertical staining jar filled with PBS for 15 minutes with slight agitation (i.e. plate rocker).

- 24 Slowly lift slide from vertical jar and allow coverslip to release from slide via gravity.
- 25 Wash slide in PBS for 5 minutes three times to remove any residual glycerol.
- 25.1 Place slides back in vertical staining jar full of PBS with slight agitation to wash.

Fluorophore I

- 26
Make a solution of 4.5% hydrogen peroxide and 24 mM sodium hydroxide made up in PBS.
- 27 Add bleaching solution to each section and incubate at room temperature for 90-120 minutes in the presence of white light (LED light).

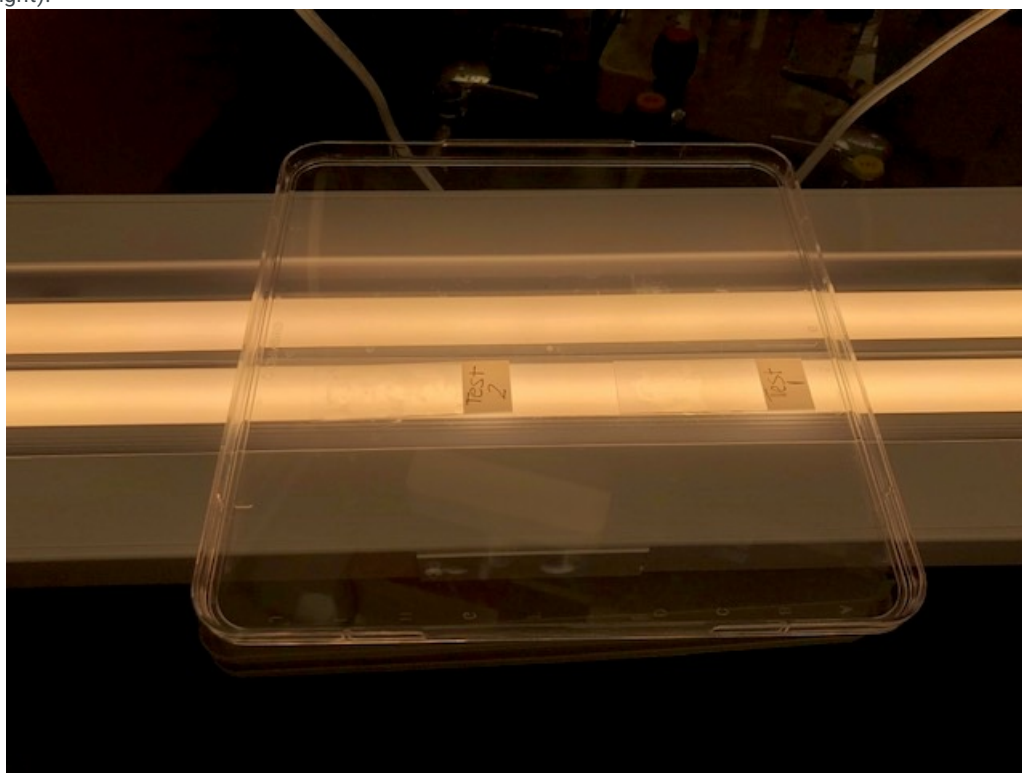


LED light system for fluorophore inactivation.

- 27.1 For this, place slides on a plastic surface on top of the LED light. We often use multiple plates for this.
- 27.2 Halfway through this incubation, the solution may be removed by pipette and replaced with fresh solution to ensure complete inactivation occurs from the LED light.
- 28 Remove solution and wash sections in PBS for 5 minutes four times.

- 28.1 After inactivating the fluorophores, slides are mounted in 70% glycerol and imaged to confirm complete fluorophore inactivation, followed by the removal of the coverslip, and three 5 minutes washes in PBS. **Hoechst stain will not bleach** and is necessary for image registration later.

Subsequent Immunofluorescence Cycles

- 29 Dilute fluorophore-conjugated primary antibodies using Antibody Diluent Reagent.
- 30 Add diluted antibodies to each section and incubate overnight at 4°C.
- 31 Remove antibody solution and wash sections with PBS for 5 minutes twice.
- 32 Incubate sections in Hoechst (1:10,000 dilution in 1X PBS) for 5 minutes (this may not be necessary, but does no harm).
- 33 Remove Hoechst, and wash sections with PBS for 5 minutes twice.
- 34 Mount slides in 70% glycerol in PBS.
- 35 Image
- 36 Store slides at 4°C in a moisture chamber.
- 37 **Sections II-IV are repeated for each remaining cycle.**

 This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited