

Dengue virus (DENV) universal MGB TaqMan 2017 version 4

Ian Mackay, Judy Northill

Abstract

This protocol was designed and developed at this laboratory.

The assay targets the capsid peptide coding region of DENV 1-4 and is designed as a qualitative screening test for human cases of DENV infection.

Citation: Ian Mackay, Judy Northill Dengue virus (DENV) universal MGB TaqMan 2017. **protocols.io**

dx.doi.org/10.17504/protocols.io.n7jdhkn

Published: 03 Apr 2018

Before start

- If using a different brand or model of real-time thermocycler, check the concentration of ROX is adequate.
- Method assumes the user is familiar with the thermocycler and software used to run the protocol.

Materials

SensiFAST™ Probe Lo-ROX Kit [BIO-84002](#) by [Bioline](#)

Protocol

Oligonucleotide sequences

Step 1.

| Name | Sequence (5'-3') |
|-------------|-----------------------------------|
| DU5-F1 | GAAYAACCAACG RAARA AGRCG |
| DU5-F2 | ATGAACCAACG RAARA AGGTGG |
| DU5-R13 | GAGAATCTCTTCGCCAACTGTG |
| DU5-R2 | TGAGAATCTCTT YGTCAR CTGYTG |
| DU5-R4 | GAGAATCTCTTCACCAACCCTTG |
| DU5-MGB2017 | 6FAM - AATATGCTGAAACGCG - MGBNFQ |

Reagents

Step 2.



REAGENTS

SuperScript™ III Platinum™ One-Step qRT-PCR Kit [11732088](#) by [Life Technologies](#)

Reaction set-up

Step 3.

The assay has been used on both a Rotor-Gene 6000 and a Rotor-Gene Q real-time thermocycler

Prepare sufficient mix for the number of reactions.

Include a suitable 'dead volume' as necessary if using a robotic dispenser.

MIX PREPARATION

| Reagent | Vol. (μl) x1 | Final reaction concentration |
|---|--------------|------------------------------|
| Nuclease-free water | 4.26 | N/A |
| DU5-F1 (200μM) | 0.09 | 900nM |
| DU5-F2 (200μM) | 0.09 | 900nM |
| DU5-R13 (200μM) | 0.03 | 300nM |
| DU5-R2 (200μM) | 0.03 | 300nM |
| DU5-R4 (200μM) | 0.03 | 300nM |
| DU5-DU5-F1 (100μM) | 0.03 | 150nM |
| 2X Reaction Mix ¹ | 10 | 1X |
| ROX reference dye (25μM) | 0.04 | 50nM |
| SuperScript™ III/Platinum™ Taq Mix ¹ | 0.4 | 1X |
| TOTAL | 15 | N/A |

¹Superscript™ III Platinum™ One-step qRT-PCR kit

- Dispense 15μL to each reaction well.
- Add 5μL of template (extracted RNA, controls or NTC [nuclease-free water]).
- Total reaction volume is 20μL

Amplification

Step 4.

| | | |
|------|------------------|-----|
| 50°C | 5min | 1X |
| 95°C | 2min | 1X |
| | | |
| 95°C | 3s | 40X |
| 60°C | 30s ¹ | |

¹Florescence acquisition step

Result analysis

Step 5.

The definition used for a satisfactory positive result from a real-time fluorogenic PCR should include each of the following:

1. A **sigmoidal curve** – the trace travels horizontally, curves upward, continues in an exponential rise and followed by a curve towards a horizontal plateau phase
2. A **suitable level of fluorescence** intensity as measured in comparison to a positive control (y-axis)
3. A **defined threshold (C_T) value** which the fluorescent curve has clearly exceeded (Fig.1 arrow), which sits early in the log-linear phase and is <40 cycles

A flat or non-sigmoidal curve or a curve that crosses the threshold with a $C_T >40$ cycles is considered a negative result. NTCs should not produce a curve

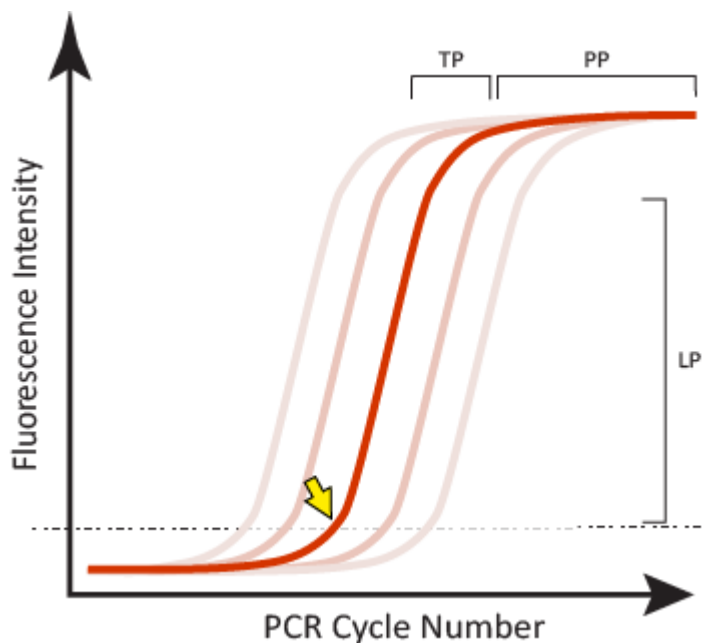


Figure 1. Examples of satisfactory sigmoidal amplification curve shape when considering an assay's fluorescent signal output. The crossing point or threshold cycle (C_T) is indicated (yellow arrow); it is the value at which fluorescence levels surpass a predefined (usually set during validation, or arbitrary) threshold level as shown in this normalized linear scale depiction. LP-log-linear phase of signal generated during the exponential part of the PCR amplification; TP-a slowing of the amplification and accompanying fluorescence signal marks the transition phase; PP-the plateau phase is reached when there is little or no increase in fluorescent signal despite continued cycling.