



Version 2

Oct 10, 2018

Working

# CRISPR Editing of Immortalized Cells with RNPs using Lipofection [↔](#) Version 2

Synthego<sup>1</sup>, Synthego<sup>1</sup><sup>1</sup>Synthego[dx.doi.org/10.17504/protocols.io.uibeuan](https://doi.org/10.17504/protocols.io.uibeuan)

Synthego

 Brittany Enzmann 

## ABSTRACT

This protocol describes how to deliver ribonucleoprotein (RNP) complexes that consist of purified Cas9 nuclease duplexed with chemically modified synthetic single guide RNA (sgRNA) to standard immortalized cell lines (adherent or suspension). Although optimized for HEK293 (human embryonic kidney 293 cells), this protocol is applicable to many other cell lines (A549, U2OS, HeLa, CHO, MCF-7). RNP delivery is accomplished using Lipofectamine™ CRISPRMAX™ Transfection Reagent. Chemically modified sgRNAs are designed to resist exonucleases and innate intracellular immune cascades that can lead to cell death. Synthego chemically modified synthetic sgRNAs are of exceptional quality and consistently drive high editing efficiencies. Data that demonstrate editing efficiency in HEK293 cells are included.

## TAGS

CRISPR

lipofection

Show tags

## EXTERNAL LINK

<https://www.synthego.com/resources/immortalized-cell-lipofection-protocol>Immortalized Cell  
Lipofection.pdf

## PROTOCOL STATUS

**Working**

Synthego uses this protocol and it is working

## GUIDELINES

This protocol is meant to serve as a starting point for lipofection of immortalized cells in a 24-well plate format. It may be necessary to experimentally optimize volumes and ratios for RNP formation for each cell type and for other culture plate formats. It is critical to add reagents in the order recommended in the steps. Prepare the RNP complexes with the Lipofectamine™ Cas9 Plus™ Reagent and Opti-MEM™ I Reduced Serum Medium in a separate tube (Tube 1) before adding diluted CRISPRMAX™ Reagent (Tube 2).

Reaction volumes are for EACH WELL and should be scaled up proportionally for the number of wells to be used.

*Abbreviations:*

CRISPR: clustered regularly interspaced short palindromic repeats Cas9: CRISPR associated protein 9

sgRNA: single guide RNA

RNP: ribonucleoprotein

PCR: polymerase chain reaction

ICE: inference of CRISPR edits

FACS: fluorescence-activated cell sorting TE: Tris EDTA

PBS: phosphate-buffered saline

GFP: green fluorescent protein

### Important Considerations

All Synthego and CRISPRMAX<sup>TM</sup> reagents should be stored according to the manufacturer's recommendations.

This protocol was optimized in HEK293 cells and can be used for other common cell lines such as A549, U2OS, HeLa, CHO, MCF-7.

Successful transfection is critically dependent on cell density. It may be necessary to optimize cell seeding densities in order to determine the most appropriate level of confluence for transfection.

For fast growing cells, seed fewer cells. Cell seeding is based on the rate of cell growth. Suggested starting cell numbers are listed in the protocol below.

In order to maximize CRISPR editing, be sure to include trypsinization ([Step 9](#)) in order to perform a reverse transfection of RNPs.

Use cells at lowest passage number possible.

Cas9 nuclease can be diluted in Opti-MEM<sup>TM</sup> I Reduced Serum Medium in order to achieve a working concentration according to the plate volume.

Synthego recommends sgRNA:Cas9 ratio of 1.3:1 for RNP formation. It may be necessary to optimize ratios for different cell lines/conditions.

RNP complexes are formed in Opti-MEM<sup>TM</sup> I Reduced Serum Medium and can be added directly to cells in culture medium irrespective of antibiotics. Following transfection, it is not necessary to remove RNP complexes or to add or change medium.

#### *Working with RNA and RNPs*

Wearing gloves and using nuclease-free tubes and reagents is recommended in order to avoid RNase contamination. • Always maintain sterile technique and use sterile filter pipette tips.

All Synthego reagents should be stored according to the manufacturer's recommendations.

Synthetic sgRNA should be dissolved in TE buffer and diluted to a working concentration using nuclease-free water. Please consult the [Synthego Quick Start Guide](#) for best practices related to dissolving and storing synthetic sgRNAs.

RNP complexes are stable at room temperature for up to 1 hour (may be stored at 4°C for up to one week, or at -20°C for up to 1 month). Note that RNPs stored at 4°C may become susceptible to contamination from microbial growth after long periods of time.

### Suggested Controls

Control	Description	Purpose
Mock	Cells transfected without Cas9 and sgRNA	Wild type sequence for comparison with experimental and other negative controls. Control toxicity from RNP, cell death from electroporation or possible viability issues associated with editing the specific gene of interest.
Negative	Cas9 complexed with a non-targeting sgRNA or no sgRNA	Ensure that there are no false positives due to contamination (no effect expected=wild type).
Positive	sgRNA with high editing efficiency (e.g., CDC42BPB, RELA)	Ensure all reagents, protocol, and equipment are functioning (effect expected).
Transfection	pMAX GFP (Lonza), GFP mRNA (SBI)	Assess transfection efficiency (without the use of RNPs).

## Timeline

Pre-Lipofection	Setup & Lipofection	Post-Lipofection		
Day 1 Seed cells Incubate (1 day)	Day 2 Assemble RNP Complexes Prepare Transfection Solution Make RNP-Transfection Solution Prepare Cell Suspension Transfect Cells Incubate (3 days)	Day 3	Day 4	Day 5 Analysis

**Note:** cell seeding may take 1-2 days and incubation after transfection may take 2-3 days.

## Additional Information















For an up-to-date list of all Synthego Protocols and other resources, please visit [synthego.com/resources](https://synthego.com/resources)

For technical assistance, contact our Scientific Support Team:

Ph: 888.611.6883

Email: [support@synthego.com](mailto:support@synthego.com)

## MATERIALS

NAME	CATALOG #	VENDOR
 PBS buffer	10010023	<a href="#">Thermo Fisher Scientific</a>
 Opti-MEM™ Reduced Serum Medium	31985062	<a href="#">Thermo Fisher Scientific</a>
 Chemically modified sgRNA	<a href="#">Chemically modified sgRNA</a>	<a href="#">Synthego</a>
 Positive control (optional); Recommended: human RELA sgRNA, CDC42BPB sgRNA		<a href="#">Synthego</a>
 Transfection control (optional); Recommended: pMAX GFP (Lonza), GFP mRNA (SBI)		Contributed by users
 TE buffer (Included with Synthego sgRNA)		<a href="#">Synthego</a>
 Nuclease-free water	R0581	<a href="#">Thermo Fisher Scientific</a>
 Cell counter		<a href="#">Thermo Fisher Scientific</a>
 Normal growth medium (Cell-type dependent)		Contributed by users
 Microcentrifuge tubes		<a href="#">Eppendorf</a>
 Tissue culture plates		<a href="#">Thermo Fisher Scientific</a>
 Trypsin		<a href="#">Thermo Fisher Scientific</a>
 Lipofectamine™ CRISPRMAX™ Cas9 Transfection Reagent (includes Cas9 Plus Reagent and CRISPRMAX Transfection Reagent)	CMAX00001	<a href="#">Thermo Fisher Scientific</a>
 Cas9 2NLS	<a href="#">Add at checkout</a>	<a href="#">Synthego</a>

## SAFETY WARNINGS

Please refer to the SDS (Safety Data Sheet) for safety warnings and hazard information.

### Pre-Lipofection - Seed cells

- Seed cells and incubate in 37°C/5% CO<sub>2</sub> incubator overnight so that they are 30-70% confluent on the day of transfection (Day 2).

 37 °C 5% CO<sub>2</sub> incubator

### Setup & Lipofection - Assemble RNP Complexes (1.3:1 sgRNA to Cas9 ratio)

- Dilute sgRNA and Cas9 to 3 µM stock concentrations (3 pmol/µl).

#### NOTE

If you followed the [Synthego QuickStart Guide](#), you may have 30 µM stock concentrations (30 pmol/µl), so be sure to dilute further to 3 µM (3 pmol/µl).

- 3 Prepare RNPs in microcentrifuge tube (Tube 1). Use the quantities (per reaction) in the table below.

RNP Preparation (Tube 1)		
Component	Molarity	Volume (per reaction)
Opti-MEM™ I Reduced Serum Medium	-	25 µl
sgRNA	3 µM (pmol/µl)	1.3 µl (3.9 pmol)
Cas9	3 µM (pmol/µl)	1 µl (3 pmol)
Lipofectamine™ Cas9 Plus Reagent	-	1 µl
Total volume	-	28.3 µl

**NOTE**

\* You may need to experimentally determine the optimum amounts of sgRNA and Cas9 nuclease. Synthego recommends a ratio of 1.3:1 sgRNA to Cas9 for RNP formation.

**NOTE**

Note: For knock-ins, see Synthego's [Nucleofection](#) or [Electroporation](#) protocol.

- 4 Incubate RNPs for 5-10 minutes at room temperature.

 00:05:00 Incubation

### Setup & Lipofection - Prepare Transfection Solution

- 5 In a separate microcentrifuge tube (Tube 2), dilute Lipofectamine™ CRISPRMAX™ Reagent in Opti-MEM™ I Reduced Serum Medium. Use the quantities (per reaction) in the table below.

Transfection Solution (Tube 2)	
Reagent	Volume (per reaction)
Opti-MEM™ I Reduced Serum Medium	25 µl
Lipofectamine™ CRISPRMAX™ Transfection Reagent	1.5 µl
Total volume	26.5 µl

 25 µl Opti-MEM I Reduced Serum Medium

 1.5 µl CRISPRMAX Reagent

- 6 Incubate transfection solution for 5 minutes at room temperature.

 00:05:00 Incubation

#### Setup & Lipofection - Make RNP-Transfection Solution



- 7 Add the transfection solution (Tube 2) directly to RNPs (Tube 1), and mix well by pipetting up and down.
- 8 Incubate for 5-10 minutes at room temperature. Do not exceed 30 minutes.

 00:05:00 Incubation

#### NOTE

*Synthego highly recommends reverse transfection (RNPs are added to wells first and cells are added second), as this method has resulted in high editing efficiencies.*

#### Setup & Lipofection - Prepare Cells

- 9
-  NOTE
- For suspension cells, resuspend in growth medium and mix well. Skip steps 9 and 10 below and proceed to step 11.*
- Wash cells with 1X PBS (enough to cover bottom of each well), then aspirate PBS.
- 10 Add trypsin (enough to cover bottom of each well), incubate for 5 minutes in a humidified 37°C/5% CO<sub>2</sub> incubator. Resuspend cells in an equivalent volume of medium to stop the trypsin reaction.
-  00:05:00 Incubation
- 11 Count cells to determine density.
- 12 Transfer  $0.42 - 1.2 \times 10^5$  cells per reaction to a microcentrifuge tube.
- 13 Centrifuge cells at 200 x g for 5 minutes.

 00:05:00 Centrifugation

- 14 Resuspend cells in 500 µl of the growth medium.

 500 µl growth medium

#### Setup & Lipofection - Transfect Cells

- 15 Add the RNP-transfection solution mixture to each well of a 24-well tissue culture plate (see table below).  
Add cell suspension to each well, and mix by pipetting (see table below).

RNP-Transfection Solution & Cell Suspension	
Reagent	Volume (per reaction)
RNP-Transfection Solution	50 µl
Cell suspension in growth medium	500 µl

 50 µl RNP-Transfection Solution

 500 µl Cell suspension in growth medium

- 16 Incubate cells for 2-3 days in a humidified 37°C/5% CO<sub>2</sub> incubator.

 37 °C Incubation

 48:00:00 Incubation

#### Post-Lipofection - Analysis


- 17 Extract DNA from cells.

- 18 Conduct analyses to determine editing efficiency: PCR, Sanger sequencing, and [ICE analysis](#). Next-Gen Sequencing, FACS, or functional tests may be conducted as alternatives.

#### NOTE

**Option:** If storing cells for future use is desired, split cells into two groups (one for analysis and one for cell culture).

## Representative Data

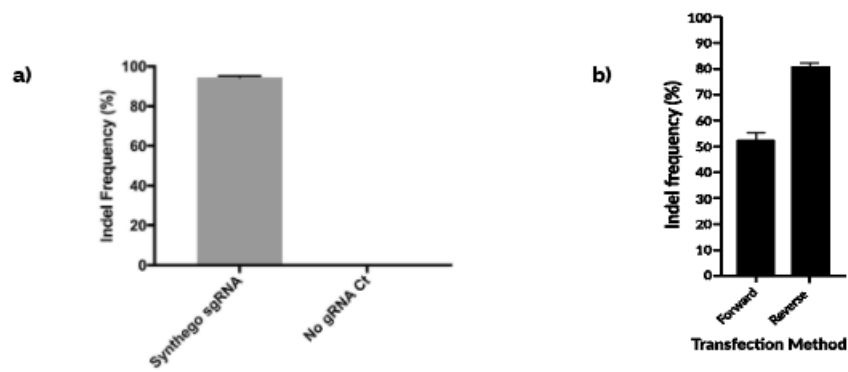

**Figure 1. a)** The percent indel frequency (mean  $\pm$ SD) of Synthego sgRNA targeted to CDC42BPB. RNPs were constructed with a sgRNA:Cas9 ratio of 3.91:3 and HEK293 cells ( $1 \times 10^5$ /reaction) were reverse-transfected using Lipofectamine™ CRISPRMAX™ Transfection Reagent. **b)** Indel frequency (mean  $\pm$ SD) in of reverse and forward-transfected HEK293 cells using Lipofectamine™ CRISPRMAX™). The indel frequency was 30% higher in reverse-transfected cells.



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited