



May 13,
2019

Working

U Michigan - Intra-Epidermal Fiber Density [↗](#)

Eva L. Feldman¹

¹University of Michigan - Ann Arbor

dx.doi.org/10.17504/protocols.io.x7kfrkw

Mouse Metabolic Phenotyping Centers
Tech. support email: info@mmpc.org

Lili Liang

ABSTRACT

Summary:

Intra-epidermal nerve fiber density (IENFD) is used as a tool to assess small fiber neuropathy. IENFD

EXTERNAL LINK

<https://mmpc.org/shared/document.aspx?id=319&docType=Protocol>

MATERIALS

NAME	CATALOG #	VENDOR	CAS NUMBER	RRID
Zamboni fix	# 1459A	Newcomer Supply		
cryomold	Cryomold	tissuetek		
sucrose	sucrose	Fisher Scientific		
PGP 9.5	14730-1-ap	Proteintech		AB_2210497
BSA		Sigma Aldrich		
Coverslip 1.5		Fisher Scientific		
slides		Fisher Scientific		
Prolong gold w/ dapi	P36931	Thermofisher		
Secondary 488 g anti mouse highly cross absorbed	A31620	Invitrogen - Thermo Fisher		

MATERIALS TEXT

Reagent Preparation:

Reagent 1:

Blocking Solution

Reagents and Materials

1) 5% BSA made in .3% triton tx100 in .1m PBS 7.4 pH

Reagent 2:

Primary antibody

Reagents and Materials

2) PGP antibody at 1/2000 made in 1% BSA in .3% triton tx100 in .1m PBS 7.4 pH

Reagent 3:

Secondary

Reagents and Materials

3) Molecular Probes anti-rabbit 488 highly cross absorbed secondary made in .3%

triton tx 100 in .1M PBS 7.4 pH.

Note:

Thermo Fisher Scientific, [RRID:SCR_008452](#)

Sigma-Aldrich [RRID:SCR_008988](#)

Proteintech Group, [RRID:SCR_008986](#)

PGP9.5/UCLH1 Antibody, Cite this (**Proteintech Group Cat# 14730-1-AP,** [RRID:AB_2210497](#))

SAFETY WARNINGS

WARNING HAZARDOUS CONDITION WARNED AGAINST. This comment describes a hazardous condition to which the technician may be exposed in the performance of this protocol. It also contains directions on how to avoid or minimize the danger. Warnings are always and only used for personnel safety, and precedes the first step that will expose the technician to the hazard.

1 Reagents: Phosphate buffer (PB, 0.1 M, pH 7.2) Sucrose

Equipment: Razor blades, embedding molds, 24 well plate

Solutions: 2% zamboni's in phosphate buffer (PB, 0.1 M, pH 7.4)
30% sucrose in PBS


Procedure:

Fresh mouse or rat:



1. Remove footpad and fix (4 hours) in 2% zamboni's fixative for 4-6 hours.
2. Place skin in 30% sucrose in 1x PBS overnight.
3. Embed tissue in OCT using scheme below.
4. Cut 10-14 sections at 30um using scheme below and place into a single well of a 24 well plate.
5. Follow IHC protocol using Proteintech rabbit anti-PGP9.5(14730-1-ap), 1:2000 followed by AlexaFluor 488 highly cross absorbed, 1:1000 (Molecular Probes). Apply coverslips with ProLong gold or platinum antifade kit with dapi (Molecular Probes)
6. Block for 1-2 hours in blocking solution.
7. Primary at RT for 1 hour, then in to cold room overnight.
8. Rinse 3x1 hour in PBS
9. Secondary at RT for 1 hour, then in to cold room overnight.
10. Rinse 3x1 hour in PBS.
11. Coverslip using dapi gold with prolong.

Imaging:

1. Three images per sample are collected on an Olympus FluoView 500 confocal microscope using a 20 X air objective at a resolution of 1024 X 1024 pixels. The optical section thickness is 3.3um. 10 images per stack are flattened using MetaMorph(version 7.7.0) arithmetic max option. Data are presented as the fibers crossing the epidermal border of PGP9.5 positive fibers per area of epidermis.

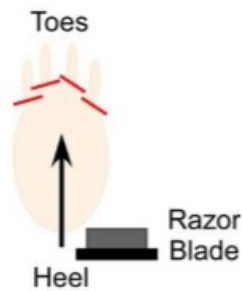

To dissect the foot pad, cut the skin at the toe joints with a sharp blade. Begin at the heel and move forward to remove the entire plantar surface of the foot. Handle the tissue by the heel only so as not to disrupt structures within the dermis and epidermis.

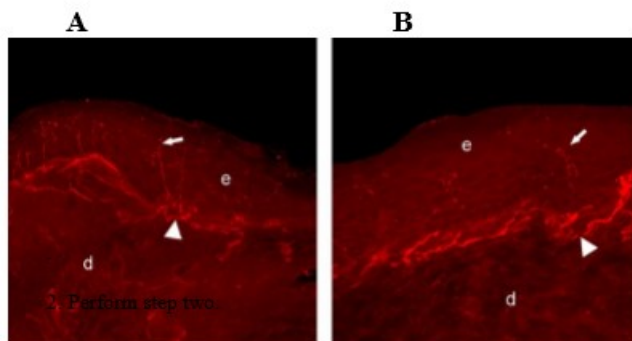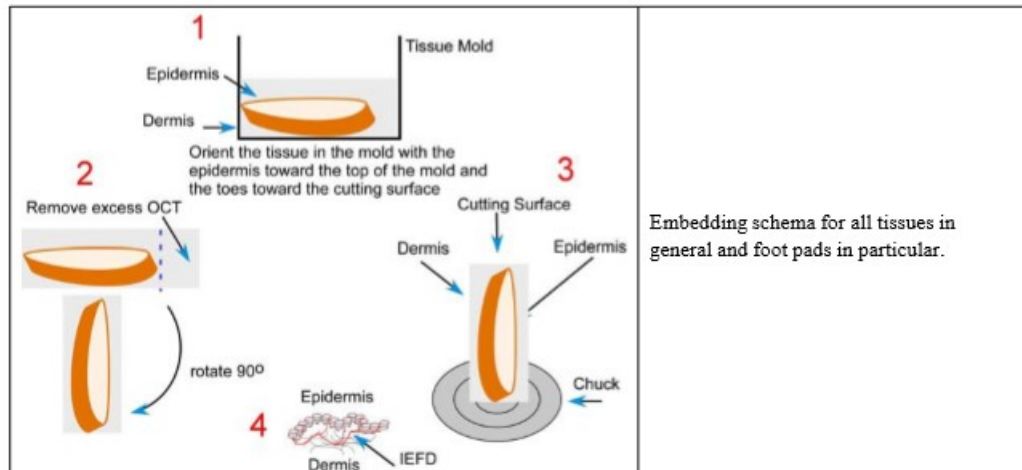

IEFD in a control (A) and diabetic (B) mouse

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited