



Jun 12,  
2019

## Agro Preparation for *Mimulus in Planta* Transformation [↗](#)

Yaowu Yuan<sup>1</sup>

<sup>1</sup>University of Connecticut

Working

[dx.doi.org/10.17504/protocols.io.3rqgm5w](https://doi.org/10.17504/protocols.io.3rqgm5w)

Mimulus



Andrea Sweigart ⚡

### ABSTRACT

This protocol is part of a [collection for \*Mimulus in planta\* transformation](#).

### EXTERNAL LINK

[http://mimubase.org/FTP/Protocols/Stable\\_Transformation/Mimulus%20in%20planta%20transformation.pdf](http://mimubase.org/FTP/Protocols/Stable_Transformation/Mimulus%20in%20planta%20transformation.pdf)



Mimulus in planta  
transformation.pdf





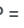







### SAFETY WARNINGS

See SDS (Safety Data Sheet) for warnings and hazards.

- 1 Inoculate 4 agro colonies into separate 5 ml tubes of LB+Kan+Gent+Rif.
- 1.1 Shake overnight at 28 °C- 30 °C.
- 2 Colony PCR from the 5 mL culture to check that it has your insert.
- 3 Inoculate correct agro colony into a flask of 300 ml LB+Kan+Gent+Rif, and cover with foil.
- 3.1 Shake for 12:00:00 to 16:00:00 at 28 °C- 30 °C.



If left longer, the bacteria will begin to die.

- 4 Make glycerol stocks by adding  **1 ml** of agro to glycerol stock tube.
- 4.1 Place in liquid nitrogen then store in  **-80 °C** freezer.
- 5 Transfer contents of flask into 500 mL centrifuge bottles and label them.
- 6 Balance the bottles (use the scale and spray bottle of water to balance).
- 7 Centrifuge to pellet the Agro.  
SPEED =  **6000 x g**  
TIME =  **00:15:00**  
TEMP =  **4 °C**
- 8 While centrifuging, make the resuspension solution:  
 **300 ml** total for each construct  
5% sucrose =  **15 g** sucrose  
 **0.1 Molarity (M)** acetosyringone = 2 mL/L =  **600 µl**
- 8.1 Make acetosyringone fresh by dissolving  **0.0196 g**/mL in methanol.
- 8.2 Dissolve sucrose.
- 8.3 Add acetosyringone.
- 8.4  
Add dH<sub>2</sub>O to bring total volume to  **300 ml**, and mix.
- 9 After centrifuging, pour off liquid into the original flasks the Agro was grown in, and bleach the flasks.
- 10 Resuspend the Agro pellet in  **300 ml resuspension solution** by adding the solution to the centrifuge bottle and shaking.
- 11 Pour into spray bottle for infiltration.

12 Add  300 µl **Silwett** to the Agro and shake to mix.



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited