



Sep 11, 2019

MojoSort™ Mouse CD4 Nanobeads Column Protocol [↗](#)Sam Li¹¹BioLegend

1

Works for me

[dx.doi.org/10.17504/protocols.io.696hh9e](https://doi.org/10.17504/protocols.io.696hh9e)

BioLegend

Sam Li
BioLegend

ABSTRACT

BioLegend MojoSort™ nanobeads work in commonly used separation columns, based on our internal research as well as validation by external testing by academic labs. This simple protocol consists of following the MojoSort™ protocol to label the cells with **pre-diluted** MojoSort™ reagents and using the columns as indicated by the manufacturer.

Note: Due to the properties of our beads, it may be possible to use far fewer beads than with other commercial suppliers. We recommend a titration to find the best dilution factor. However, as a general rule, dilutions ranging from 1:3 to 1:20 for the Nanobeads can be used. Please contact BioLegend Technical Service (tech@biolegend.com) if further assistance is needed.

EXTERNAL LINK

<https://www.biolegend.com/protocols/mojosort-mouse-cd4-nanobeads-column-protocol/4766/>

GUIDELINES

MojoSort™ magnetic particles can be used with other commercially available magnetic separators, both free standing magnets and column-based systems. Because MojoSort™ protocols are optimized for the MojoSort™ separator, the protocols may need to be adjusted for other systems. Please contact BioLegend Technical Service (tech@biolegend.com) for more information and guidance. We do not recommend using MojoSort™ particles for BD's iMag™ or Life Technologies' DynaMag™.

MATERIALS

NAME ▾	CATALOG # ▾	VENDOR ▾
MojoSort™ Buffer	480017	BioLegend
TruStain FcX™ PLUS (anti-mouse CD16/32) Antibody	156603	BioLegend
MojoSort™ Mouse CD4 Nanobeads	480069	

MATERIALS TEXT

Additional reagents:

- commercially available cell separation columns
- 5 mL polypropylene tubes
- 70 µm cell strainer


- 1 Prepare cells from your tissue of interest or blood without lysing erythrocytes.
- 2 In the final wash of your sample preparation, resuspend the cells in MojoSort™ Buffer by adding up to 4 mL in a 5 mL (12 x 75 mm) polypropylene tube.
Note: Keep MojoSort™ Buffer on ice throughout the procedure.

- 3 Filter the cells with a 70 μm cell strainer, centrifuge at 300xg for 5 minutes, and resuspend in an appropriate volume of MojoSort™ Buffer. Count and adjust the cell concentration to 1×10^8 cells/mL by adding MojoSort™ Buffer. 5m
- 4 Aliquot 100 μL (10^7 cells) into a new tube. **Add 5 μL of Mouse TruStain FcX™** mix well and **incubate at room temperature for 10 minutes**. Scale up the volume accordingly if separating more cells. For example, if the volume of Mouse TruStain FcX™ for 1×10^7 cells is 5 μL , add 50 μL for 1×10^8 cells. When working with less than 10^7 cells, use indicated volumes for 10^7 cells. 10m
- 5 Vortex the antibody-conjugated Nanobeads (to resuspend) at max speed, 5 touches, and prepare the dilutions to test. **Add 10 μL of pre-diluted conjugated Nanobeads**. Mix well and **incubate on ice for 15 minutes**. Scale up the volume accordingly if separating more cells. For example, add 100 μL of pre-diluted Nanobeads for separating 1×10^8 cells in 1 ml of MojoSort™ Buffer. When working with less than 10^7 cells, use indicated volumes for 10^7 cells. 15m
- 6 Add the appropriate amount of MojoSort™ Buffer and proceed to separation. At least 500 μL is needed for column separation. **Note:** There are several types of commercially available columns, depending on your application. Choose the one that fits best your experiment:

	Max. number of labeled cells	Max. number of total cells	Cell suspension volume	Column rinse volume	Cell wash volume	Elution volume
Small Capacity	1×10^7	2×10^8	500 μL for up to 10^8 cells	1ml	1 ml	1 ml
Medium Capacity	1×10^8	2×10^9	500 μL for up to 10^9 cells	3ml	3 ml	5 ml
Large Capacity	1×10^9	2×10^{10}	500 μL for up to 10^{10} cells	20-50ml	30 ml	20 ml

Example of magnetic separation with medium capacity columns:

- 7 Place the column in a magnetic separator that fits the column.
- 8 Rinse the column with 3 mL of cell separation buffer.
- 9 Add the labeled cell suspension to the column through a 30 μm filter and collect the fraction containing the unlabeled cells.
- 10 Wash the cells in the column **3 times** with 3 mL of buffer and collect the fraction containing the unlabeled cells. Combine with the collected fraction from step 3. These cells may be useful as controls, to monitor purity/yield, or other purposes.
- 11 Take away the column from the magnet and place it on a tube. Then add 5 mL of buffer and flush out the magnetically labeled fraction with a plunger or supplied device. These are the positively isolated cells of interest; do not discard. To increase the purity of the magnetically labeled fraction repeat the isolation process with a new, freshly prepared column.

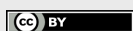

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited