



Working

# Protocol for use with NEBNext rRNA Depletion Kit (Human/Mouse/Rat) E6310

Lauren Moran<sup>1</sup>, New England Biolabs<sup>1</sup>

<sup>1</sup>New England Biolabs

dx.doi.org/10.17504/protocols.io.sxuefnw

New England Biolabs (NEB)









#### ABSTRACT

The NEBNext rRNA Depletion Kit (Human/Mouse/Rat) depletes both cytoplasmic (5S rRNA, 5.8S rRNA, 18S rRNA and 28S rRNA) and mitochodrial ribosomal RNA (12S rRNA and 16S rRNA) from human, mouse and rat total RNA preparations. This product is suitable for both intact and degraded RNA (e.g. FFPE RNA). The resulting rRNA-depleted RNA is suitable for RNA-Seq, random $primed\ cDNA\ synthesis, or\ other\ downstream\ RNA\ analysis\ applications.$ 

**PROTOCOL STATUS** 

#### Working

	NAME V NEBNext RNase H	CATALOG #	VENDOR V
	NERNAYT PNaca H		
_	INLDINEXT INVASE II	E6318	New England Biolabs
	RNase H Reaction Buffer	E6312	New England Biolabs
	NEBNext rRNA Depletion Solution	E6313	New England Biolabs
	NEBNext Probe Hybridization Buffer	E6314	New England Biolabs
	DNase I (RNase-free)	E6316	New England Biolabs
	DNase I Reaction Buffer	E6315	New England Biolabs
	Nuclease-free Water	E6317	New England Biolabs
	NEBNext RNA Sample Purification Beads	E6315	New England Biolabs
	Magnetic Rack	View	Contributed by users
	80% Ethanol (freshly prepared)	View	Contributed by users
	Thermal cycler	View	Contributed by users
	Agencourt RNAClean XP Beads	A63987	Beckman Coulter
	DNase I (e.g., NEB #M0303) and DNase I Cleanup Reagants or Kit for Removal of DNA Prior to Depletion	View	New England Biolabs
	Random Primers	E7422	New England Biolabs
TEPS M	MATERIALS		
	NAME Y	CATALOG #	VENDOR ~
	NEBNext RNase H	E6318	New England Biolabs
	RNase H Reaction Buffer	E6312	New England Biolabs
	Nuclease-free Water	E6317	New England Biolabs
	DNase I Reaction Buffer	E6315	New England Biolabs

NAME ~	CATALOG #	VENDOR >
DNase I (RNase-free)	E6316	New England Biolabs
Nuclease-free Water	E6317	New England Biolabs
NEBNext Sample Purification Beads	E6315	New England Biolabs
Nuclease-free Water	E6317	New England Biolabs
First Strand Synthesis Reaction Buffer	E7421	New England Biolabs
Random Primers	E7422	New England Biolabs
NEBNext rRNA Depletion Solution	E6313	New England Biolabs
NEBNext Probe Hybridization Buffer	E6314	New England Biolabs
80% Ethanol (freshly prepared)		Contributed by users
NEBNext Strand Specificity Reagent	E7766	New England Biolabs
NEBNext First Strand Synthesis Enzyme Mix	E7761	New England Biolabs
NEBNext Second Strand Synthesis Reaction Buffer with dUTP Mix	E7426	New England Biolabs
NEBNext Second Strand Synthesis Enzyme Mix	E7425	New England Biolabs
Nuclease-free Water	E7764	New England Biolabs
NEBNext Sample Purification Beads	E6315	New England Biolabs
Fresh 80% Ethanol		Contributed by users
(0.1X) TEBuffer	E7763	New England Biolabs
NEBNext Ultra II End Prep Reaction Buffer	E7647	New England Biolabs
NEBNext Ultra II End Prep Enzyme Mix	E7646	New England Biolabs
NEBNext Ligation Enhancer	E7374	New England Biolabs
NEBNext Ultra II Ligation Master Mix	E7648	New England Biolabs
NEBNext USER Enzyme	E7458	New England Biolabs
NEBNext Sample Purification Beads	E7767	New England Biolabs
80% Ethanol (freshly prepared)		Contributed by users
NEBNext Sample Purification Beads	E6315	New England Biolabs
80% Ethanol (freshly prepared)		Contributed by users

BEFORE STARTING

# **RNA Sample Recommendations**

**RNA Integrity:** Assess the quality of the input RNA by running the RNA sample on an Agilent Bioanalyzer RNA 6000 Nano/Pico Chip to determine the RNA Integrity Number (RIN). RNA with different RIN values require different fragmentation times or no fragmentation at all.

For intact (RIN > 7) or partially degraded RNA samples (RIN = 2 to 7) follow the library preparation protocol in Chapter 2 (current chapter). See Table 2.5.1. for the recommended fragmentation times, based on RIN.

For highly degraded samples (RIN = 1 to 2) (e.g. FFPE), which do not require fragmentation, follow the library preparation protocol in Chapter 3.

**RNA Purity:** Samples should be free of DNA. The RNA sample should be free of salts (e.g., Mg2+, or guanidinium salts), divalent cation chelating agents (e.g., EDTA, EGTA, citrate), or organics (e.g., phenol and ethanol).

\_\_\_\_\_\_

--

### **Input Amount Requirement**

 $5 \text{ ng} - 1 \mu \text{g}$  total RNA (DNA-free) in up to 12  $\mu$ l of Nuclease-free Water, quantified by Qubit Fluorometer and quality checked by Bioanalyzer.

The protocol is optimized for approximately 200 nt RNA inserts. To generate libraries with longer RNA insert sizes, refer to Appendix A (Chapter 6) for recommended fragmentation times and size selection conditions.

Keep all of the buffers on ice, unless otherwise indicated.

# Probe Hybridization to RNA

- 1 Dilute the total RNA with Nuclease-free Water to a final volume of 12 μl in a PCR tube. Keep the RNA on ice.
- 2 Prepare a RNA/Probe master mix as follows:

RNA Probe Master Mix	Volume
NEBNext rRNA Depletion Solution	1 μΙ
Probe Hybridization Buffer	2 μΙ
Total Volume	3 μΙ





- 3 Add 3  $\mu$ l of the above mix to 12  $\mu$ l total RNA (from Step 1), resulting in a total volume of 15  $\mu$ l.
- 4 Mix by pipetting up and down at least 10 times.
- 5 Spin down briefly in a microcentrifuge.
- 6 Place samples in a thermocycler with a heated lid set to approximately 105°C, and run the following program, which will take

approximately 15-20 minutes to complete:

Temp	Time
95°C	2 min
95-22°C	0.1°C/sec
22°C	5 min hold

7 Spin down the samples in a microcentrifuge, and place on ice. Proceed immediately to RNase H Digestion Step.

# RNase H Digestion

Assemble the RNAse H master mix on ice as follows.

RNase H Master Mix	Volume
NEBNext RNase H	2 μΙ
RNase H Reaction Buffer	2 μΙ
Nuclease-free Water	1 μΙ
Total Volume	5 μΙ







- 9 Mix thoroughly by pipetting up and down at least 10 times.
- $10 \quad \text{ Briefly spin down the samples in a microcentrifuge}.$
- Add 5  $\mu$ l of the RNase H master mix to the RNA sample from Step 7, resulting in a total volume of 20  $\mu$ l.
- 12 Mix thoroughly by pipetting up and down at least 10 times.

Incubate the sample in a thermocycler for 30 minutes at 37°C with the lid set to 40°C (or on).

- 13 **© 00:30:00** Incubate
- Briefly spin down the samples in a microcentrifuge, and place on ice. Proceed immediately to DNase I Digestion to prevent non-specific degradation of RNA.

# DNase I Digestion

 $15 \qquad \text{Assemble the DNase I master mix} \, \textbf{on ice} \, \, \text{in a nuclease-free tube}.$ 

DNase I Master Mix	Volume
DNase I Reaction Buffer	5 μΙ
DNase I (RNase-free)	2.5 μΙ
Nuclease-free Water	22.5μΙ
Total Volume	30 μΙ







- 16 Mix thoroughly by pipetting up and down at least 10 times.
- 17 Briefly spin down the sample in a microcentrifuge.
- 18 Add 30 μl of DNase I master mix to 20 μl RNA sample from Step 14, resulting in a total volume of 50 μl.
- 19 Mix thoroughly by pipetting up and down 10 times.
- 20 Incubate the sample in a thermocycler for 30 minutes at 37°C with the heated lid set to 40°C (or on).

**©00:30:00** Incubate

21 Briefly spin down the sample in a microcentrifuge, and place on ice. Proceed immediately to RNA Purification.

RNA Purification after rRNA Depletion Using Agencourt RNAClean XP Beads or NEBNext RNA Sample Purification Beads

- 27 Vortex Agencourt RNAClean XP Beads or RNA Sample Purification Beads to resuspend.
- Add 110 µl (2.2X) beads to the RNA sample from Step 21 and mix thoroughly by pipetting up and down at least 10 times.



24 Incubate the sample for 15 minutes on ice to bind RNA to the beads.

**© 00:15:00 Incubate** 

- Place the tube on a magnetic rack to separate beads from the supernatant. After the solution is clear, carefully remove and discard the supernatant. Be careful not to disturb the beads, which contain RNA.
- Add 200 µl of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for 30 seconds, and then carefully remove and discard the supernatant. Be careful not to disturb the beads that contain the RNA.



27 Repeat Step 26 once for a total of 2 washing steps.

28 Completely remove residual ethanol, and air dry the beads for up to 5 minutes while the tube is on the magnetic rack with the lid open.

Caution: Do not over-dry the beads. This may result in lower recovery of RNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.

Remove the tube from the magnet. Elute the RNA from the beads by adding 7 µl Nuclease-free Water. Mix well by pipetting up and down at least 10 times and briefly spin the tube.



30 Incubate for 2 minutes at room temperature. Place the tube in the magnet until the solution is clear (~2 minutes).

**©00:02:00** Incubate

- 31 Remove 5  $\mu$ l of the supernatant containing RNA and transfer to a nuclease-free tube.
- 32 Place the sample on ice and proceed to RNA Fragmentation and Priming.

# RNA Fragment at ion and Priming

ARNA fragmentation is only required for intact or partially degraded RNA. Recommended fragmentation times can be found in Table 2.5.1.

2.5.1. Assemble the following fragmentation and priming reaction on ice:

Fragmentation and Priming Reaction	Volume
Ribosomal RNA Depleted Sample (Step 32)	5 µl
NEBNext First Strand Synthesis Reaction Buffer	4 µl
Random Primers	1 µl
Total Volume	10 μΙ





- 34 Mix thoroughly by pipetting up and down 10 times.
- Place the sample on a thermocycler and incubate the sample at 94°C following the recommendations in Table 35 below for libraries with inserts ~200 nt.

Table 35: Suggested fragmentation times based on RIN value of RNA input.

	RNA Type	RIN	Frag. Time
--	----------	-----	------------

Intact RNA	> 7	15 min. @ 94°C
Partially Degraded DNA	2-6	7-8 min. @ 94°C

#### **NOTE**

Refer to Appendix A (Chapter 6) for fragmentation conditions if you are preparing libraries with large inserts (> 200 bp). Conditions in Appendix A only apply for intact RNA.

36 Immediately transfer the tube to ice and proceed to First Strand cDNA Synthesis.

# First Strand cDNA Synthesis

37 Assemble the first strand synthesis reaction **on ice** by adding the following components to the fragmented and primed RNA from Step 36:

First Strand Synthesis Reaction	Volume
Fragmented and Primed RNA (Step 36)	10 µl
NEBNext Strand Specificity Reagent	8 µl
NEBNext First Strand Synthesis Enzyme Mix	2 μΙ
Total Volume	20 µl





NEBNext Strand

Specificity Reagent by New England Biolabs

Catalog #: E7766

# REAGENT



NEBNext First Strand

Synthesis Enzyme Mix

by New England Biolabs

Catalog #: E7761

- 38 Mix thoroughly by pipetting up and down 10 times.

Note: If you are following recommendations in Appendix A (Chapter 6), for libraries with longer inserts (> 200 bases), increase the incubation at 42°C from 15 minutes to 50 minutes at Step 2 below.

Step 1: 10 minutes at 25°C

Step 2: 15 minutes at 42°C

Step 3: 15 minutes at 70°C

Step 4: Hold at 4°C

40 Proceed directly to Second Strand cDNA Synthesis.

# Second Strand cDNA Synthesis

41 Assemble the second strand cDNA synthesis reaction on ice by adding the following components into the first strand synthesis

product from Step 40.

Second Strand Synthesis Reaction	Volume
First Strand Synthesis Product (Step 40)	20 μΙ
NEBNext Second Strand Synthesis Reaction Buffer with dUTP (10X)	8 μΙ
NEBNext Second Strand Synthesis Enzyme Mix	4 µl
Nuclease-free Water	48 µl
Total Volume	80 µl

#### **REAGENT**



NEBNext Second Strand

Synthesis Reaction Buffer with dUTP Mix

by New England Biolabs

Catalog #: E7426

#### **REAGENT**



**NEBNext Second Strand** 


Synthesis Enzyme Mix

by New England Biolabs

Catalog #: E7425



#### **REAGENT**

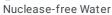

by New England Biolabs

Catalog #: E7764

- 42 Keeping the tube on ice, mix thoroughly by pipetting up and down at least 10 times.
- 43 Incubate in a thermocycler for 1 hour at 16°C with the heated lid set at  $\leq$  40°C (or off).

**© 01:00:00** Incubate

Purification of Double-stranded cDNA Using SPRIselect Beads or NEBNext Sample Purification Beads

- Vortex SPRIselect Beads or NEBNext Sample Purification Beads to resuspend.
- Add 144  $\mu$ l (1.8 $\chi$ ) of resuspended beads to the second strand synthesis reaction (~80  $\mu$ l). Mix well on a vortex mixer or by pipetting up and down at least 10 times.





**NEBNext Sample** 

Purification Beads


by New England Biolabs

Catalog #: E6315

- 46 Incubate for 5 minutes at room temperature. © 00:05:00 Incubate
- Briefly spin the tube in a microcentrifuge to collect any sample on the sides of the tube. Place the tube on a magnet to separate beads from the supernatant. After the solution is clear, carefully remove and discard the supernatant. Be careful not to disturb the beads, which contain DNA.

Caution: do not discard beads.

Add 200 µl of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for 30 seconds, and then carefully remove and discard the supernatant.

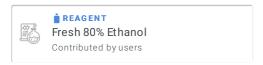

49 Repeat Step 48 once for a total of 2 washing steps.

☼ go to step #48 Repeat Step

 $\,\,50\,\,$  Air dry the beads for up to 5 minutes while the tube is on the magnetic rack with lid open.

Caution: Do not over-dry the beads. This may result in lower recovery of DNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.

Factor of the season of the s



52  $\;$  Remove 50  $\mu l$  of the supernatant and transfer to a clean nuclease-free PCR tube.

**NOTE** 

If you need to stop at this point in the protocol samples can be stored at  $-20\,^{\circ}\text{C}$ .

#### End Prep of cDNA Library

Assemble the end prep reaction on ice by adding the following components to the second strand synthesis product from Step 52.

End Prep Reaction	Volume
Second Strand Synthesis Product (Step 52)	50 μl
NEBNext Ultra II End Prep Reaction Buffer	7 µl
NEBNext Ultra II End Prep Enzyme Mix	3 μΙ
Total Volume	60 µl

If a master mix is made, add 10  $\mu$ l of master mix to 50  $\mu$ l of cDNA for the End Prep reaction.

# **REAGENT**



NEBNext Ultra II End Prep

Reaction Buffer by New England Biolabs

Catalog #: E7647

#### **REAGENT**





Enzyme Mix

by New England Biolabs

Catalog #: E7646

Set a 100 µl or 200 µl pipette to 50 µl and then pipette the entire volume up and down at least 10 times to mix thoroughly. Perform a 54 quick spin to collect all liquid from the sides of the tube.

It is important to mix well. The presence of a small amount of bubbles will not interfere with performance.

Incubate the sample in a thermocycler with the heated lid set at  $\geq 75$ °C as follows. 55

30 minutes at 20°C

30 minutes at 65°C

Hold at 4°C.

Proceed immediately to Adaptor Ligation. 56

# Adapt or Ligation

57 Dilute the NEBNext Adaptor\* prior to setting up the ligation reaction in ice-cold Adaptor Dilution Buffer and keep the adaptor on ice.

Total RNA Input	Dilution Required
1,000 ng-101 ng	5-fold dilution in Adaptor Dilution Buffer
100 ng-10 ng	25-fold dilution in Adaptor Dilution Buffer
5 ng	200-fold dilution in Adaptor Dilution Buffer

<sup>\*</sup>The NEBNext adaptor is provided in NEBNext oligos kit. NEB has several oligo kit options, which are supplied separately from the library prep kit.

58 Assemble the ligation reaction on ice by adding the following components, in the order given, to the end prep reaction product from Step 56.

Ligation Reaction	Volume
End Prepped DNA (Step 56)	60 µl
Diluted Adaptor (Step 57)	2.5 µl
NEBNext Ligation Enhancer	1 µl
NEBNext Ultra II Ligation Master Mix	30 µl
Total Volume	93.5 µl

#### **NOTE**

The Ligation Master Mix and Ligation Enhancer can be mixed ahead of time and is stable for at least 8 hours @ 4°C. We do not recommend premixing the Ligation Master Mix, Ligation Enhancer and adaptor prior to use in the Adaptor Ligation Step.



NEBNext Ligation Enhancer

by New England Biolabs

Catalog #: E7374

#### **REAGENT**

NEBNext Ultra II Ligation



Master Mix

by New England Biolabs
Catalog #: E7648

59 Set a 100 μl or 200 μl pipette to 80 μl and then pipette the entire volume up and down at least 10 times to mix thoroughly. Perform a quick spin to collect all liquid from the sides of the tube.

Caution: The NEBNext Ultra II Ligation Master Mix is very viscous. Care should be taken to ensure adequate mixing of the ligation reaction, as incomplete mixing will result in reduced ligation efficiency. The presence of a small amount of bubbles will not interfere with performance.

60 Incubate 15 minutes at 20°C in a thermocycler.

**© 00:15:00** Incubate

61 Add 3 μl (blue) USER Enzyme to the ligation mixture from Step 60, resulting in total volume of 96.5 μl.



62 Mix well and incubate at 37°C for 15 minutes with the heated lid set to  $\geq$  45°C.

**© 00:15:00 Incubate** 

63 Proceed immediately to Purification of the Ligation Reaction.

Purification of the Ligation Reaction Using SPRIselect Beads or NEBNext Sample Purification Beads

64 Mote: If you are selecting for libraries with larger insert size (> 200 nt) follow the size selection recommendations in Appendix A, Chapter 6.

Add 87  $\mu$ l (0.9X) resuspended SPRIselect Beads or NEBNext Sample Purification Beads and mix well on a vortex mixer or by pipetting up and down at least 10 times.



65 Incubate for 10 minutes at room temperature.

**© 00:10:00** Incubate

Quickly spin the tube in a microcentrifuge and place the tube on an appropriate magnetic rack to separate beads from the supernatant. After the solution is clear (~ 5 minutes), discard the supernatant that contains unwanted fragments.

Caution: do not discard beads.

Add 200 µl of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for 30 seconds, and then carefully remove and discard the supernatant.



Repeat Step 67 once for a total of 2 washing steps.

☆ go to step #67 Repeat Step

- Briefly spin the tube, and put the tube back in the magnetic rack.
- 70 Completely remove the residual ethanol, and air dry beads until the beads are dry for up to 5 minutes while the tube is on the magnetic rack with the lid open.

Caution: Do not over-dry the beads. This may result in lower recovery of DNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.

Remove the tube from the magnetic rack. Elute DNA target from the beads by adding 17 μl 0.1X TE (provided) to the beads. Mix well on a vortex mixer or by pipetting up and down. Quickly spin the tube and incubate for 2 minutes at room temperature. Put the tube in the magnet until the solution is clear.

**© 00:02:00** Incubate

72 Without disturbing the bead pellet, transfer 15 µl of the supernatant to a clean PCR tube and proceed to PCR enrichment.

**NOTE** 

If you need to stop at this point in the protocol samples can be stored at  $-20\,^{\circ}$ C.

# PCR Enrichment of Adaptor Ligated DNA

73 — Check and verify that the concentration of your oligos is 10 μM on the label.

Luse Option A for any NEBNext oligos kit where index primers are supplied in tubes. These kits have the forward and reverse primers supplied in separate tubes.

Use Option B for any NEBNext oligos kit where index primers are supplied in a 96-well plate format. These kits have the forward and reverse (i7 and i5) primers combined.

Set up the PCR reaction as described below based on the type of oligos (PCR primers) used.

#### 74 Option A:

Component	Volume Per One Library	
Adaptor Ligated DNA (Step 72)	15 μΙ	
NEBNext Ultra II Q5 Master Mix	25 μΙ	
Universal PCR Primer/i5 Primer*,**	5 µl	
Index (X) Primer/i7 Primer*,**	5 μl	
Total Volume	50 μΙ	

## Option B:

Component	Volume Per One	
	Library	
Adaptor ligated DNA (Step 72)	15 µl	
NEBNext Ultra II Q5 Master Mix	25 μΙ	
Index (X)/Universal Primer Mix*	10 μΙ	
Total Volume	50 μl	

<sup>\*</sup> NEBNext Oligos must be purchased separately from the library prep kit. Refer to the corresponding NEBNext Oligo kit manual for determining valid barcode combinations.

75 Mix well by gently pipetting up and down 10 times. Quickly spin the tube in a microcentrifuge.

76 Place the tube on a thermocycler with the heated lid set to 105°C and perform PCR amplification using the following PCR cycling conditions (refer to Table 76.A and Table 76.B):

Table 76.A:

Cycle Step	Temp	Time	Cycles
Initial Denaturation	98°C	30 seconds	1
Denaturation	98°C	10 seconds	7-16*,**
Annealing/Extension	65°C	75 seconds	
Final Extension	65°C	5 minutes	1
Hold	4°C	∞	

<sup>\*</sup>The number of PCR cycles should be adjusted based on RNA input.

Table 76.B: Recommended PCR cycles based on total RNA input amount:

Total RNA Input	Recommended PCR Cycles
1,000 ng	7–8
100 ng	11-12
10 ng	14-15
5 ng	15-16

https://www.neb.com/-/media/catalog/datacards-or-manuals/manuale7760.pdf

# Purification of the PCR Reaction using SPRIselect Beads or NEBNext Sample Purification Beads

77 Vortex SPRIselect Beads or NEBNext Sample Purification Beads to resuspend.

Add  $45 \,\mu\text{l}$  (0.9X) of resuspended beads to the PCR reaction (~  $50 \,\mu$ l). Mix well on a vortex mixer or by pipetting up and down at least 10 times.

<sup>\*\*</sup> Use only one i7 primer/index primer per sample. Use only one i5 primer (or the universal primer for single index kits) per sample


<sup>\*\*</sup> It is important to limit the number of PCR cycles to avoid overamplification. If overamplification occurs, a second peak ~ 1,000 bp will appear on the Bioanalyzer trace (see Figure 7.2 on page 79 in manual).



- 79 Incubate for 5 minutes at room temperature. © 00:05:00 Incubate
- Quickly spin the tube in a microcentrifuge and place the tube on an appropriate magnetic rack to separate beads from the supernatant. After the solution is clear (about 5 minutes), carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

Caution: Do not discard beads.

Add 200 µl of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for 30 seconds, and then carefully remove and discard the supernatant.

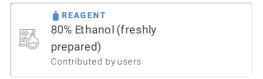

- 82 Repeat Step 81 once for a total of 2 washing steps. 5 go to step #81 Repeat Step
- Air dry the beads for up to 5 minutes while the tube is on the magnetic rack with the lid open.

  Caution: Do not over-dry the beads. This may result in lower recovery of DNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.
- Remove the tube from the magnetic rack. Elute the DNA target from the beads by adding 23 µl 0.1X TE (provided) to the beads. Mix well on a vortex mixer or by pipetting up and down ten times. Quickly spin the tube in a microcentrifuge and incubate for 2 minutes at room temperature. Place the tube in the magnetic rack until the solution is clear.
- Transfer 20  $\mu$ l of the supernatant to a clean PCR tube, and store at -20 °C.

# Assess Library Quality on an Agilent Bioanalyzer DNA Chip


- Run 1 µl library on a DNA 1000 chip. If the library yield is too low to quantify on this chip, please run the samples on a DNA High Sensitivity chip. A dilution may be necessary for running on a Bioanalyzer High Sensitivity DNA Chip.
- 87 Check that the electropherogram shows a narrow distribution with a peak size approximately 300 bp.

**■**NOTE

Note: If a peak at ~ 80 bp (primers) or 128 bp (adaptor-dimer) is visible in the bioanalyzer traces,

bring up the sample volume (from Step 85) to 50  $\mu$ l with 0.1X TE buffer and repeat the SPRIselect Bead or NEBNext Sample Purification Bead Cleanup Step (Section "Purification of the PCR Reaction using SPRIselect Beads or NEBNext Sample Purification Beads").

88 Figure 88: Example of RNA library size distribution on a Bioanalyzer.

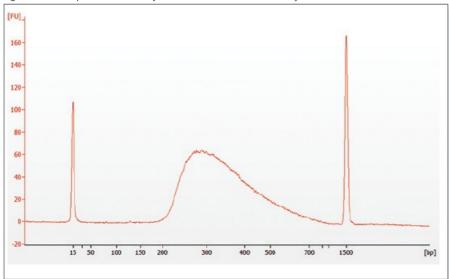

This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited