

## Isolated astrocyte culture - protocol 2

Yasmin Bar El<sup>1</sup>

<sup>1</sup>School of Physics and Astronomy, Tel-Aviv University, Tel-Aviv, Israel

Sep 21, 2018 dx.doi.org/10.17504/protocols.io.tqcemsw

Working



Yasmin Bar El ⚡

### ABSTRACT

An alternative protocol for isolated astrocyte cell culture preparation. Dissociation, preparation and plating of mice cortex neurons and glia cells on flasks. Isolation process of the culture from neurons, to obtain an isolated astrocyte culture. Plating isolated astrocytes on MEA.

### PROTOCOL STATUS


#### Working

We use this protocol in our group and it is working

#### Sample collection

- 1 Dissect and place on ice cortices from pups (post natal day 0 or 1 mice).

#### Cell lysis

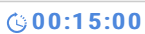
- 2
  1. Chopp with scissors in a papain-based dissociation buffer (2.5 mM CaCl<sub>2</sub>, 0.83 mM EDTA, 137 U papain (Sigma-Aldrich)), 100 µl DNase (Sigma-Aldrich), 3–5 crystals of L-Cysteine (Sigma-Aldrich), HBSS with 20 mM HEPES (pH 7.4);
  2. Place on a rotating shaker for 15 min  00:15:00 at room temperature.

#### Cell preparation

- 3
  1. After centrifuging, discard the supernatant, resuspend the pellet in modified essential medium (MEM) without L-glutamine with essential amino acids (Beit Haemek, 06-1025-01-1A), 5% heat-inactivated fetal calf serum (Biological Industries), heat-inactivated 5% horse serum (Beit Haemek, 04-004-1), 2 mM glutamine (Beit Haemek, 03-020-1c), 3 mg/ml glucose, 2% B-27 (Gibco, 17504-044), 0.5% Pen/Strep (100 U/ml penicillin, 100 µg/ml streptomycin; Beit Haemek, 03-031-1B);
  2. Triturate seven times.

#### Cell culture

- 4
  1. Plate the cells on poly-D-lysine (PDL, Sigma-Aldrich, P7405-5MG) coated T75 flasks;
  2. Maintain cultures at 37 °C with 5% CO<sub>2</sub>. Replace growth medium every 3–4 days to astrocyte enriched growth medium: MEM-EAGLE (without L-glutamine, with essential amino acids) (Beit Haemek), 3 mg/ml glucose, 10% heat-inactivated fetal calf serum (Biological Industries), 0.8% GlutaMAX (100X; Gibco), 0.5% Pen/Strep, 2 mM glutamine (Beit Haemek);
  3. After 20 DIV detach using trypsin/EDTA;
  4. Re-plate the cells on poly-D-lysine (PDL, Sigma-Aldrich, P7405-5MG) coated micro-electrode arrays (MEAs; 200/30iR-Ti-gr and 500/30iR-Ti-pr; Multichannel Systems) with a cell density of 2000–2500 cells/mm<sup>2</sup> (~106 cells per dish).

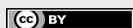

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited