

HMW DNA extraction for insects

Tom Harrop

Abstract

Adapted from Qiagen's [genomic DNA handbook](#) and the [user protocol for mosquito DNA extraction](#).

Citation: Tom Harrop HMW DNA extraction for insects. **protocols.io**

dx.doi.org/10.17504/protocols.io.pnwdmfe

Published: 30 Apr 2018

Protocol

Tissue homogenization

Step 1.

Place insects in a 1.7 mL tube with two 3.2 mm steel balls (RNase-free).

Snap freeze in liquid nitrogen.



REAGENTS

✓ 3.2 mm stainless steel beads, RNase free NEXSSB32-RNA by Contributed by users

Tissue homogenization

Step 2.

The goal of this step is to pulverise the insects without allowing them to thaw. Grinding time requires optimisation to balance DNA yield against physical damage to the DNA.

Grind insects using a Retsch mixer mill MM 400 / Qiagen TissueLyser or similar with **blocks pre-cooled in liquid N₂** for 90 seconds at 20 Hz.

After grinding, immediately return the tubes to liquid nitrogen. Ground samples can be stored at -80°C until you're ready to continue.

Lysis

Step 3.

Prepare the lysis buffer by adding 1.5 µL of RNase A (100 mg/mL) to 1438.5 µL of Buffer G2 per

sample.



REAGENTS

RNase A [19101](#) by [Qiagen](#)

Buffer G2 [1014636](#) by [Qiagen](#)

Lysis

Step 4.

Remove the powdered sample from liquid nitrogen or the -80°C freezer and immediately add 1440 µL of lysis buffer from step 3.

Lysis

Step 5.

Incubate for 30 minutes at 37°C with inversion.



TEMPERATURE

37 °C Additional info:

Lysis

Step 6.

Add 60 µL of Proteinase K (20 mg/mL) and incubate with inversion at 50°C for 2 hours.



TEMPERATURE

50 °C Additional info:



REAGENTS

Proteinase K [19131](#) by [Qiagen](#)

Pellet debris and recover supernatant

Step 7.

Centrifuge at max speed for 20 minutes to pellet the debris. Using a P1000, **slowly** pipette 1200 µL of the supernatant into a 15 mL Falcon tube.

Pellet debris and recover supernatant

Step 8.

Dilute the lysate to around 3 mL. The volume doesn't matter, but if the sample is more dilute it will run through the column faster.



ANNOTATIONS

Tom Harrop 07 May 2018

Dilute the sample by adding buffer G2 to the required volume.

DNA binding and washing

Step 9.

Equilibrate a QIAGEN Genomic-tip 20/G with 1 mL Buffer QBT. Wait for the Genomic-tip to drain by

gravity. You can do this during step 7.



REAGENTS

✓ Buffer QBT by Contributed by users

DNA binding and washing

Step 10.

Carefully apply the lysate to the Genomic-tip. Don't vortex the sample. Gently invert it a few times to mix it, then pour as much as you can from the 15 mL tube into the Genomic-tip. Pulse spin the 15 mL tube to collect any remaining lysate, and **slowly** pipette it into the Genomic-tip with a P1000. Wait for the Genomic-tip to drain by gravity. All of the draining steps can take a long time, and this step is particularly slow if the sample is viscous (30–60 minutes), but resist the urge to apply positive pressure.

DNA binding and washing

Step 11.

Wash the QIAGEN Genomic-tip with 1 mL Buffer QC.

Repeat this step 3 times for a total of 4 mL of Buffer QC.



REAGENTS

✓ Buffer QC by Contributed by users

DNA recovery

Step 12.

Elute the DNA into a new 15 mL tube with 2 mL of Buffer QF. Wait patiently for the tip to drain.



REAGENTS

✓ Buffer QF by Contributed by users

DNA recovery

Step 13.

After elution I split the sample into three because our lab doesn't have a fast centrifuge for 15 mL tubes.

Slowly pipette 667 μ L of the eluate into three clean 1.7 mL tubes with a P1000. Add 467 μ L of room-temperature isopropanol to each tube, mix by inversion about 10 times, and centrifuge for 20 min at 15,000g at 4°C to pellet DNA



REAGENTS

2-Propanol (IsoPropanol) [PC8601.SIZE.4L](#) by [Bio Basic Inc.](#)

DNA recovery

Step 14.

The pellet can be hard to see at this point. **Slowly** pour off the supernatant. I usually keep the supernatant from this step and the ethanol wash until I've tested the final yield of my preps. Pulse spin the tubes, and remove any remaining traces of supernatant with a P10.

DNA recovery

Step 15.

Add about 1 mL of ice-cold, 70% ethanol. Invert a few times, and centrifuge at 15,000g for 10 min at 4°C. Pour off the ethanol, pulse spin the tubes, and remove any remaining traces of ethanol.

Repeat this step once for a second wash with 70% ethanol.



REAGENTS



Ethanol by Contributed by users

DNA recovery

Step 16.

Air dry the pellet **briefly**. It can be difficult to resuspend, so don't leave it too long. I try to get all the drops of ethanol with a P10, and then watch the tubes until the traces evaporate.

DNA recovery

Step 17.

Resuspend each pellet in 55 µL of Qiagen EB (Tris-Cl pH 8). Don't pipette the sample up and down to resuspend it. Leave the tubes overnight at room temperature with inversion to resuspend.

Resuspension volume is up to you, but if you're doing a MinION run, this leaves you around 5 µL for QC.



REAGENTS

Buffer EB [19086](#) by [Qiagen](#)

Sample QC

Step 18.

Visualise 1 µL of sample to estimate the molecular weight. We don't have PFGE or a TapeStation / Fragment Analyzer, so I use 0.7% Agarose.


Quantify your sample with a Qubit. If you're using the HS kit, you may have to dilute the sample 1:10.

Analyse 1 µL in a UV spectrophotometer (e.g. Nanodrop).

Results

Step 19.

This gel shows 10 samples extracted with this method. Each sample is one third of a single sample that was run through a Genomic-tip and split into three at the precipitation stage. Each original sample was a pool of 20 adult [Argentine stem weevils](#) (Coleoptera: Curculionidae). DNA extractions from these insects with other methods have resulted in low 260/230 ratios.

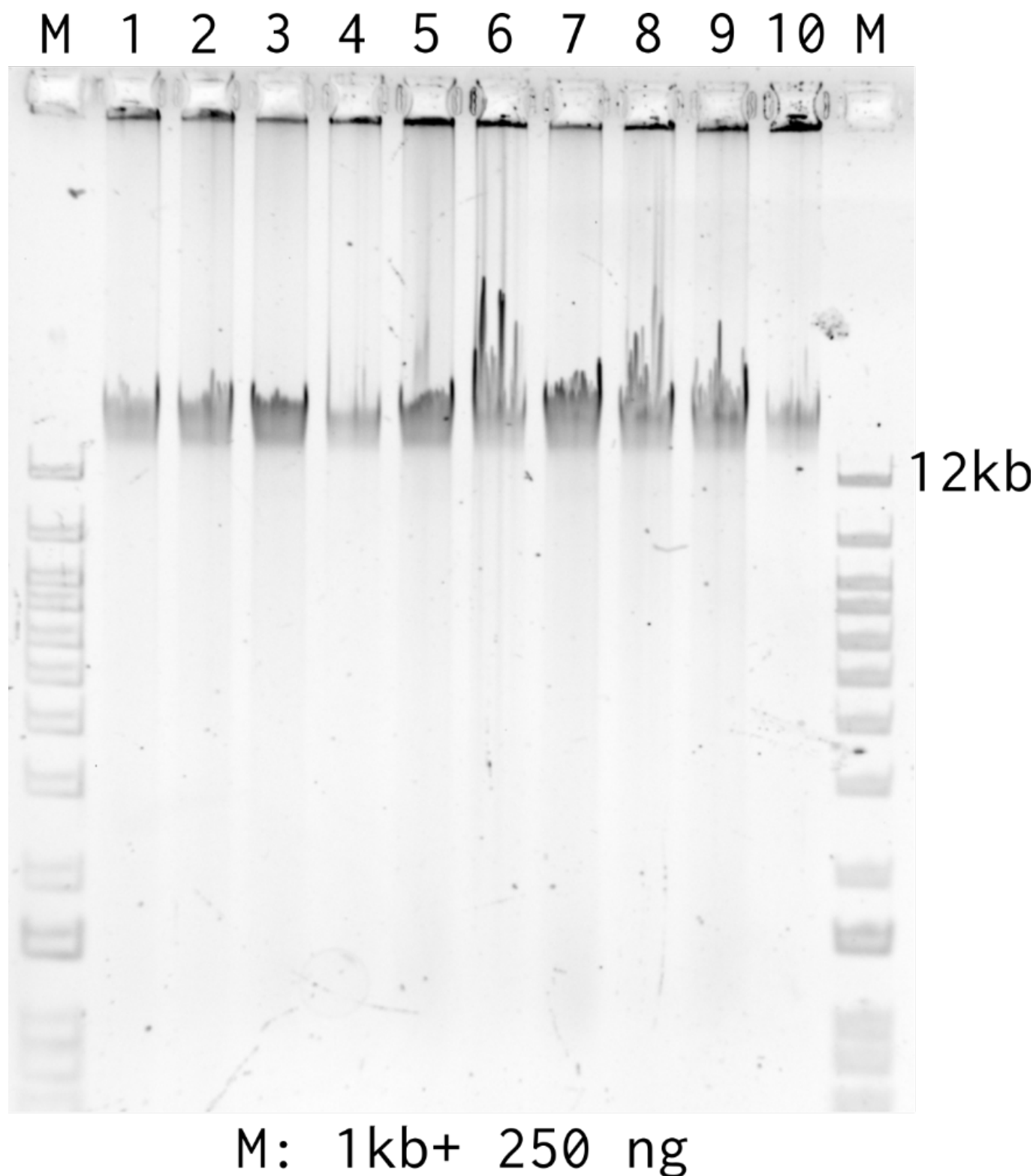

0.7 % Agarose in TBE, run at 80 V for 60 minutes and post-stained for 30 minutes in 0.1 mg / mL EtBr. We ran 1 μ L of each sample.

Results

Step 20.

These are the qubit and nanodrop results for the same 10 samples.

ID	Qubit (ng/μL)	Nanodrop (ng/μL)	260/280	260/230
1	52	35.4	1.87	1.96
2	56	113.1	1.85	2.24
3	59	122.7	1.88	2.28
4	30.4	61.4	1.8	2.1
5	53	54	1.78	2.19
6	53.3	104.2	1.81	2.28
7	79.6	214.6	1.78	2.18
8	132	133.7	1.88	2.41
9	210	127.7	1.85	2.36
10	114	233.8	1.8	2.18

Results

Step 21.

We used all of sample 8 for a library prep with the Oxford Nanopore SQK-LSK108 kit. At the end of the prep we recovered 1536 ng of DNA. We loaded the whole library onto a R9.4.1 flowcell and ran it for 40 hours.

Here are some metrics from the basecalling:

GB data called	GB >Q7	N50 KB >Q7	Mean KB >Q7	Median KB >Q7	Max KB >Q7	GB >20KB >Q7	GB >50KB >Q7	GB >100KB >Q7	
5.5	5.2	21.4	8.6	3.7		216	2.8	0.8	0.1

Here is a weighted read length histogram:

