

Yeast Colony PCR: It doesn't get any easier than this!

Joe Horecka and Angela M. Chu

Abstract

A protocol for yeast colony PCR, starting with intact cells. Can be used to analyze genomic and plasmid DNA, with PCR products up to 2 kb being no problem. It doesn't get easier than this: pick a small amount of cells from a plate into water, lyse cells at 99C for 5 min. Use the lysed cells as template by adding to 1/10th final volume in a PCR. We use this protocol for diagnostics (Phire HS II enzyme) and for Sanger sequencing (Phusion HS II or NEB Q5 enzyme). The main reasons why others might not get this protocol to work could be due to using too many cells and/or not using pure water for the cell lysis. Good luck, and have fun!


See the 'BEFORE START' section of this protocol for photographs that show how to pick cells from a colony and make a suspension. See the 'WARNINGS' section for important points and trouble shooting.

Citation: Joe Horecka and Angela M. Chu Yeast Colony PCR: It doesn't get any easier than this!. **protocols.io**
dx.doi.org/10.17504/protocols.io.gzwbx7e


Published: 31 Jan 2017

Before start

Below are a few photographs of picking cells from a colony and making cell suspensions. Illustrated in the last photo are three cell suspension densities, one of them being just right. We believe a problem of our protocol failing for some researchers is using too many cells, perhaps with the idea of 'the more the better.' If you are having problems, try using less cells.

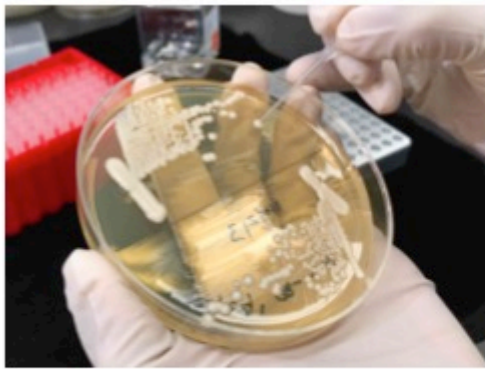

Picking cells from a colony

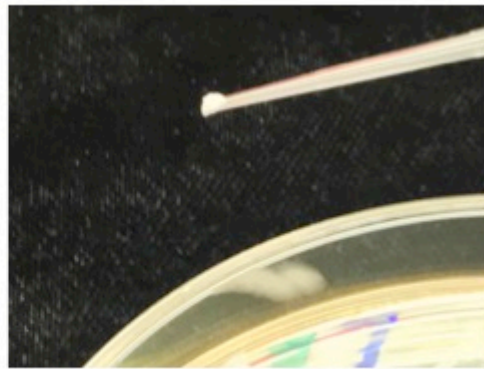

Cells on a small tip



Suspend cells in water



Cell suspensions:
1 = water only; 2 = just right;
3 = too few; 4 = too many.

Protocol

Preparation of heat-popped yeast cells

Step 1.

1. Aliquot 50 μ L water per each 0.2 ml PCR tube.
2. Use a small pipet tip to pick a small amount of yeast (less than the size of sesame seed). Avoid touching the agar. Transfer the yeast to the water, swirling to get cells off the tip. Cap tubes. Gently vortex to suspend cells.
3. Place the tubes in a PCR machine with program: 99C, 5 min; 4C, hold. Remove tubes from machine, quick-spin, place tubes on ice. The heat-popped yeast are ready to use or can be frozen at -20C for later with no loss of activity. *Important:* gently vortex cells into suspension just before adding to PCR reactions. DO NOT use supernatant from pelleted yeast cells. The DNA template is in the cells. Use the cell suspension.

PCR using heat-popped yeast cells - an example with 20 μ L PCRs and Phire Hot Start II enzyme

Step 2.

From here forward, PCR is essentially the same as with any template (e.g., plasmid, pure genomic DNA). Make a PCR master mix following the PCR enzyme manufacturer. Below are the final concentrations of reagents in a master mix for 20 μ L PCRs using Phire Hot Start II (Thermo Scientific):

PCR Master Mix (20 μ L reactions)

1X Phire HS II buffer

0.2 mM each dNTP

0.4 μ L Phire HS II enzyme/20 μ L PCR

0.5 μ M each primer (can be added later to individual tubes)

The master mix can be made with or without primers. Aliquot master mix to new 0.2 ml PCR tubes or to a 96-well plate. If primers are included in the mix, aliquot 18 μ L of the mix. If primers are added later, reduce volume/tube accordingly. In any case, the typical tube-filling order is: 1) master mix, 2) primers (if not in the master mix), and 3) 2 μ L vortexed, heat-popped yeast. Quick-spin samples using a mini centrifuge (PCR tubes) or a centrifuge with a plate rotor (PCR plates).

Run PCR (Phire HS II extension time is 10-15 sec/kb PCR product):

98C, 30 sec	denature
98C, 5 sec	
60C, 5 sec	35 cycles PCR
72C, X sec	
72C, 1 min	final extension
4C, hold	end

Remove tubes from PCR machine and quick-spin. If using an E-Gel, add 20 μ L water to each tube, vortex, and load 20 μ L to each E-Gel lane. For other gels, add concentrated loading dye and run on the gel.

Warnings

Surprisingly, some researchers can't get this protocol to work. The three most common mistakes are

likely:

1) Using too many cells. See photographs in 'BEFORE START' for cell suspension density. If you're having problems, try using less cells.

2) Suspending cells in something other than molecular biology grade water. Don't use zymolyase, NaOH, etc. Use water only!

3) Centrifuging the heat-popped cells and using the clear supernatant. This protocol lyses cell membranes with water and heat. It's not a nucleic acid extraction. THE DNA IS IN THE SUSPENDED, LYSED CELLS. Use the cell suspension.