

## Protocol for use with Purified mRNA or rRNA Depleted RNA and NEBNext Ultra II Directional RNA Library Prep Kit for Illumina (E7760, E7765)

New England Biolabs<sup>1</sup>

<sup>1</sup>New England Biolabs

dx.doi.org/10.17504/protocols.io.t8xerxn

New England Biolabs (NEB)

Tech. support phone: +1(800)632-7799 email: [info@neb.com](mailto:info@neb.com)

### ABSTRACT

The NEBNext Ultra II Directional RNA Library Prep Kit for Illumina contains the enzymes and buffers required to convert a broad range of input amounts of RNA into high quality directional (strand-specific) libraries for next-generation sequencing on the Illumina platform. The fast, user-friendly workflow has minimal hands-on time and is compatible with poly(A) mRNA enrichment and rRNA depletion methods.

### PROTOCOL STATUS

**Working**

### GUIDELINES

#### Section 4

#### RNA Sample Requirements

This Section can be used for libraries without any enrichment or depletion of total RNA with RIN scores > 7.

#### RNA Integrity:

RNA Integrity Number (RIN) is computed using ribosomal RNA (rRNA) amount in the sample. If rRNA is removed by any method, the RIN value should not be used to evaluate the integrity of the RNA sample. In this case, we recommend that the fragmentation time is empirically determined if the RNA sample is suspected to be low quality. **The following recommendation apply to the total RNA samples only.**

Assess the quality of the input RNA by running the RNA sample on an Agilent Bioanalyzer RNA 6000 Nano/Pico Chip to determine the

RNA Integrity Number (RIN). RNA with different RIN values require different fragmentation times or no fragmentation at all.

For intact (RIN > 7) or partially degraded RNA samples (RIN = 2 to 7) follow the library preparation protocol in Section 4 (current Section). See Table 4.1.1 for the recommended the fragmentation times.









For highly degraded samples (RIN = 1 to 2) (e.g. FFPE), which do not require fragmentation, follow the library preparation protocol in Section 5.

#### RNA Sample Requirements:














The RNA sample should be free of salts (e.g.  $Mg^{2+}$ , or guanidinium salts, divalent cation chelating agents (e.g. EDTA or EGTA) or organics (e.g. phenol or ethanol). RNA must be free of DNA. gDNA is a common contaminant from RNA preps. It may be carried over from the interphase of organic extractions or when the silica matrix of solid phase RNA purification methods is overloaded. If the total RNA sample may contain gDNA contamination, treat the sample with DNase I to remove all traces of DNA (DNase is not provided in this kit). After treatment with DNase I the enzyme should be removed from the sample. DNase I can be removed from the extraction using phenol/chloroform extraction and ethanol precipitation.

### MATERIALS

| NAME ▾                  | CATALOG # ▾ | VENDOR ▾  |
|-------------------------|-------------|---|
| NEBNext RNase H         | E6318       |  New England Biolabs |
| RNase H Reaction Buffer | E6312       |  New England Biolabs |

| NAME ▾   | CATALOG # ▾ | VENDOR ▾  |
|--|-------------|---|
| NEBNext rRNA Depletion Solution  | E6313       |  New England Biolabs |
| NEBNext Probe Hybridization Buffer   | E6314       |  New England Biolabs |
| DNase I (RNase-free)   | E6316       |  New England Biolabs |
| DNase I Reaction Buffer  | E6315       |  New England Biolabs |
| Nuclease-free Water  | E6317       |  New England Biolabs |
| NEBNext RNA Sample Purification Beads  | E6315       |  New England Biolabs |
| Magnetic Rack  | View        |   |
| 80% Ethanol (freshly prepared)   | View        |   |
| Thermal cycler   | View        |   |
| Agencourt RNAClean XP Beads  | A63987      | Beckman Coulter   |
| DNase I (e.g., NEB #M0303) and DNase I Cleanup Reagents or Kit for Removal of DNA Prior to Depletion | View        |  New England Biolabs |
| Random Primers   | E7422       |  New England Biolabs |

#### STEPS MATERIALS

| NAME ▾  | CATALOG # ▾ | VENDOR ▾  |
|---|-------------|---|
| NEBNext Strand Specificity Reagent                            | E7766       |  New England Biolabs |
| NEBNext First Strand Synthesis Enzyme Mix                     | E7761       |  New England Biolabs |
| NEBNext Second Strand Synthesis Reaction Buffer with dUTP Mix | E7426       |  New England Biolabs |
| NEBNext Second Strand Synthesis Enzyme Mix                    | E7425       |  New England Biolabs |
| Nuclease-free Water   | E7764       |  New England Biolabs |
| NEBNext Sample Purification Beads                             | E7767       |  New England Biolabs |
| Fresh 80% Ethanol   |             |   |
| (0.1X) TE Buffer  | E7763       |  New England Biolabs |
| NEBNext Ultra II End Prep Reaction Buffer                     | E7647       |  New England Biolabs |
| NEBNext Ultra II End Prep Enzyme Mix                          | E7646       |  New England Biolabs |
| NEBNext Ligation Enhancer                                     | E7374       |  New England Biolabs |
| NEBNext Ultra II Ligation Master Mix                          | E7648       |  New England Biolabs |
| NEBNext USER Enzyme   | E7458       |  New England Biolabs |
| NEBNext Sample Purification Beads                             | E7767       |  New England Biolabs |
| 80% Ethanol (freshly prepared)                                |             |   |

| NAME ▾                            | CATALOG # ▾ | VENDOR ▾  |
|-----------------------------------|-------------|---|
| NEBNext Sample Purification Beads | E6315       |  New England Biolabs |
| 80% Ethanol (freshly prepared)    |             |   |

#### SAFETY WARNINGS

Please refer to the SDS (Safety Data Sheet) for safety warnings and hazard information.

#### BEFORE STARTING

##### Input Amount Requirements

1 ng – 100 ng purified mRNA or rRNA depleted RNA that is **quantified after the purification**. RNA should be DNA free in up to 5 µl of Nuclease-free Water, quantified by Qubit Fluorometer and quality checked by Bioanalyzer.

The protocol is optimized for approximately 200 nt RNA inserts. To generate libraries with longer RNA insert sizes, refer to Appendix A, (Section 6) for recommended fragmentation times and size selection conditions.


This protocol has been optimized using Universal Human Reference Total RNA.

#### RNA Fragmentation and Priming

- [ ! ] RNA fragmentation is only required for intact or partially degraded RNA. Recommended fragmentation times can be found in Table 1 (step 3).**

Assemble the fragmentation and priming reaction **on ice** in a nuclease-free tube by adding the following components:

| Fragmentation and Priming Mix                  | Volume |
|--|--------|
| Purified mRNA or rRNA Depleted RNA             | 5 µl   |
| NEBNext First Strand Synthesis Reaction Buffer | 4 µl   |
| Random Primers                                 | 1 µl   |
| Total Volume                                   | 10 µl  |

- Mix thoroughly by pipetting up and down 10 times.
- Place the sample in a thermocycler and incubate the sample at  **94 °C** following the recommendations in Table 1 below for fragment sizes ~200 nt.

**Table 1.** Suggested fragmentation times based on RIN value of RNA input.

| RNA Type               | RIN | Frag. Time     |
|------------------------|-----|----------------|
| Intact RNA             | > 7 | 15 min @ 94°C  |
| Partially Degraded RNA | 2–6 | 7–8 min @ 94°C |

#### NOTE

**Refer to Appendix A, (Chapter 6) for fragmentation conditions if you are preparing libraries with large inserts (> 200 bp). Conditions in Appendix A, (Chapter 6) only apply for intact RNA.**

- Immediately transfer the tube to ice and proceed to First Strand cDNA Synthesis.

## First Strand cDNA Synthesis

- 5 Assemble the first strand synthesis reaction **on ice** by adding the following components to the fragmented and primed RNA from Step 4:

| First Strand Synthesis Reaction           | Volume       |
|---|--------------|
| Fragmented and Primed RNA (Step 4)        | 10 µl        |
| NEBNext Strand Specificity Reagent        | 8 µl         |
| NEBNext First Strand Synthesis Enzyme Mix | 2 µl         |
| <b>Total Volume</b>                       | <b>20 µl</b> |



NEBNext Strand Specificity Reagent

by [New England Biolabs](#)

Catalog #: [E7766](#)



NEBNext First Strand Synthesis Enzyme

Mix

by [New England Biolabs](#)

Catalog #: [E7761](#)

- 6 Mix thoroughly by pipetting up and down 10 times.

- 7 [!] Incubate the sample in a preheated thermocycler with the heated lid set at  $\geq 80^{\circ}\text{C}$  as follows:

**Note:** If you are following recommendations in Appendix A (Chapter 6), for libraries with longer inserts (> 200 bases), increase the incubation at  $42^{\circ}\text{C}$  from 15 minutes to 50 minutes at Step 2 below.

Step 1: 00:10:00 at  $25^{\circ}\text{C}$

Step 2: 00:15:00 at  $42^{\circ}\text{C}$

Step 3: 00:15:00 at  $70^{\circ}\text{C}$

Step 4: Hold at  $4^{\circ}\text{C}$

- 8 Proceed directly to Second Strand cDNA Synthesis.

## Second Strand cDNA Synthesis

- 9 Assemble the second strand cDNA synthesis reaction **on ice** by adding the following components into the first strand synthesis product from Step 8).

| Second Strand Synthesis Reaction                                | Volume       |
|---|--------------|
| First Strand Synthesis Product (Step 8)                         | 20 µl        |
| NEBNext Second Strand Synthesis Reaction Buffer with dUTP (10X) | 8 µl         |
| NEBNext Second Strand Synthesis Enzyme Mix                      | 4 µl         |
| Nuclease-free Water   | 48 µl        |
| <b>Total Volume</b>   | <b>80 µl</b> |




NEBNext Second Strand Synthesis


Reaction Buffer with dUTP Mix

by [New England Biolabs](#)

Catalog #: [E7426](#)



**NEBNext Second Strand Synthesis Enzyme Mix**  
by [New England Biolabs](#)  
Catalog #: [E7425](#)




**Nuclease-free Water**  
by [New England Biolabs](#)  
Catalog #: [E7764](#)

10 Keeping the tube on ice, mix thoroughly by pipetting up and down at least 10 times.



11 Incubate in a thermocycler for  **01:00:00** at  **16 °C** with the heated lid set at ≤  **40 °C** (or off).

#### Purification of Double-stranded cDNA Using SPRIselect Beads or NEBNext Sample Purification Beads

12 Vortex SPRIselect Beads or NEBNext Sample Purification Beads to resuspend.





**NEBNext Sample Purification Beads**  
by [New England Biolabs](#)  
Catalog #: [E7767](#)


13 Add  **144 µl** (1.8X) of resuspended beads to the second strand synthesis reaction (~  **80 µl**). Mix well on a vortex mixer or by pipetting up and down at least 10 times.

14 Incubate for  **00:05:00** at room temperature.

15 Briefly spin the tube in a microcentrifuge to collect any sample on the sides of the tube. Place the tube on a magnet to separate beads from the supernatant. After the solution is clear, carefully remove and discard the supernatant. Be careful not to disturb the beads, which contain DNA.

**Caution: Do not discard beads.**

16 Add  **200 µl** of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for  **00:00:30**, and then carefully remove and discard the supernatant.





**Fresh 80% Ethanol**


17 Repeat Step 16 once for a total of 2 washing steps.


 **go to step #16 Repeat Step**

- 18 Air dry the beads for up to 5 minutes while the tube is on the magnetic rack with lid open.

**Caution: Do not over-dry the beads. This may result in lower recovery of DNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.**

- 19 Remove the tube from the magnetic rack. Elute the DNA from the beads by adding  53 µl 0.1X TE Buffer (provided) to the beads. Mix well on a vortex mixer or by pipetting up and down at least 10 times. Quickly spin the tube and incubate for  00:02:00 at room temperature. Place the tube on the magnetic rack until the solution is clear.

 (0.1X) TE Buffer  
by New England Biolabs  
Catalog #: E7763

- 20 Remove  50 µl of the supernatant and transfer to a clean nuclease-free PCR tube.



**NOTE**

**If you need to stop at this point in the protocol samples can be stored at -20°C.**

## End Prep of cDNA Library

- 21 Assemble the end prep reaction on ice by adding the following components to the second strand synthesis product from Step 20.

| End Prep Reaction                         | Volume       |
|---|--------------|
| Second Strand Synthesis Product (Step 20) | 50 µl        |
| NEBNext Ultra II End Prep Reaction Buffer | 7 µl         |
| NEBNext Ultra II End Prep Enzyme Mix      | 3 µl         |
| <b>Total Volume</b>                       | <b>60 µl</b> |

If a master mix is made, add  0 µl of master mix to  50 µl of cDNA for the End Prep reaction.

 NEBNext Ultra II End Prep Reaction Buffer  
by New England Biolabs  
Catalog #: E7647

 NEBNext Ultra II End Prep Enzyme Mix  
by New England Biolabs  
Catalog #: E7646

- 22 Set a 100 µl or 200 µl pipette to 50 µl and then pipette the entire volume up and down at least 10 times to mix thoroughly. Perform a quick spin to collect all liquid from the sides of the tube.

**NOTE**

**It is important to mix well. The presence of a small amount of bubbles will not interfere with performance.**

23 Incubate the sample in a thermocycler with the heated lid set at  $\geq 75^{\circ}\text{C}$  as follows.

00:30:00 at  $20^{\circ}\text{C}$

00:30:00 at  $65^{\circ}\text{C}$

Hold at  $4^{\circ}\text{C}$

24 Proceed immediately to Adaptor Ligation.

### Adaptor Ligation

25 **[ ! ] Dilute the red NEBNext Adaptor\* prior to setting up the ligation reaction in ice-cold Adaptor Dilution Buffer and keep the adaptor on ice.**

| Purified RNA | Dilution Required                           |
|--------------|---|
| 100 ng–11 ng | 5–fold dilution in Adaptor Dilution Buffer  |
| 10 ng–1 ng   | 25–fold dilution in Adaptor Dilution Buffer |

\*The NEBNext adaptor is provided in NEBNext oligos kit. NEB has several oligo kit options, which are supplied separately from the library prep kit.

26 Assemble the ligation reaction **on ice** by adding the following components, in the order given, to the end prep reaction product from Step 24.

| Ligation Reaction                    | Volume Per One Library               |
|--------------------------------------|--------------------------------------|
| End Prepped DNA (Step 24)            | 60 $\mu\text{l}$                     |
| Diluted Adaptor (Step 25)            | 2.5 $\mu\text{l}$                    |
| NEBNext Ligation Enhancer            | 1 $\mu\text{l}$                      |
| NEBNext Ultra II Ligation Master Mix | 30 $\mu\text{l}$                     |
| <b>Total Volume</b>                  | <b>93.5 <math>\mu\text{l}</math></b> |

#### NOTE

**The Ligation Master Mix and Ligation Enhancer can be mixed ahead of time and is stable for at least 8 hours @  $4^{\circ}\text{C}$ . We do not recommend premixing the Ligation Master Mix, Ligation Enhancer and adaptor prior to use in the Adaptor Ligation Step.**



NEBNext Ligation Enhancer  
by [New England Biolabs](#)  
Catalog #: [E7374](#)






NEBNext Ultra II Ligation Master Mix  
by [New England Biolabs](#)  
Catalog #: [E7648](#)

27 Set a 100  $\mu\text{l}$  or 200  $\mu\text{l}$  pipette to 80  $\mu\text{l}$  and then pipette the entire volume up and down at least 10 times to mix thoroughly. Perform a quick spin to collect all liquid from the sides of the tube.

**[ ! ] Caution: The NEBNext Ultra II Ligation Master Mix is very viscous. Care should be taken to ensure adequate mixing of the ligation reaction, as incomplete mixing will result in reduced ligation efficiency. The presence of a small amount of bubbles will not interfere with performance.**

28 Incubate 00:15:00 at  $20^{\circ}\text{C}$  in a thermocycler.

29 Add  3 µl (blue) USER™ Enzyme to the ligation mixture from Step 60, resulting in total volume of  96.5 µl




NEBNext USER Enzyme  
by [New England Biolabs](#)  
Catalog #: [E7458](#)


30 Mix well and incubate at  37 °C for  00:15:00 with the heated lid set to ≥  45 °C

31 Proceed immediately to Purification of the Ligation Reaction.

#### Purification of the Ligation Reaction Using SPRIselect Beads or NEBNext Sample Purification Beads


32 **[ ! ] Note: If you are selecting for libraries with larger insert size (> 200 nt) follow the size selection recommendations in Appendix A, Chapter 6.**

Add  87 µl (0.9X) resuspended SPRIselect Beads or NEBNext Sample Purification Beads and mix well on a vortex mixer or by pipetting up and down at least 10 times.






NEBNext Sample Purification Beads  
by [New England Biolabs](#)  
Catalog #: [E7767](#)

33 Incubate for  00:10:00 at room temperature.

34 Quickly spin the tube in a microcentrifuge and place the tube on an appropriate magnetic rack to separate beads from the supernatant. After the solution is clear (~  00:05:00 ), discard the supernatant that contains unwanted fragments.

**Caution: Do not discard beads.**

35 Add  200 µl of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for  00:00:30 , and then carefully remove and discard the supernatant.



80% Ethanol (freshly prepared)

36 Repeat Step 35 once for a total of 2 washing steps.


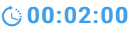
 [go to step #35 Repeat Step](#)


37 Briefly spin the tube, and put the tube back in the magnetic rack.



- 38 Completely remove the residual ethanol, and air dry beads until the beads are dry for up to 5 minutes while the tube is on the magnetic rack with the lid open.

**Caution: Do not over-dry the beads. This may result in lower recovery of DNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.**

- 39 Remove the tube from the magnetic rack. Elute DNA target from the beads by adding  17 µl 0.1X TE (provided) to the beads. Mix well on a vortex mixer or by pipetting up and down. Quickly spin the tube and incubate for  00:02:00 at room temperature. Put the tube in the magnet until the solution is clear.

- 40 Without disturbing the bead pellet, transfer  15 µl of the supernatant to a clean PCR tube and proceed to PCR enrichment.

**NOTE**

**If you need to stop at this point in the protocol, samples can be stored at –20°C.**

### PCR Enrichment of Adaptor Ligated DNA

- 41 [!] Check and verify that the concentration of your oligos is 10 µM on the label.

[!] Use **Option A** for any NEBNext oligos kit where index primers are supplied in tubes. These kits have the forward and reverse primers supplied in **separate** tubes.

Use **Option B** for any NEBNext oligos kit where index primers are supplied in a 96-well plate format. These kits have the forward and reverse (i7 and i5) primers **combined**.

- 42 Set up the PCR reaction as described below based on the type of oligos (PCR primers) used.

#### Option A: Forward and Reverse Primers Separate:

| Component                          | Volume Per One Library |
|------------------------------------|------------------------|
| Adaptor Ligated DNA (Step 40)      | 15 µl                  |
| NEBNext Ultra II Q5 Master Mix     | 25 µl                  |
| Universal PCR Primer/i5 Primer*,** | 5 µl                   |
| Index (X) Primer/i7 Primer*,**     | 5 µl                   |
| <b>Total Volume</b>                | <b>50 µl</b>           |

#### Option B: Forward and Reverse Primers Combined:

| Component                      | Volume Per One Library |
|--------------------------------|------------------------|
| Adaptor ligated DNA (Step 40)  | 15 µl                  |
| NEBNext Ultra II Q5 Master Mix | 25 µl                  |
| Index (X)/i7 Primer Mix*       | 10 µl                  |
| <b>Total Volume</b>            | <b>50 µl</b>           |

\* NEBNext Oligos must be purchased separately from the library prep kit. Refer to the corresponding NEBNext Oligo kit manual for determining valid barcode combinations.

\*\* Use only one i7 primer/ index primer per sample. Use only one i5 primer (or the universal primer for single index kits) per sample

- 43 Mix well by gently pipetting up and down 10 times. Quickly spin the tube in a microcentrifuge.

- 44 Place the tube on a thermocycler with the heated lid set to  105 °C and perform PCR amplification using the following PCR cycling conditions (refer to Table 44.A and Table 44.B):

**Table 44.A:**

| Cycle Step | Temp | Time | Cycles |
|------------|------|------|--------|
|------------|------|------|--------|

|                      |      |            |          |
|----------------------|------|------------|----------|
| Initial Denaturation | 98°C | 30 seconds | 1        |
| Denaturation         | 98°C | 10 seconds | 6–13*,** |
| Annealing/Extension  | 65°C | 75 seconds |          |
| Final Extension      | 65°C | 5 minutes  | 1        |
| Hold                 | 4°C  | ∞          |          |

\* The number of PCR cycles should be adjusted based on RNA input.

\*\* It is important to limit the number of PCR cycles to avoid overamplification. If overamplification occurs, a second peak ~ 1,000 bp will appear on the Bioanalyzer trace (see Figure 7.2 in manual).

**Table 44.B:** Recommended PCR cycles based on total RNA input amount:



| Total RNA Input | Recommended PCR Cycles |
|-----------------|------------------------|
| 1,000 ng        | 6–7                    |
| 50 ng           | 7–8                    |
| 10 ng           | 9–10                   |
| 1 ng            | 12–13                  |

**NOTE**

**PCR cycles are recommended based on high quality Universal Human Reference Total RNA. It may require optimization based on the sample quality to prevent PCR over-amplification.**

## Purification of the PCR Reaction using SPRIselect Beads or NEBNext Sample Purification Beads

45 Vortex SPRIselect Beads or NEBNext Sample Purification Beads to resuspend.

46 Add  **45 µl** (0.9X) of resuspended beads to the PCR reaction (~  **50 µl** ). Mix well on a vortex mixer or by pipetting up and down at least 10 times.




**NEBNext Sample Purification Beads**



by [New England Biolabs](#)

Catalog #: [E6315](#)

47 Incubate for  **00:05:00** at room temperature.

48 Quickly spin the tube in a microcentrifuge and place the tube on an appropriate magnetic rack to separate beads from the supernatant. After the solution is clear (~  **00:05:00** ), carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

**Caution: Do not discard beads.**

49 Add  **200 µl** of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for  **00:00:30** , and then carefully remove and discard the supernatant.





80% Ethanol (freshly prepared)

50 Repeat Step 49 once for a total of 2 washing steps.

[go to step #49 Repeat Step](#)


51 Air dry the beads for up to 5 minutes while the tube is on the magnetic rack with the lid open.

**Caution: Do not over-dry the beads. This may result in lower recovery of DNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.**

52 Remove the tube from the magnetic rack. Elute the DNA target from the beads by adding  23 µl 0.1X TE (provided) to the beads. Mix well on a vortex mixer or by pipetting up and down ten times. Quickly spin the tube in a microcentrifuge and incubate for  00:02:00 at room temperature. Place the tube in the magnetic rack until the solution is clear.

53 Transfer  20 µl of the supernatant to a clean PCR tube, and store at -20°C.

#### Assess Library Quality on an Agilent Bioanalyzer DNA Chip


54 Run  1 µl library on a DNA 1000 chip. If the library yield is too low to quantify on this chip, please run the samples on a DNA High Sensitivity chip. A dilution may be necessary for running on a Bioanalyzer High Sensitivity DNA Chip.

55 Check that the electropherogram shows a narrow distribution with a peak size approximately 300 bp.

#### NOTE

**If a peak at ~ 80 bp (primers) or 128 bp (adaptor-dimer) is visible in the bioanalyzer traces, bring up the sample volume (from Step 53) to 50 µl with 0.1X TE buffer and repeat the SPRIselect Bead or NEBNext Sample Purification Bead Cleanup Step (Section "Purification of the PCR Reaction using SPRIselect Beads or NEBNext Sample Purification Beads").**

56 Figure 56: Example of RNA library size distribution on a Bioanalyzer.

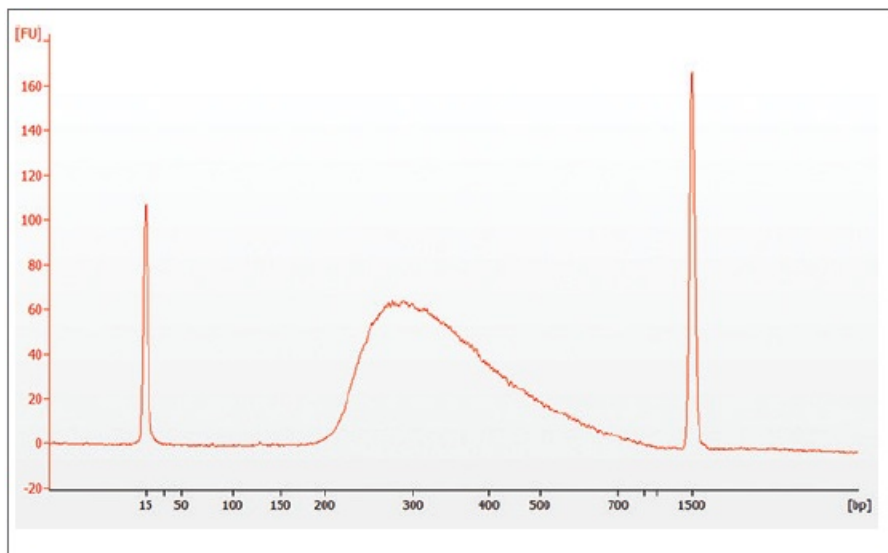

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited