

Fluorometric real-time quantification of in vivo fluorescence in *E. coli*

Dennis Dienst, Alice Pawlowski, Ilka Maria Axmann

Abstract

Citation: Dennis Dienst, Alice Pawlowski, Ilka Maria Axmann Fluorometric real-time quantification of in vivo fluorescence in *E. coli*. **protocols.io**

dx.doi.org/10.17504/protocols.io.nrxdd7n

Published: 31 May 2018

Materials

- ✓ Microplate Reader Synergy Mx [View](#) by Contributed by users
- ✓ Costar 96 clear bottom black side plates by Contributed by users
- LB Broth (Lennox) vegetal [0155](#) by [Carl Roth](#)
- Mineral oil for Molecular Biology [HP50](#) by [Carl Roth](#)

Protocol

Overnight culture preparation

Step 1.

Prepare o/n culture of each 3 independent clones

- add 5 mL of LB medium to a sterile culture tube
- add the appropriate amount of antibiotic(s)
- use a sterile inoculation needle, pipette tip or toothpick to pick a colony from the Agar plate
- gently whirl the tip in the medium (or drop it)
- incubate at 37 °C and 230 rpm

🔧 TEMPERATURE

37 °C Additional info: 230 rpm

🔗 LINK:

<https://www.addgene.org/protocols/inoculate-bacterial-culture/>



Inoculation of 96-well plates

Step 2.

1:100 dilution of o/n cultures in microtiter (96-well) plate

- pipet 99 μL fresh LB medium (incl. the appropriate antibiotic) into each well of the 96-well plate
- add 1 μL of each o/n culture to the medium
- add 100 μL LB medium to 6 wells as blank
- overlay each well with 40 μL of mineral oil (to prevent evaporation)

96-well plate:
clear bottom/ black side

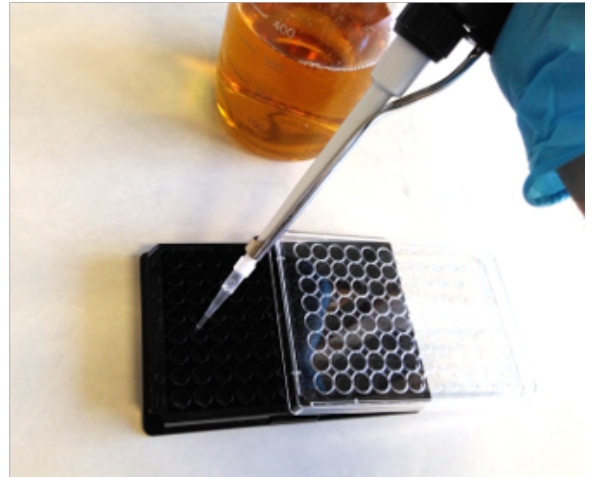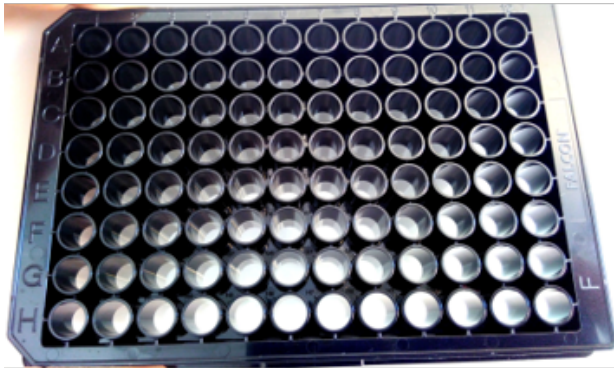

📌 NOTES

Dennis Dienst 13 Apr 2018

- pipet in a sterile or semi-sterile environment, i.e. in sterile hood, since also contaminants would have the chance to grow during the following o/n incubation
- due to the autofluorescence of LB medium, you might prefer to use a mineral medium like [M9 medium](#) - note that the choice of media will significantly influence the growth characteristics and might require longer cultivation periods

Plate reader settings and measurement

Step 3.

- remove the lid from the microtiter plate
- insert microtiter plate into plate reader device (e.g. Synergy Mx, Biotek)

Settings

Gain: 50

Absorbance: **OD₆₀₀** (8 measurements per data point)

Fluorescence **GFP**: Excitation λ - **480/ 9 nm**

Emission λ - **520/ 20 nm**

Fluorescence **YFP**: Excitation λ - **500/ 9 nm**

Emission λ - **535/ 20 nm**

Frequency of measurements: **5 min** (241 reads, 10 measurements per data point)

Temperature: **37 °C**

Agitation: **continuous** 'medium mode'

Time: **20 h**

TEMPERATURE

37 °C Additional info: 'medium mode'

Data collection and basic analysis

Step 4.

- export data to .csv/.xlsx file
- subtract blank (LB) values from both the fluorescence and the OD₆₀₀ values
- estimate relative fluorescence by dividing GFP/YFP-fluorescence by the corresponding OD₆₀₀
- the figure below depicts exemplary curves of relative YFP (mVenus) fluorescence accumulation over time in *E. coli* strain W3110
 - the unit of measurement is

relative mVenus fluorescence [Ex.:Em. 500:535nm * OD₆₀₀⁻¹]

- as a quantitative proxy for the expression dynamics calculate the slopes of relative fluorescence within the linear phase of relative fluorophor accumulation

- divide single values of relative fluorescence by the corresponding value of the control group at a given time point, or a given OD₆₀₀

