

Nov 28,
2019

Marchantia Cryopreservation of Gemmae

Forked from [Marchantia cryopreservation of gemmae](#)

[Marius Rebmann](#)¹, [marta tomaselli](#)², [Linda Silvestri](#)³, [Susana Sauret-Gueto](#)¹, Michelle Lim⁴

¹Plant Sciences, University of Cambridge, OpenPlant, ²University of Cambridge, Open Plant, ³University of Cambridge, ⁴Plant Sciences, University of Cambridge, OpenPlant,

1 Works for me [dx.doi.org/10.17504/protocols.io.9vxh67n](https://doi.org/10.17504/protocols.io.9vxh67n)

OpenPlant Project



Marius Rebmann
Plant Sciences, University of Cambridge, OpenPlant



ABSTRACT

Simplified cryopreservation protocol for *Marchantia polymorpha* gemmae, based on ([Tanaka et al. 2016, Plant and Cell Physiology](#)). Enables long term storage of gemmae at -80°C.

MATERIALS

NAME	CATALOG #	VENDOR
Agar	A20021	Melford
Abcisic acid	A1049	Sigma Aldrich
Sodium Alginate	S1320	Duchefa Biochemie
Calcium chloride	22328.262	VWR Scientific
Sucrose	10634932	Fisher Scientific
Glycerol	G5516	Sigma Aldrich
Gamborg B5 Medium with Vitamins	G0210	Duchefa Biochemie

MATERIALS TEXT

- Petridish, 9cm/4.5cm or 6/12-well transparent multi-well plates
- Liquid nitrogen (N₂)

Reagent setup:

Preculture plates

1/2 strength Gamborg B5, pH 5.8 + 1.2% agar + 0.3 M sucrose + 10µM ABA

Autoclave media without ABA, add ABA from a filter sterilised stock solution just before pouring plates. 9cm/4.5cm petridishes or 6/12-well transparent multi-well plates can be used for plates, depending on desired throughput.

Alginate Solution

1/2 strength Gamborg B5 + 3% Sodium alginate

CaCl₂ Solution

1/2 strength Gamborg B5 + 0.1M CaCl₂

Loading solution

1/2 strength Gamborg B5 + 2M glycerol, + 1M sucrose

Thawing solution

1/2 strength Gamborg B5 + 1M sucrose

Rinse solution


1/2 strength Gamborg B5 + 1% sucrose

Preculture

- 1 Collect fresh gemmae and plate on **preculture plates**, incubate 1-3 days under normal growth conditions (e.g. 21°C constant light)


Encapsulation

- 2 Place small drops (approx. 40µL) of **alginate solution** on an empty 4.5cm petridish/multi-well plate to form beads.

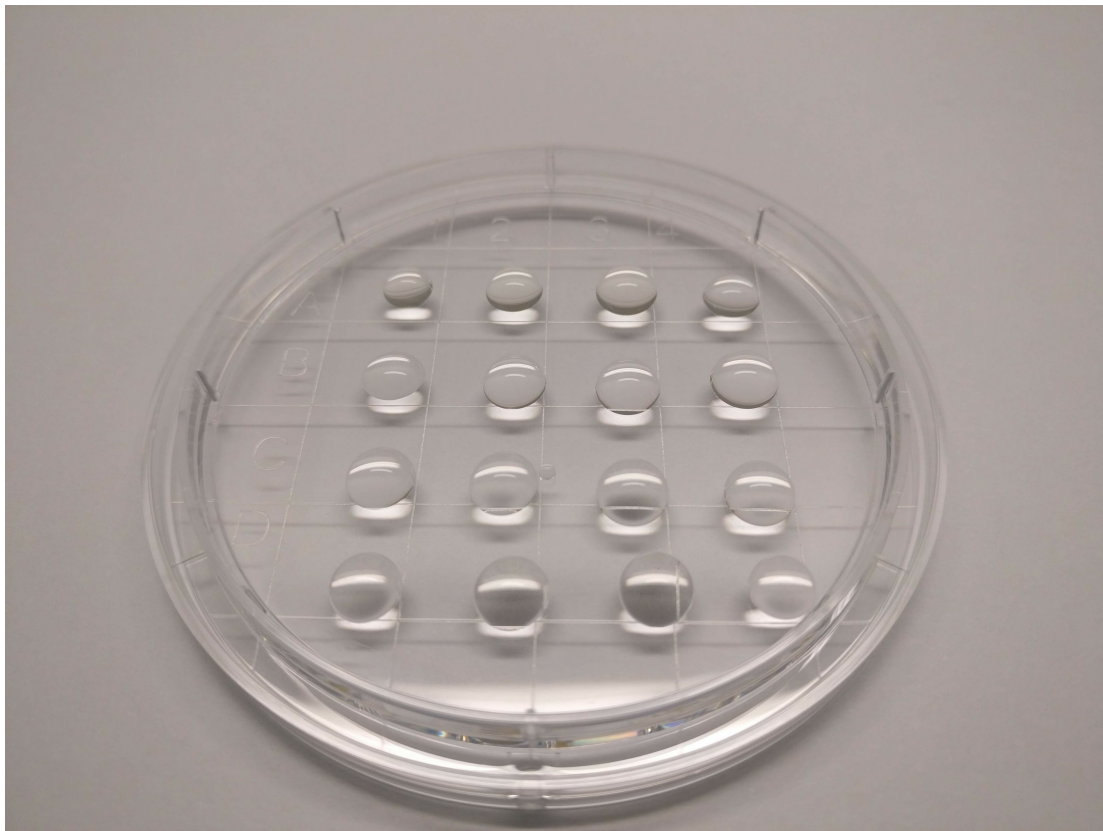


- 3 Use forceps to transfer gemmae from preculture plates, place 1-5 gemmae inside each bead.

Proceed quickly to step 4 to avoid gemmae sinking to the bottom or edge of the bead.

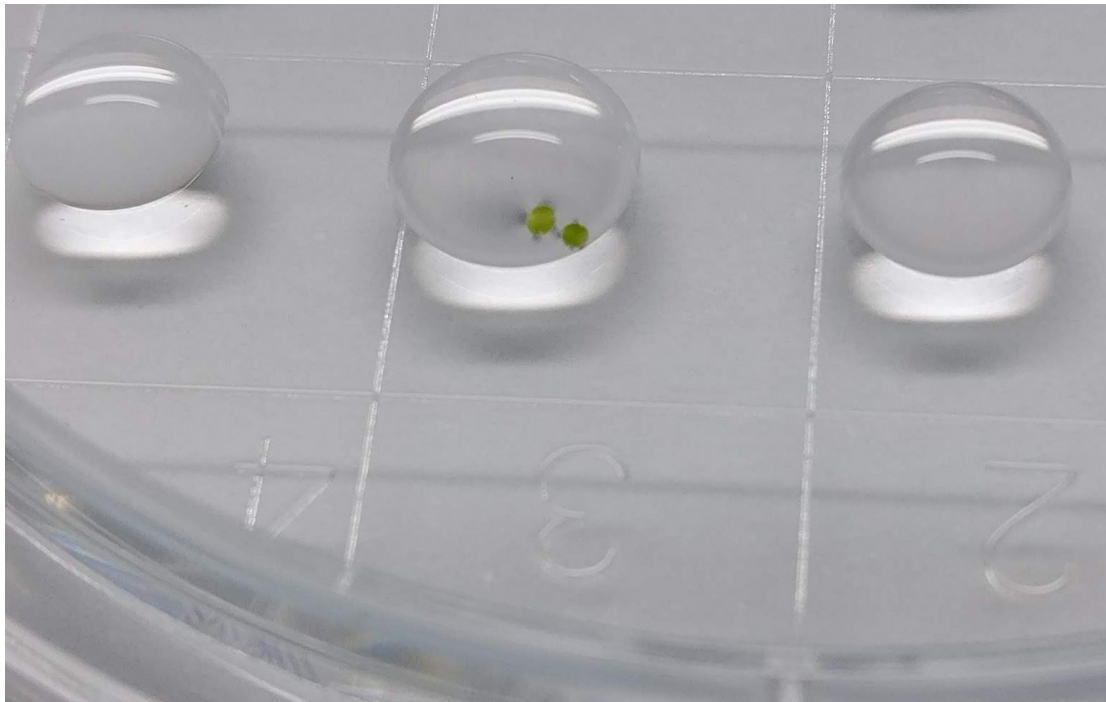

- 4 Add a drop of **CaCl₂ solution** to each bead, make sure to not displace gemmae from the bead.


- 4.1 Allow beads to solidify for 10min.

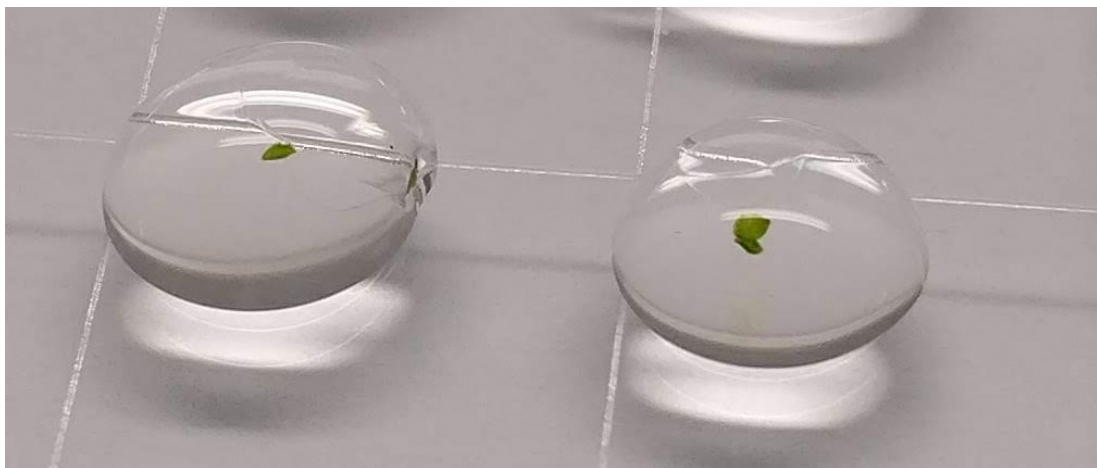

- 5 Use a serological pipette to fill the plate with **loading solution**, submerging the beads. Soak beads for 30min.

6 Use a serological pipette to remove the loading solution without disrupting the beads.

6.1 Remove lid and air dry beads in flowhood for >2h.

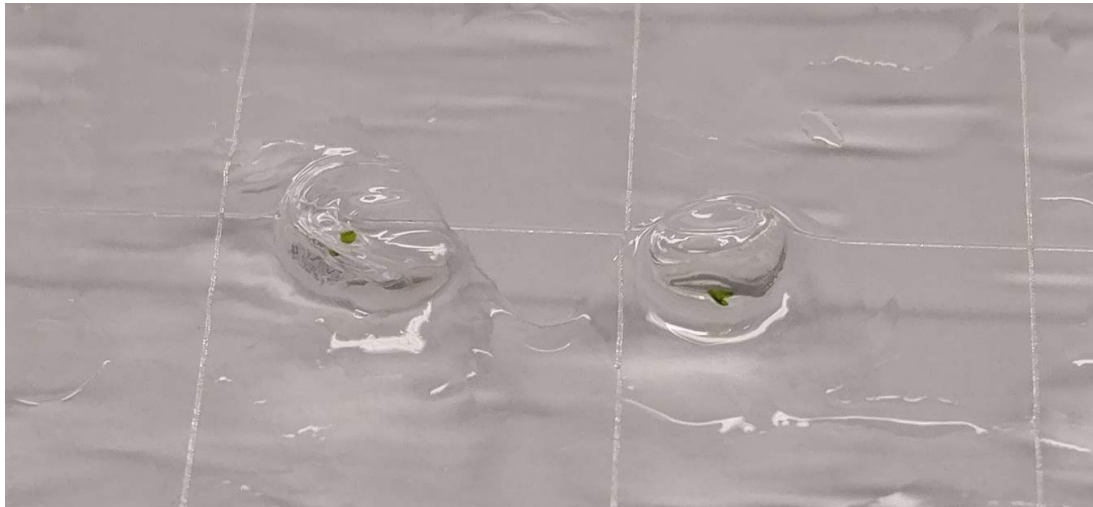

7 Use forceps or scalpels to transfer beads into 1.5mL Eppendorf tubes (max. 10 beads/tube).

7.1 Flash freeze tubes in liquid nitrogen for >2 minutes.

7.2 Move tubes to -80°C for long term storage.

Thawing and Recovery

8 Remove tubes from freezer and immediately place in a 37°C water bath or heat block for 2 minutes.

9 Add 1.5mL **thawing solution** to each tube, incubate at room temperature for 10 minutes.

10 Replace thawing solution with **rinse solution**. Incubate at room temperature for 10 minutes.

11 Remove rinse solution, transfer beads on 1/2 Gamborg plates, incubate under normal growth conditions (e.g. 21°C constant light) until thalli grow



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited