



Sep 20, 2019

Agarose Gel Electrophoresis

¹¹Northeast Forest University**1** *Works for me* dx.doi.org/10.17504/protocols.io.7h3hj8n**2019 iGEM NEFU_China**Tech. support email: shengyiyanywork@gmail.com

ABSTRACT

The concentration of the agarose gel is 1%.


MATERIALS TEXT

Agarose
50X TAE
Gelred

SAFETY WARNINGS

Please wear gloves for the experiment, Be careful when heating, protect your hands.

- 1 Weigh 0.20 g of agarose.
- 2 Add 20 ml of TAE 1X.
Need to prepare 1X TAE with 50X TAE.
- 3 Heat up until the solution is homogeneous.
- 4 Add 1 μ l Gelred into solution and mix well.
- 5 Put the solution into bed for polymerize, make sure "comb" is well placed and the solution is balanced.
- 6 20 min later, use the gel.

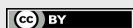

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited