



## ITS2 metabarcoding protocol

Forked from [Pollen metabarcoding](#)

[Tomasz Suchan](#)<sup>1</sup>

<sup>1</sup>W. Szafer Institute of Botany, Polish Academy of Sciences

[dx.doi.org/10.17504/protocols.io.sn9edh6](https://doi.org/10.17504/protocols.io.sn9edh6)

Molecular Biogeography Group



Tomasz Suchan

W. Szafer Institute of Botany,...



### PROTOCOL STATUS

#### Working

We use this protocol in our group and it is working

### GUIDELINES

Perform reactions in small batches until you are confident that there is no cross-contamination among the samples. Including isolation blanks and PCR blanks is crucial for the quality control.

#### ITS2 primers used in the 1st PCR:

ITS2-4R GTGACTGGAGTT CAGACGTGTGCTCTTCCGATCTNNNNNNTCCTCCGCTTATTGATATGC

ITS2-S2F ACACTCTTCCCTACACGACGCTCTTCCGATCTNNNNNNATGCGATACTTGGTGTGAAT

#### Indexing primers used in the 2nd PCR (xxxxxxx - index):

AATGATACGGCGACCACCGAGATCTACACxxxxxxxACACTCTTCCCTACACGACGC

CAAGCAGAAGACGGCATACGAGATxxxxxxxGTGACTGGAGTT CAGACGTGTGC

PCR\_F\_A501: AATGATACGGCGACCACCGAGATCTACACT**GAACCTT**ACACTCTTCCCTACACGACGC

PCR\_F\_A502: AATGATACGGCGACCACCGAGATCTACACT**GCTAAGT**ACACTCTTCCCTACACGACGC

PCR\_F\_A503: AATGATACGGCGACCACCGAGATCTACACT**GTTCTCT**ACACTCTTCCCTACACGACGC

PCR\_F\_A504: AATGATACGGCGACCACCGAGATCTACACT**AAGACAC**ACACTCTTCCCTACACGACGC

PCR\_F\_A505: AATGATACGGCGACCACCGAGATCTACAC**CTAATCGA**ACACTCTTCCCTACACGACGC

PCR\_F\_A506: AATGATACGGCGACCACCGAGATCTACAC**CTAGAACA**ACACTCTTCCCTACACGACGC

PCR\_F\_A507: AATGATACGGCGACCACCGAGATCTACACT**AAGTTCC**ACACTCTTCCCTACACGACGC

PCR\_F\_A508: AATGATACGGCGACCACCGAGATCTACACT**AGACCTA**ACACTCTTCCCTACACGACGC

PCR\_F\_D501: AATGATACGGCGACCACCGAGATCTACACT**ATAGCCT**ACACTCTTCCCTACACGACGC

PCR\_F\_D502: AATGATACGGCGACCACCGAGATCTACAC**ATAGAGGC**ACACTCTTCCCTACACGACGC

PCR\_F\_D503: AATGATACGGCGACCACCGAGATCTACAC**CTATCCT**ACACTCTTCCCTACACGACGC

PCR\_F\_D504: AATGATACGGCGACCACCGAGATCTACAC**GGCTCTGA**ACACTCTTCCCTACACGACGC

PCR\_F\_D505: AATGATACGGCGACCACCGAGATCTACAC**AGGCGAAG**ACACTCTTCCCTACACGACGC

PCR\_F\_D506: AATGATACGGCGACCACCGAGATCTACACT**AATCTTA**ACACTCTTCCCTACACGACGC

PCR\_F\_D507: AATGATACGGCGACCACCGAGATCTACAC**CAGGACGT**ACACTCTTCCCTACACGACGC

PCR\_F\_D508: AATGATACGGCGACCACCGAGATCTACAC**GTACTGAC**ACACTCTTCCCTACACGACGC

PCR\_R\_A701: CAAGCAGAAGACGGCATACGAGAT**GTCTGTGAT**GTGACTGGAGTT CAGACGTGTGC

PCR\_R\_A702: CAAGCAGAAGACGGCATACGAGAT**ACCCTGT**GTGACTGGAGTT CAGACGTGTGC

PCR\_R\_A703: CAAGCAGAAGACGGCATACGAGAT**TGGATCTG**GTGACTGGAGTT CAGACGTGTGC

PCR\_R\_A704: CAAGCAGAAGACGGCATACGAGAT**CCGTTTGT**GTGACTGGAGTT CAGACGTGTGC

PCR\_R\_A705: CAAGCAGAAGACGGCATACGAGAT**TGCTGGGT**GTGACTGGAGTT CAGACGTGTGC

PCR\_R\_A706: CAAGCAGAAGACGGCATACGAGAT**GAGGGTT**GTGACTGGAGTT CAGACGTGTGC

PCR\_R\_A707: CAAGCAGAAGACGGCATACGAGAT**AGGTTGGG**GTGACTGGAGTT CAGACGTGTGC


PCR\_R\_A708: CAAGCAGAAGACGGCATACGAGAT**GTGTGGTGGT**GTGACTGGAGTT CAGACGTGTGC

PCR\_R\_A709: CAAGCAGAAGACGGCATACGAGAT**TGGGTTTC**GTGACTGGAGTT CAGACGTGTGC

PCR\_R\_A710: CAAGCAGAAGACGGCATACGAGATT**GGT CACAGT**GACTGGAGTTCAGACGTGTGC  
PCR\_R\_A711: CAAGCAGAAGACGGCATACGAGATT**TGACCCT**GTGACTGGAGTTCAGACGTGTGC  
PCR\_R\_A712: CAAGCAGAAGACGGCATACGAGAT**CCACT CCT**GTGACTGGAGTTCAGACGTGTGC

PCR\_R\_D701: CAAGCAGAAGACGGCATACGAGAT**CGAGT AAT**GTGACTGGAGTTCAGACGTGTGC  
PCR\_R\_D702: CAAGCAGAAGACGGCATACGAGATT**CT CCGGAGT**GACTGGAGTTCAGACGTGTGC  
PCR\_R\_D703: CAAGCAGAAGACGGCATACGAGAT**AATGAGCGG**TGACTGGAGTTCAGACGTGTGC  
PCR\_R\_D704: CAAGCAGAAGACGGCATACGAGAT**GGAAT CTC**GTGACTGGAGTTCAGACGTGTGC  
PCR\_R\_D705: CAAGCAGAAGACGGCATACGAGATT**TTCTGAAT**GTGACTGGAGTTCAGACGTGTGC  
PCR\_R\_D706: CAAGCAGAAGACGGCATACGAGAT**ACGAATT C**GTGACTGGAGTTCAGACGTGTGC  
PCR\_R\_D707: CAAGCAGAAGACGGCATACGAGAT**AGCTT CAG**GTGACTGGAGTTCAGACGTGTGC  
PCR\_R\_D708: CAAGCAGAAGACGGCATACGAGAT**GCGCATT AGT**GACTGGAGTTCAGACGTGTGC  
PCR\_R\_D709: CAAGCAGAAGACGGCATACGAGAT**CATAGCCG**GTGACTGGAGTTCAGACGTGTGC  
PCR\_R\_D710: CAAGCAGAAGACGGCATACGAGATT**TCGCGGAGT**GACTGGAGTTCAGACGTGTGC  
PCR\_R\_D711: CAAGCAGAAGACGGCATACGAGAT**GCGCGAGAGT**GACTGGAGTTCAGACGTGTGC  
PCR\_R\_D712: CAAGCAGAAGACGGCATACGAGAT**CTATCGCT**GTGACTGGAGTTCAGACGTGTGC

MATERIALS

NAME	CATALOG #	VENDOR
Q5 Hot Start High-Fidelity DNA Polymerase - 100 units	M0493S	 New England Biolabs
Water, nuclease free		
dNTP mix (25 mM of each)		
Phire Plant Direct PCR Kit	F130WH	Thermo Fisher Scientific

BEFORE STARTING  
Prepare 5 µM primer solutions.

1st PCR

- 1 Prepare the mix:

11.64 µl molecular grade water

4 µl Q5 polymerase buffer

1 µl ITS2-S2F primer

1 µl ITS2-4R primer

0.2 µl Q5 Hot Start polymerase
- 2 Add 2 µl of the sample to 18 µl of the mix. Use water instead of the sample for the blanks.

**SAFETY INFORMATION**

Add one blank sample at this step (= "PCR blank") and also use the blank from the DNA extraction step (= "extraction blank").

<http://null>
- 3 Run the PCR program: initial denaturation at 98°C for 5 min; 20 cycles of denaturation at 98°C for 40 s, annealing at 49°C for 40 s and elongation at 72°C for 40 s; followed by a final extension step at 72°C for 5 min.

## Purification

4 Perform AMPure purification with ratio 1x. Elute in 10 µl.

### ✓ **PROTOCOL**



AMPure purification protocol

by **Tomasz Suchan**,

**W. Szafer Institute of Botany, Polish Academy of Sciences**

PREVIEW

START EXPERIMENT

4.1 Remove beads from the fridge, equilibrate to the room temperature, mix thoroughly but do not vortex.

4.2 Add desired ratio of AMPure beads to the DNA to be purified and mix well by pipetting.

4.3 Incubate 5 minutes.

4.4 Place on the magnetic rack.

4.5 Let it stand for 5 minutes on the rack, aspirate and discard supernatant.

4.6 Add 200 µl of freshly prepared 70% EtOH, incubate 30 seconds, aspirate and discard EtOH.

4.7 Repeat the wash: add 200 µl of freshly prepared 70% EtOH, incubate 30 seconds, aspirate and discard EtOH.

4.8 Wait until the pellet dries out completely. DO NOT remove the plate/tubes from the magnetic rack, as the dried AMPure might pop out from the tubes!

4.9 Add desired volume of Tris 10 mM or water, wetting the dried AMPure pellets (add 1 µl to the final volume to avoid pipetting out the beads)

4.10 Remove from the magnetic rack.

- 4.11 Resuspend by pipetting or vortexing.
- 4.12 Incubate 10 minutes, incubating in 37°C can improve DNA yield.
- 4.13 Place on the magnetic rack.
- 4.14 Let it stand for 5 minutes, pipette out and save supernatant. The eluted DNA is in the supernatant, do not discard it!

## 2nd PCR

- 5 Prepare the mix:

 4.82 µl molecular grade water

 2 µl Q5 reaction buffer

 0.08 µl dNTPs (25 mM each)

 0.1 µl Q5 Hot Start polymerase


- 6 Add 1 µl of the template to 7 µl of the mix.
- 7 Add 1 µl of each 5 uM primer (forward and reverse).

### ⚠ SAFETY INFORMATION

Use different pair of indexed primers for each sample.

- 8 Run the PCR program: 30 s denaturation at 98°C; 12 cycles of denaturation at 98 °C for 10 s, combined annealing and extension at 72°C for 30 s (shuttle PCR); the final extension at 72°C for 5 min.
- 9 Check profiles on a gel.

### EXPECTED RESULT

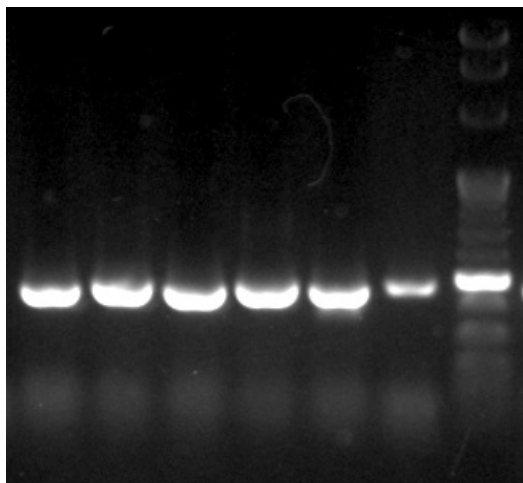

## Pooling

- 10 Pool all the samples together.

### NOTE

The sample concentrations can be normalized at this step, depending on the experimental questions. Purify the samples separately using AMPure, quantify and pool in the equimolar proportions. Alternatively you can use bead-based normalization, e.g. doi: [10.1186/1471-2164-15-645](https://doi.org/10.1186/1471-2164-15-645).

- 11 Perform AMPure purification with ratio 1x. Elute in 100 µl.

### PROTOCOL



AMPure purification protocol

by **Tomasz Suchan**,

**W. Szafer Institute of Botany, Polish Academy of Sciences**

PREVIEW

START EXPERIMENT

- 11.1 Remove beads from the fridge, equilibrate to the room temperature, mix thoroughly but do not vortex.

- 11.2 Add desired ratio of AMPure beads to the DNA to be purified and mix well by pipetting.

- 11.3 Incubate 5 minutes.

- 11.4 Place on the magnetic rack.

- 11.5 Let it stand for 5 minutes on the rack, aspirate and discard supernatant.

- 11.6 Add 200 µl of freshly prepared 70% EtOH, incubate 30 seconds, aspirate and discard EtOH.
- 11.7 Repeat the wash: add 200 µl of freshly prepared 70% EtOH, incubate 30 seconds, aspirate and discard EtOH.
- 11.8 Wait until the pellet dries out completely. DO NOT remove the plate/tubes from the magnetic rack, as the dried AMPure might pop out from the tubes!
- 11.9 Add desired volume of Tris 10 mM or water, wetting the dried AMPure pellets (add 1 µl to the final volume to avoid pipetting out the beads)
- 1.10 Remove from the magnetic rack.
- 1.11 Resuspend by pipetting or vortexing.
- 1.12 Incubate 10 minutes, incubating in 37°C can improve DNA yield.
- 1.13 Place on the magnetic rack.
- 1.14 Let it stand for 5 minutes, pipette out and save supernatant. The eluted DNA is in the supernatant, do not discard it!
- 12 Check the concentration using Qubit and the profile using TapeStation/Fragment Analyzer. Calculate molarity from that and proceed to the sequencing. Add 15% PhiX to the sequencing run.



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited