



May 09,  
2019

Working

## U Cinn - Triglyceride Assay [↗](#)

Patrick Tso<sup>1</sup>, Dana Lee<sup>1</sup>

<sup>1</sup>University of Cincinnati

dx.doi.org/10.17504/protocols.io.xm4fk8w

Mouse Metabolic Phenotyping Centers

Tech. support email: [info@mmpc.org](mailto:info@mmpc.org)

Lili Liang

### ABSTRACT

#### Summary:

Determinations of triglycerides in plasma/serum/lymph will be made using a Randox Triglycerides colorimetric kit. The triglycerides are determined after enzymatic hydrolysis with lipases.

### EXTERNAL LINK

<https://mmpc.org/shared/document.aspx?id=200&docType=Protocol>

### MATERIALS

NAME	CATALOG #	VENDOR	CAS NUMBER	RRID
Triglyceride Assay Kit	TR213	<a href="#">Randox</a>		

### MATERIALS TEXT

#### Working Reagent:

*Reagents and Materials:*

**Enzyme Reagent R1b**

**Buffer R1a**


#### Procedure:

Prepare **Working Reagent** by reconstituting one vial of **Enzyme Reagent R1b** with a portion of **Buffer R1a** and then transfer entire contents to bottle **R1a**, rinsing vial **R1b** several times.

#### Note:

**Randox Life Sciences**, [RRID:SCR\\_005525](#)

- 1 Prepare **Working Standards** by making a serial dilution of the stock 200mg/dl standard. (See diagram below) \*Stock standard is included in kit

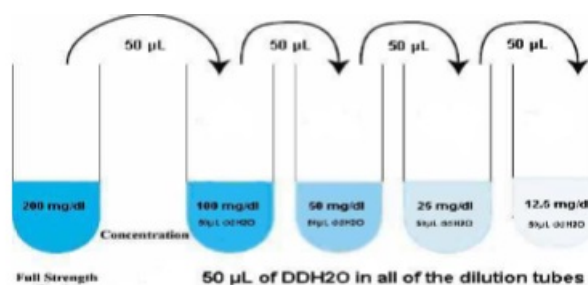

- 2 Prepare **Working Reagent** by reconstituting one vial of **Enzyme Reagent R1b** with a portion of **Buffer R1a** and then transfer entire contents to bottle **R1a**, rinsing vial **R1b** several times.
- 3 Using a 96 well flat bottom plate, into separate wells, pipette 2µL of deionized water, standard, or sample to be assayed.
- 4 Add 200µL of reconstituted **Working Reagent** to all wells.
- 5 Incubate plate for 5 minutes at 37°C
- 6 Determine the absorbance (abs) of the standards and of each unknown at 500nm.
- 7 Calculate values of unknowns from the standard curve.

**Specimen:** Serum or Plasma. Specimen stable for 7 days at 2-8°C or 3 months at -20°C.

**Assay Linearity:** 1172 mg/dl

**Working Reagent Stability After Reconstitution:** 21 days stored at 2-8°C. PROTECT FROM LIGHT

**Stability of Final Reaction:** 60 minutes



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited