



Oct 03, 2019

MojoSort™ Human anti-APC Nanobeads Protocol [↗](#)Sam Li¹¹BioLegend

1

Works for me

[dx.doi.org/10.17504/protocols.io.7x2hpgq](https://doi.org/10.17504/protocols.io.7x2hpgq)

BioLegend

Sam Li
BioLegend

ABSTRACT

Product description and procedure summary:

Target cells are positively selected or depleted by incubating the sample with an anti-human APC conjugated antibody, followed by incubation with magnetic anti-APC Nanobeads. The magnetically labeled fraction is retained by the use of a magnetic separator. These are the APC+ cells, do not discard them if those are the cells of interest. Some of the downstream applications include functional assays, gene expression, phenotypic characterization, etc.

Note: This procedure is optimized for the isolation of 10^7 to 2×10^8 cells per tube. If working with fewer than 10^7 cells, keep volumes as indicated for 10^7 cells. For best results, optimize the conditions to your specific cell number and tissue. Prepare fresh MojoSort™ Buffer solution by diluting the 5X concentrate with sterile distilled water. *Scale up volumes if using 14 mL tubes and Magnet, and place the tube in the magnet for 10 minutes.*

EXTERNAL LINK

<https://www.biolegend.com/protocols/mojosort-human-anti-apc-nanobeads-protocol/4755/>

GUIDELINES

MojoSort™ magnetic particles can be used with other commercially available magnetic separators, both free standing magnets and column-based systems. Because MojoSort™ protocols are optimized for the MojoSort™ separator, the protocols may need to be adjusted for other systems. Please contact BioLegend Technical Service (tech@biolegend.com) for more information and guidance. We do not recommend using MojoSort™ particles for BD's IMag™ or Life Technologies' DynaMag™.

Application notes: To use this product in magnetic separation columns, a titration of the Nanobeads should be performed. Optimal concentration for magnetic separation columns is lot-specific. Please contact BioLegend Technical Service (tech@biolegend.com) for further assistance on how to use MojoSort™ Nanobeads in magnetic separation columns.

MATERIALS

NAME ▾	CATALOG # ▾	VENDOR ▾
MojoSort™ Magnet	480019	BioLegend
MojoSort™ Human anti-APC Nanobeads	480089	BioLegend
MojoSort™ Buffer	480017	BioLegend

MATERIALS TEXT

- Adjustable pipettes
- 70µm filters (one per sample)
- 5mL (12 x 75mm) or 14mL (17 x 100mm) polypropylene tubes
- Reagents for sample preparation
- Reagents and instruments (Flow cytometer) to determine yield and purity


- 1 Prepare cells from your tissue of interest. Kits for human samples have been optimized for PBMCs, please prepare the cells using a suitable method.
- 2 In the final wash of your sample preparation, resuspend the cells in MojoSort™ Buffer by adding up to 4 mL in a 5 mL (12 x 75 mm) polypropylene tube.
Note: Keep MojoSort™ Buffer on ice throughout the procedure.
- 3 Filter the cells with a 70 µm cell strainer, centrifuge at 300xg for 5 minutes, and resuspend in a small volume of MojoSort™ Buffer. Count and adjust the cell concentration to 1 x 10⁸ cells/mL by adding MojoSort™ Buffer.
- 4 Aliquot 100 µL of cell suspension (10⁷ cells) into a new tube. **Add 5 µL of Human TruStain FcX™ (Fc Receptor Blocking Solution)**, mix well and **incubate at room temperature for 10 minutes**. Scale up the volume accordingly if separating more cells. For example, if the volume of Human TruStain FcX™ for 1x10⁷ cells is 5 µL, add 50 µL for 1 x 10⁸ cells. When working with less than 10⁷ cells, use indicated volumes for 10⁷ cells.
- 5 Check the recommended usage for flow cytometric staining of the APC-conjugated antibody indicated in the antibody technical datasheet. Calculate the volume to stain 10⁷ cells (or desired amount of cells). **Add the appropriate volume of APC-conjugated antibody** to the cell suspension, mix well and **incubate on ice for 15 minutes**.
Optional: Take an aliquot before adding the antibody to monitor purity and yield.
- 6 Wash the cells by adding MojoSort™ Buffer up to 4 mL; centrifuge the cells at 300xg for 5 minutes.
- 7 Discard supernatant and resuspend in 100 µL of MojoSort™ Buffer.
- 8 Resuspend the beads by vortexing, maximum speed, 5 touches. **Add the appropriate volume of Human anti-APC Nanobeads** to the cell suspension, mix well and **incubate on ice for 15 minutes**. The volume of MojoSort™ Human anti-APC Nanobeads should be adjusted depending on starting percentage of cell type to be isolated. For 1x10⁷ cells in 100 µl of buffer, use the following volumes:

Cell Type	Cell Frequency in PBMCs	Volume of MojoSort™ anti-APC Nanobeads (for 1x10 ⁷ cells in 100 µl of buffer)
CD3	50-60%	20 µL
CD4 (including monocytes)	50-55%	20 µL
CD8a (including low expressing cells)	15-20%	10 µL
CD19	5-10%	10 µL
TCR Vδ2	1-3%	5 µL

Scale up the volume accordingly if separating more cells. For example, if the volume of Nanobeads for 1x10⁷ cells is 10 µL, add 100 µL for 1 x 10⁸ cells. When working with less than 10⁷ cells, use indicated volumes for 10⁷ cells.

- 9 Wash the cells by adding MojoSort™ Buffer up to 4 mL; centrifuge the cells at 300xg for 5 minutes.
- 10 Discard the supernatant.

- 11 Add 2.5 mL of MojoSort™ Buffer.
Note: If you observe aggregates, filter the suspension. To maximize yield, you can disrupt the aggregates by pipetting the solution up and down.
- 12 Place the tube in the magnet for 5 minutes.
Optional: Take a small aliquot before placing the tube in the magnet to monitor purity and yield. Keep unused cells to be used as control or other applications if needed.
- 13 Pour out the unlabeled fraction. If these are your cells of interest, **DO NOT DISCARD**. Resuspend the labeled cells in 2.5 mL MojoSort™ Buffer.
- 14 Repeat steps 11-13 on the labeled fraction twice more for a total of **3 separations**. Pool the unlabeled fractions and keep the labeled cells. The fraction that is not of interest may be useful as staining controls, to monitor purity/yield, or other purposes.
Optional: Take a small aliquot to monitor purity and yield.

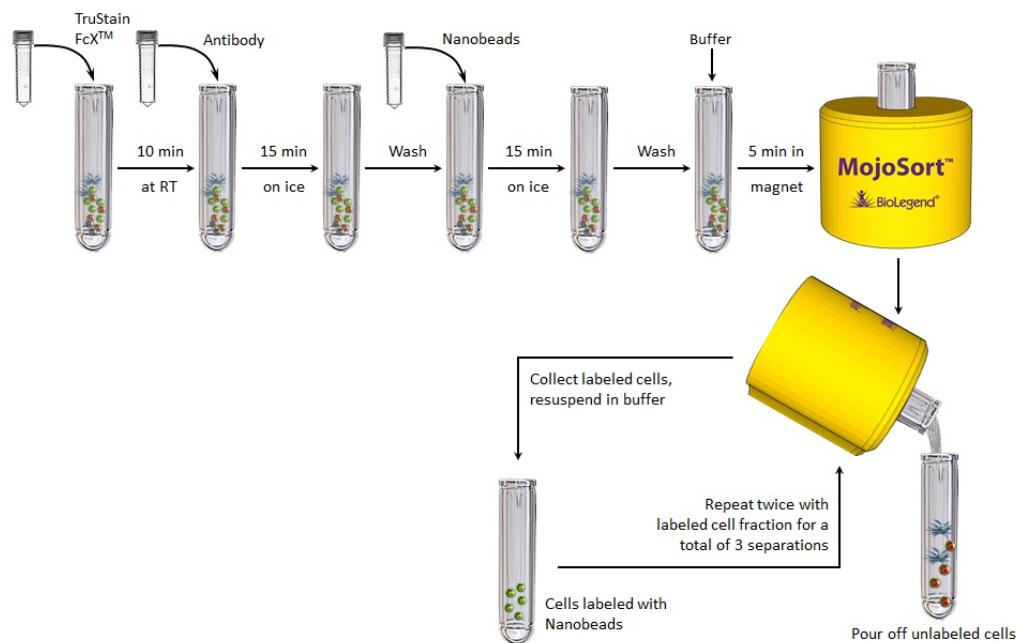

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited