



Dec 18, 2019

# Attachment 1: Preparation of Fluorescein Diacetate and Propidium Iodide Stock Solutions (FDA/PI)

Integrated Islet Distribution Program<sup>1</sup>, Carol J Swanson<sup>1</sup>, Barbara Olack<sup>1</sup>
<sup>1</sup>[Integrated Islet Distribution Program, City of Hope, Duarte, CA]

1 Works for me

dx.doi.org/10.17504/protocols.io.5bgg2jw



#### ABSTRACT

This Standard Operating Procedure is adapted from the work of the 'National Institutes of Health-Sponsored Clinical Islet

Transplantation Consortium Phase 3 Trial: Manufacture of a Complex Cellular Product at Eight Processing Facilities" following the SOP cited in the document 'Purified Human Pancreatic Islet - Viability Estimation of Islet Using Fluorescent Dyes (FDA/PI): Standard Operating Procedure of the NIH Clinical Islet Transplantation Consortium'

This SOP defines the procedure for assessment of viability of human isolated islet preparations, which include endocrine and exocrine tissue, for use in the National Institute of Diabetes and Digestive and Kidney Diseases (NIDDK) sponsored research in the Integrated Islet Distribution Program (IIDP). This protocol is written to assist the participating islet isolation centers and investigators who are part of this program.

Fluorescein Diacetate/ Propidium Iodide (FDA)/(PI) Viability Assay is a rapid fluorometric method to test the integrity of the plasma membrane simultaneously using inclusion and exclusion dyes; the assay differentiates between viable and nonviable cells and is, consequently, used for determination of viability of islet preparations.

#### **GUIDELINES**

Fluorescein Diacetate/ Propidium Iodide (FDA)/(PI) Viability Assay is a rapid fluorometric method to test the integrity of the plasma membrane simultaneously using inclusion and exclusion dyes; the assay differentiates between viable and nonviable cells and is, consequently, used for determination of viability of islet preparations.

- The inclusion dye is Fluorescein Diacetate (FDA) and the exclusion dye is Propidium Iodide (PI). The final concentrations are as follows:
  - ♦ FDA: 0.46 µM
  - ♦ PI: 14.34 μM
- Fluorescein Diacetate is a nonpolar ester, which passes through plasma membranes and is hydrolyzed by intracellular esterases to
  produce free fluorescein. The polar fluorescein is confined within cells with an intact plasma membrane and can be observed under
  appropriate excitation conditions. FDA functions as an inclusion dye, i.e., viable cells will appear bright green fluorescent using FDA.
- Propidium iodide functions as an exclusion dye that cannot penetrate living cells but readily enters dead or dying cells. Once PI
  penetrates through the cell membrane, it binds to nucleic acids and causes them to fluoresce bright orange/red. PI absorbs in green
  light and fluoresces orange/red.
- Both of the fluorescent dyes used in this assay are light sensitive and must be kept in the dark, covered with aluminum foil.
- The fluorescent dyes are temperature sensitive and must be stored as follows:
  - ♦ FDA: ≤ -20°C
  - ♦ PI: 2-8° C

**MATERIALS** 

NAME VENDOR VENDOR VENDOR Sigma Aldrich



NAME Y	CATALOG #	VENDOR V
Propidium Iodide	P4170	Sigma Aldrich
Acetone	179124	Sigma Aldrich
Thermo Scientific™ Nunc™ Cell Culture/Petri Dishes	174926	Fisher Scientific
Corning® Dulbecco's Phosphate-Buffered Saline 1X without calcium and magnesium	21-031-CM	Corning


#### MATERIALS TEXT



Snap Cap Microcentrifuge Tube or equivalent
Polypropylene Microcentrifuge Tube
Corning Costar Snap Cap Microcentrifuge
Tube
2 mL snap cap polypropylene micro tube

07200210 👄



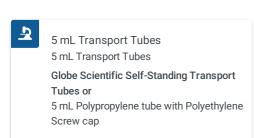

22-010-1227



### SAFETY WARNINGS

## Always wear gloves and observe standard chemical procedures:

Fluorescein Diacetate: FDA MSDSAction.pdf

Protect from light. Avoid contact and inhalation. Nitrile gloves are recommended in the MSDS when handling FDA.

Propidium Iodide: PI MSDSAction.pdf

 Use personal protective equipment. Product may be toxic if inhaled, swallowed, or splashed on skin. Avoid dust formation. Avoid breathing vapors, mist, or gas. Ensure adequate ventilation. Wear gloves and observe Safety Data Sheet. Suspected of causing genetic defects.

Acetone: AcetoneMSDSAction.pdf

Solvent/Flammable. Keep away from heat, spark, and open flame. Keep container tightly closed. Use with adequate ventilation. Avoid contact with eyes. Avoid prolonged or repeated breathing of vapor. Handle with gloves. Gloves must be inspected prior to use. Use proper glove removal technique (without touching glove's outer surface) to avoid skin contact with this product. Dispose of contaminated gloves after use in accordance with applicable laws and good laboratory practices. Wash and dry hands. Store aliquots in tightly sealed in glass or polypropylene tubes with polypropylene or polyethylene closure.

	Preparation	of Stock Flu	orescein Diacetate	(FDA)	) Solution
--	-------------	--------------	--------------------	-------	------------

1	Prepare [M] 24 Micromolar (µM) Fluorescein Diacetate (FDA) according to the formulation in sub-steps below.
	The expiration date, for both PI and FDA stains, is six months from the date of preparation.
	Both of the fluorescent dyes used in this assay are light sensitive and must be kept in the dark, covered with aluminum foil.
1.1	Calculation:
	FDA FW = 416.4 Stock Concentration = 24 µM Volume required = 200 mL
	Formula: FW X Concentration X Volume = $(416.4) \times (24 \times 10^{-6}) \times (200 \times 10^{-3}) = 0.00199 \text{ g FDA}$ Sigma Aldrich has a Mass Molarity Calculator that can be used to determine smaller quantities: https://www.sigmaaldrich.com/chemistry/stockroom-reagents/learning-center/technical-library/mass-molarity-
1.2	Dissolve 0.00199 g of FDA in 200 mL of acetone in a glass bottle and cover with aluminum foil.
	■ 0.00199 g FDA  FDA Supplier  FDA Lot #  Expiration Date
	■ 200 ml Acetone  Acetone Supplier  Acetone Lot #  Expiration Date
1.3	■ Store Stock FDA solution tightly sealed in 4 mL or smaller aliquots in polypropylene tubes or glass tubes at ≤ § -20 °C , protected from light, for up to six months.

• Record the following on each aliquot tube:

## Preparation of Stock Propidium Iodide (PI)

- Prepare [M] **750 Micromolar (μM)** Propidium Iodide (PI) according to the formulation in sub-steps below.
- 2.1 <u>Calculation</u>:

PI FW = 668.4 Stock Concentration = 750  $\mu$ M Volume required = 25 mL

FW X Concentration X Volume =  $(668.4) \times (750 \times 10^{-6}) \times (25 \times 10^{-3}) = 0.0125 \text{ g PI}$ 

**i** protocols.io 4 12/18/2019

	■0.0125 g Propidium Iodide	
	PI Supplier	
	■ PI Lot #	
	Expiration Date	
	■25 ml DPBS	
	<ul> <li>PBS Supplier</li> </ul>	
	■ PBS Lot #	-
	Expiration Date	_
2.3	<ul> <li>Store Stock PI solution tightly sealed in 0.5 mL or sm tubes at &amp; 2 °C to &amp; 8 °C , protected from lig</li> <li>Record the following on each aliquot tube:         <ul> <li>PI (750 µM) Stock Solution Assigned Lot number:_ Expiration Date:</li></ul></li></ul>	
	This is an open access protocol distributed under the tened use, distribution, and reproduction in any medium, provided	ns of the Creative Commons Attribution License, which permits the original author and source are credited

Dissolve 0.0125 g of PI in 25 mL of DPBS and cover with aluminum foil.

2.2

**i** protocols.io 5 12/18/2019