

In vitro transcription of guide RNAs

Mark Dewitt & Julia Wong

Abstract

In vitro T7 template synthesis and transcription
sgRNA protocol with SPRI beads

Citation: Mark Dewitt & Julia Wong In vitro transcription of guide RNAs. **protocols.io**

dx.doi.org/10.17504/protocols.io.dm749m

Published: 13 Aug 2015

Guidelines

Assembly Oligos:

T7FwdVar:

GGATCCTAATACGACTCACTATAG---guide-sequence---GTTTTAGAGCTAGAA


T7RevLong:

AAAAAAGCACCGACTCGGTGCCACTTTTTCAAGTTGATAACGGACTAGCCTATTTTAACTTGCTATTTCTAGCTC
TAAAC

Amplification Primers:

T7FwdAmp: GGATCCTAATACGACTCACTATAG

T7RevAmp: AAAAAAGCACCGACTCGG

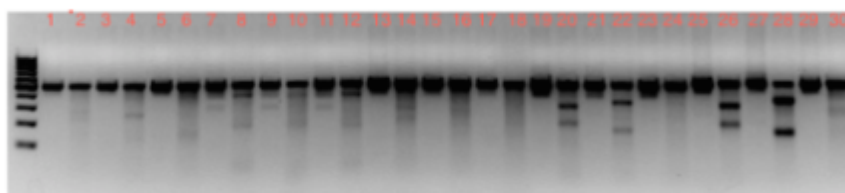

Odd numbers: just PCR product

Even numbers: with T7E1 treatment

Here's a gel that shows some good guides and some that didn't work| examples of more promising looking guides would be #'s 20, 22, 26, and 28.

Protocol

Making the Template (for 1 sgRNA)

Step 1.

see guidelines for the primers

Set up the following 50 µl reaction

■ AMOUNT

50 µl Additional info:

. T7 PCR Template for 1sgRNA Mixture

CONTACT: [Jacob Corn](#)

Step 1.1.

35.5 µl H₂O

Step 1.2.

10 µl 5x Phusion HF Buffer

Step 1.3.

1 µl 10 mM dNTPs

Step 1.4.

1 µl T7FwdVar (1 µM) ← THIS IS THE PART WE DESIGN

Step 1.5.

1 µl T7RevLong (1 µM)

Step 1.6.

0.5 µl T7FwdAmp (100 µM)

Step 1.7.

0.5 µl T7RevAmp (100 µM)

Step 1.8.

0.5 µl Phusion HF DNA polymerase (2u/µl)

PCR

Step 2.

Run PCR:

95° 30 sec

95° 10 sec

57° 10 sec

72° 10 sec

30x steps 2-4

72° 2 min

4° hold

No PCR cleanup necessary at this point

Making template

Step 3.

Mix:	Final conc.
2 µl 10x Buffer	1x
2 µl ATP (100 mM)	10 mM
2 µl GTP (100 mM)	10 mM
2 µl CTP (100 mM)	10 mM
2 µl UTP (100 mM)	10 mM
8 µl DNA template (85 ng/µl) ←PCR PDTS 25 ng/µ	
12 µl T7 RNA polymerase mix	
20 µl total	

Incubated this mix for 18 hours at 37° in a thermalcycler. Added 1 µl of RNase-free DNase; incubated

20 min, room T.

 **DURATION**

00:18:00

Run PCR

Step 4.

Run PCR

95C 30s

95C 10s

57C 10s

72C 10s

steps 2-4, 30 cycles

72 2min

4C hold

No PCR cleanup necessary at this point

T7 transcription

Step 5.

Make 20 µl total T7 transcription mix

 **AMOUNT**

20 µl Additional info:

 **PROTOCOL**

. [T7 transcription mix](#)

CONTACT: [Jacob Corn](#)

T7 transcription

Step 6.

Incubate transcription mix for 18 hours at 37° in a thermocycler

 **DURATION**

18:00:00

T7 transcription

Step 7.

Add 1 µl of RNase-free DNase; incubated 20 min, room T

 **DURATION**

00:20:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 8.

IVT sgRNAs to 20uL (or 10uL)

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 9.

Bring volume to 150 uL with 100% EtOH (this helps binding of small fragments)

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 10.

Add 5X SPRI (we use homemade SeraPure beads for RNA binding)

5*10 (IVT sgRNA)= 50 uL of SPRI Beads

5*20 (IVT sgRNA)= 100 uL SPRI Beads

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 11.

Pipette to mix 10 times

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 12.

Incubate 5 minutes at room temperature

 DURATION

00:05:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 13.

Place on magnetic stand, 5 min

 DURATION

00:05:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 14.

Discard supernatant

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 15.

Wash#1 Add 200 uL, 80% EtOH. Wait 2 min. Remove EtOH.

 DURATION

00:02:00

Step 16.

8b)Wash #2: Add 200 uL, 80% EtOH. Wait 2 min. Remove EtOH.

 DURATION

00:02:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 17.

Air dry 5-10 min (pellet will change from a glossy/wet to matte/dry.)

 DURATION

00:10:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 18.

Elute 20 uL of water or TE. Pipette mix 10 times.

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 19.

Incubate 2 minutes at room temperature

 DURATION

00:02:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 20.

Place on magnetic stand, 5 min

 DURATION

00:05:00

SPRI Beads clean-up of sgRNAs (96-well plate):

Step 21.

Keep Supernatant. Transfer to a new plate / tubes.