



May 08,
2019

In devel.

iGEM Calibration Protocol - Flow Cytometry Cell Size

Jacob Beal¹, Cheryl Telmer¹, Richard Tennant¹, Paul Rutten¹

¹iGEM Measurement Committee

dx.doi.org/10.17504/protocols.io.2pegdje

Richard Tennant

ABSTRACT

Conversion from forward scatter to E μ m is not a linear function, so data cannot be converted simply by multiplying with a scaling factor. We thus recommend use of software tools for processing data with size calibration.

This protocol can be combined with bead-based fluorescence calibration.

MATERIALS

NAME	CATALOG #	VENDOR	CAS NUMBER	RRID
Size calibration beads	PPS-6K	SpheroTech		

MATERIALS TEXT


- Experimental samples: transformed per your desired experimental procedure
- Flow cytometer with a forward scatter channel (typically a 488nm excitation and a 488nm/10nm emission filter). Measurements on this channel will be converted to equivalent μ m diameter (E μ m).
- [TASBE Flow Analytics](#) - Analysis software

BEFORE STARTING

Read through this entire protocol carefully before you start your experiment and prepare any materials you may need. See the "Results" section for an example of a completed data analysis spreadsheet.


- 1 Prepare and culture experimental samples according to your desired experimental procedure.
- 2 Prepare experimental samples as needed for running through your flow cytometer.
- 3 Follow SpheroTech directions for preparation for PPS-6K bead sample.
- 4 Measure all samples in flow cytometer
 - Using the bead sample, adjust FSC channel voltage so that the top bead peak is below detector saturation.
 - Using the negative process control sample, adjust forward-scatter and side-scatter voltages to place the strong cell peak as close to the center of the detector range as possible (see example illustrated below).

Center FSC and SSC

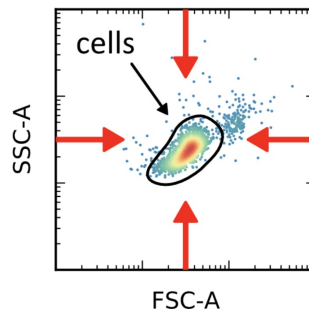

- Instrument gating should be set to ensure that no cell events are discarded (see example illustrated below).

Threshold Debris

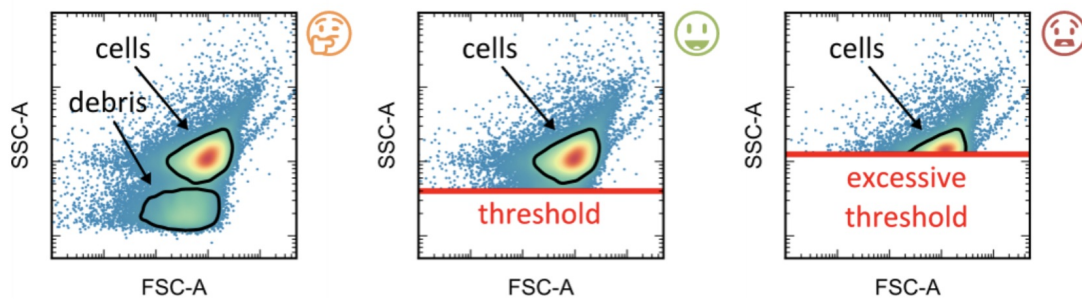

- Collect at least 10,000 events per sample.

5 Compute the $E_{\mu m/a.u.}$ scaling function from bead sample using [TASBE Flow Analytics](#) to create a TASBE color model.

6 Convert experimental data to $E_{\mu m}$ using [TASBE Flow Analytics](#) analysis functions.



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited