



Oct 02, 2019

## MojoSort™ Streptavidin Nanobeads Protocol - Positive Selection [↗](#)

Sam Li<sup>1</sup><sup>1</sup>BioLegend

1

Works for me

[dx.doi.org/10.17504/protocols.io.7wghpbw](https://doi.org/10.17504/protocols.io.7wghpbw)

BioLegend

Sam Li  
BioLegend

### ABSTRACT

#### Product description and procedure summary:

**If your target cells are the labeled cells (the positive fraction), use the Streptavidin Nanobeads Protocol – Positive Selection. If your target cells are the unlabeled cells (negative fraction), use the Streptavidin Nanobeads Protocol - Negative Selection.**

Target cells are isolated by incubating your sample with a Biotin-conjugated antibody or antibody cocktail, followed by incubation with magnetic Streptavidin Nanobeads (Cat. No. 480015/480016). The magnetically labeled fraction is retained by the use of a magnetic separator. The untouched cells are collected. These are the cells of interest; do not discard the liquid. Some of the downstream applications include functional assays, gene expression, phenotypic characterization, etc.

**Note:** This procedure is optimized for the isolation of  $10^7$  to  $2 \times 10^8$  cells per tube. If working with fewer than  $10^7$  cells, keep volumes as indicated for  $10^7$  cells. For best results, optimize the conditions to your specific cell number and tissue. Prepare fresh MojoSort™ Buffer solution by diluting the 5X concentrate with sterile distilled water.

*Scale up volumes if using 14 mL tubes and Magnet, and place the tube in the magnet for 10 minutes.*

### EXTERNAL LINK

<https://www.biolegend.com/protocols/mojosort-streptavidin-nanobeads-protocol-positive-selection/4748/>

### GUIDELINES

MojoSort™ magnetic particles can be used with other commercially available magnetic separators, both free standing magnets and column-based systems. Because MojoSort™ protocols are optimized for the MojoSort™ separator, the protocols may need to be adjusted for other systems. Please contact BioLegend Technical Service ([tech@biolegend.com](mailto:tech@biolegend.com)) for more information and guidance. We do not recommend using MojoSort™ particles for BD's IMag™ or Life Technologies' DynaMag™.

**Application notes:** To use this product in magnetic separation columns, a titration of the Nanobeads should be performed. Optimal concentration for magnetic separation columns is lot-specific. Please contact BioLegend Technical Service ([tech@biolegend.com](mailto:tech@biolegend.com)) for further assistance on how to use MojoSort™ Nanobeads in magnetic separation columns.

## General Tips & FAQ

- How much of Biotinylated antibodies should I use?
- For the biotinylated antibodies, we recommend to do a titration to determine the optimal concentration. In general, if using only a single biotin-conjugated antibody (and not a cocktail), then using the flow cytometry recommended concentration may work okay as long as you are washing after antibody incubation. If you are using a cocktail or plan to remove the wash after antibody incubation, you should perform titrations and optimizations.
- 
- How much of Streptavidin Nanobeads should I use?
- The amount of beads to use always depends on the frequency of the target, among a few other factors. The standard volume to use is 10  $\mu$ L of beads per  $10^7$  cells. You can benefit from titrating the beads, using 2, 5, 10, and 20  $\mu$ L of beads (for example). In some cases less than 10  $\mu$ L of beads is enough to achieve the purity and you can save on the reagent. Please note that Streptavidin Nanobeads volume should not exceed more than 20% of the total volume.
-

**Note:** The Antibody/Nanobead volume to add should not exceed more than 20% of the 100µL cell suspension volume. Thus, for 100µL of cell suspension do not add more than 20µL of Antibody/Nanobeads. If you need to add more than 20µL, resuspend the cells in step 3 or step 6 at a higher concentration. For example, to add 50µL of Antibody/Nanobeads, resuspend the cells to a final concentration of  $2 \times 10^8$  cells/mL. You can then add 50µL of Antibody/Nanobeads to 50µL of cells. Always keep the total volume around 100µL.

#### MATERIALS


NAME ▾	CATALOG # ▾	VENDOR ▾
MojoSort™ Buffer	480017	BioLegend
MojoSort™ Magnet	480019	BioLegend
MojoSort™ Streptavidin Nanobeads	480015, 480016	BioLegend

#### MATERIALS TEXT

- Adjustable pipettes
- 70µm filters (one per sample)
- 5mL (12 x 75mm) or 14mL (17 x 100mm) polypropylene tubes
- Reagents for sample preparation
- Reagents and instruments (Flow cytometer) to determine yield and purity

- 1 Prepare cells from your tissue of interest or blood without lysing erythrocytes.
- 2 In the final wash of your sample preparation, resuspend the cells in MojoSort™ Buffer by adding up to 4 mL in a 5 mL (12 x 75 mm) polypropylene tube.  
**Note:** Keep MojoSort™ Buffer on ice throughout the procedure.
- 3 Filter the cells with a 70µm cell strainer, centrifuge at 300xg for 5 minutes, and resuspend in an appropriate volume of MojoSort™ Buffer. Count and adjust the cell concentration to  $1 \times 10^8$  cells/mL.
- 4 Aliquot 100µL of cell suspension ( $10^7$  cells) into a new tube. Check the recommended usage for flow cytometric staining of the Biotin-conjugated antibody indicated in the antibody technical datasheet. Calculate the volume to stain  $10^7$  cells (or desired amount of cells). **Add the appropriate volume (see General Tips and FAQ below) of Biotin-conjugated antibody** to the cell suspension, mix well and **incubate on ice for 15 minutes**.  
*Optional: Take an aliquot before adding the cocktail to monitor purity and yield.*
- 5 Wash the cells by adding MojoSort™ Buffer up to 4mL. Centrifuge the cells at 300xg for 5 minutes.
- 6 Discard supernatant and resuspend in 100µL of MojoSort™ Buffer.
- 7 Resuspend the beads by vortexing, maximum speed, 5 touches. **Add the appropriate volume (see General Tips and FAQ below) of Streptavidin Nanobeads**. Mix well and **incubate on ice for 15 minutes**. Scale up the volume accordingly if separating more cells. For example, if the volume of Nanobeads for  $1 \times 10^7$  cells is 10µL, add 100µL for  $1 \times 10^8$  cells. When working with less than  $10^7$  cells, use indicated volumes for  $10^7$  cells.
- 8 Wash the cells by adding MojoSort™ Buffer up to 4mL. Centrifuge the cells at 300xg for 5 minutes.

- 9 Discard the supernatant.
- 10 Add 2.5mL of MojoSort™ Buffer.  
**Note:** If you observe aggregates, filter the suspension. To maximize yield, you can disrupt the aggregates by pipetting the solution up and down.
- 11 Place the tube in the magnet for 5 minutes.  
*Optional: Take a small aliquot before placing the tube in the magnet to monitor purity and yield. Keep unused cells to be used as control or other applications if needed.*
- 12 Pour out the unlabeled fraction. If these are your cells of interest, **DO NOT DISCARD**. Resuspend the labeled cells in 2.5mL MojoSort™ Buffer.
- 13 Repeat steps 10-12 on the labeled fraction twice more for a total of **3 separations**. Pool the unlabeled fractions and keep the labeled cells. The fraction that is not of interest may be useful as staining controls, to monitor purity/yield, or other purposes.  
*Optional: Take a small aliquot to monitor purity and yield.*

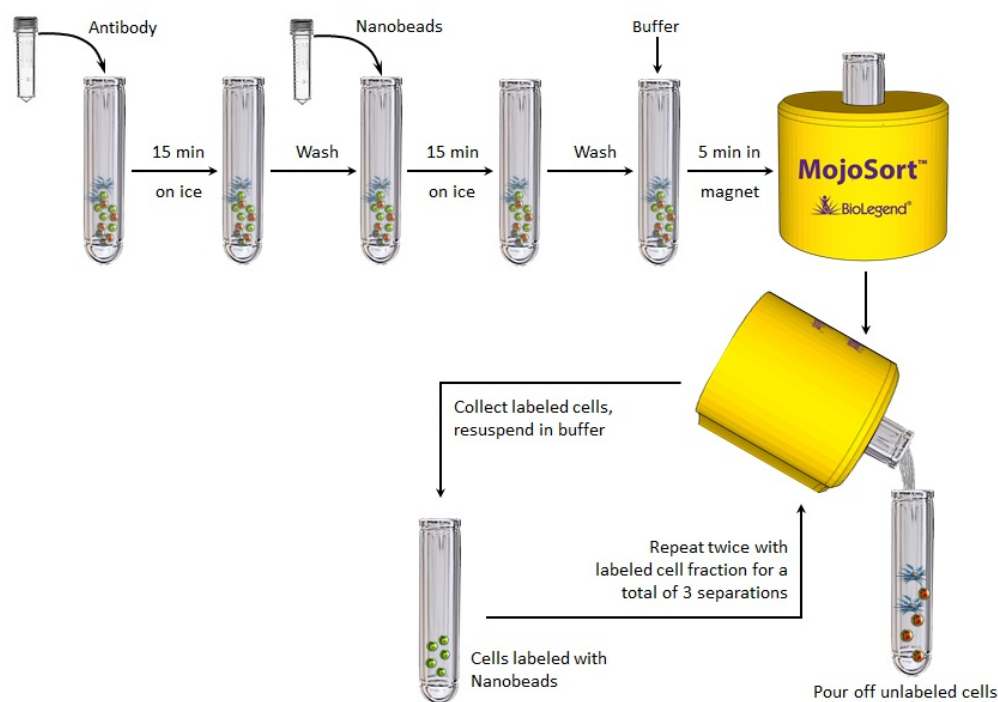

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited