

May 09, 2019

Working

U Cinn - Phospholipids Assay 👄

Patrick Tso1, Dana Lee1

¹University of Cincinnati

dx.doi.org/10.17504/protocols.io.xmvfk66

Mouse Metabolic Phenotyping Centers

Tech. support email: info@mmpc.org



Lili Liang 69

ABSTRACT

Summary:

Determinations of phospholipids in plasma/serum/lymph will be made using the Wako Phospholipids C enzymatic assay. In this assay, phospholipids in the sample are hydrolyzed ultimately producing a blue pigment. The amount of phospholipids in the sample is determined by measuring the absorbance of the blue color.

EXTERNAL LINK

https://mmpc.org/shared/document.aspx?id=198&docType=Protocol

MATERIALS

NAME CATALOG # **VENDOR CAS NUMBER** RRID FUJIFILM Wako Pure 433-36201 Phospholipid C assay kits standard included in kit **Chemical Corporation**

MATERIALS TEXT

Working Reagent:

Reagents and Materials:

Color Reagent

Buffer


Procedure:

Reconstitute one vial of Color Reagent with a portion of Buffer then transferring entire contents to Buffer bottle, rinsing Color vial several times.

Note:

FUJIFILM Wako RRID:SCR_013651

Prepare working standards by making a serial dilution of the stock 300mg/dl standard. (See diagram below) *Stock standard included in kit.

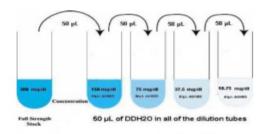

Prepare Working Reagent by reconstituting one vial of Color Reagent with a portion of Buffer then transferring entire contents to Buffer bottle, rinsing Color Reagent vial several times.
Using a 96 well flat bottom plate, into separate wells, pipette 2μL of deionized water, standard, or sample to be assayed.
Add 300μL of Working Reagent to all wells.
Incubate plate for 5 minutes at 37°C.
Determine the absorbance (abs) of the standards and of each unknown at 600nm.
Calculate values of unknowns from the standard curve.
Specimen: Serum or Plasma. Specimen stable for 7 days at 2-8°C or 3 months at -20°C.
Assay Linearity: 1000 mg/dl

This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Reagent Stability: 7 days at 2-8°C Stability of Final Reaction: 60 minutes