



Jul 21, 2019

Calibration Protocol - Particle Standard Curve with Microspheres

Forked from [Calibration Protocol - Particle Standard Curve with Microspheres](#)

Jacob Beal¹, Traci Haddock-Angelli², Markus Gershater³, Vishal Sanchania³, Russell Buckley-Taylor¹, Geoff Baldwin⁴, Natalie Farny¹, [Richard Tennant](#)¹, [Paul Rutten](#)¹

¹iGEM Measurement Committee, ²iGEM, ³Synthace, ⁴Imperial College London

1 Works for me dx.doi.org/10.17504/protocols.io.5n5g5g6

iGEM Measurement

Tech. support email: pauljrutten@gmail.com

Paul Rutten
The University of Oxford



ABSTRACT

You will prepare a dilution series of monodisperse silica microspheres and measure the Abs₆₀₀ in your plate reader.

The size and optical characteristics of these microspheres are similar to cells, and there is a known amount of particles per volume. This measurement will allow you to construct a standard curve of particle concentration which can be used to convert 600 nm absorbance measurements into an estimated equivalent number of cells.

EXTERNAL LINK

<https://2019.igem.org/Measurement>



iGEM Data Analysis
Template - Particle
Standard Curve - v1.xlsx

MATERIALS

| NAME | CATALOG # | VENDOR |
|--------------------------------|-----------|--------|
| 96 well plate | | |
| double distilled water (ddH2O) | | |
| 300µl Silica beads | | |

STEPS MATERIALS

| NAME | CATALOG # | VENDOR |
|--------------------|-----------|--------|
| 300µl Silica beads | | |
| ddH2O | | |

MATERIALS TEXT


300 µL Silica beads are provided in the iGEM Measurement Kit. The 96-well plate should preferably be black with a clear flat bottom.

BEFORE STARTING

Read through this entire protocol carefully before you start your experiment and prepare any materials you may need.

Prepare the Microsphere Stock Solution

- 1 Obtain the tube labeled "Silica Beads" from the Measurement Kit and vortex vigorously for 30 seconds.

 300µl Silica beads



Microspheres should NOT be stored at 0°C or below, as freezing affects the properties of the microspheres. If you believe your microspheres may have been frozen, please contact the iGEM Measurement Committee for a replacement (measurement@igem.org).

- 2 Immediately pipet 100 µL microspheres into a 1.5 mL eppendorf tube

- 3 Add 900 µL of ddH₂O to the microspheres

 ddH₂O


- 4 Vortex well. This is your Microsphere Stock Solution

Prepare the serial dilution of microspheres

- 5 Accurate pipetting is essential. Serial dilutions will be performed across columns 1-11. **Column 12 must contain ddH₂O only.**

Initially you will setup the plate with the microsphere stock solution in column 1 and an equal volume of 1x ddH₂O in columns 2 to 12.

You will perform a serial dilution by consecutively transferring 100 µl from column to column with good mixing.

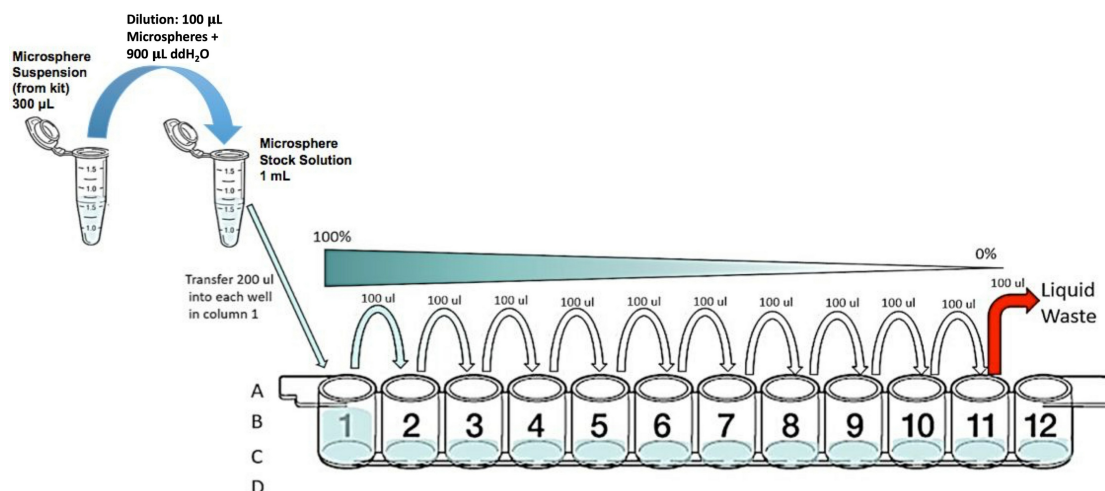

- 6 Add 100 µl of ddH₂O into wells A2, B2, C2, D2....A12, B12, C12, D12

- 7 Vortex the tube containing the stock solution of microspheres vigorously for 10 seconds
- 8 Immediately add 200 μ l of microspheres stock solution into A1
- 9 Transfer 100 μ l of microsphere stock solution from A1 into A2
- 10 Mix A2 by pipetting up and down 3x and transfer 100 μ l into A3
- 11 Mix A3 by pipetting up and down 3x and transfer 100 μ l into A4
- 12 Mix A4 by pipetting up and down 3x and transfer 100 μ l into A5
- 13 Mix A5 by pipetting up and down 3x and transfer 100 μ l into A6
- 14 Mix A6 by pipetting up and down 3x and transfer 100 μ l into A7
- 15 Mix A7 by pipetting up and down 3x and transfer 100 μ l into A8
- 16 Mix A8 by pipetting up and down 3x and transfer 100 μ l into A9
- 17 Mix A9 by pipetting up and down 3x and transfer 100 μ l into A10
- 18 Mix A10 by pipetting up and down 3x and transfer 100 μ l into A11
- 19 Mix A11 by pipetting up and down 3x and transfer 100 μ l into liquid waste



Take care not to continue serial dilution into column 12

- 20 Repeat dilution series for rows B, C, D
- 21 **IMPORTANT!**

Re-Mix (pipette up and down) each row of your plate **immediately before** putting in the plate reader! (This is important because the beads begin to settle to the bottom of the wells within about 10 minutes, which will affect the measurements.)



Take care to mix gently and avoid creating bubbles on the surface of the liquid

Measure OD

22 Measure OD₆₀₀ of all samples in instrument

23 Record the data in your notebook

24 Import data into this Excel sheet:



[iGEM Data Analysis Template - Particle Standard Curve - v1.xlsx](#)

Congratulations!

25 You have now completed this calibration protocol



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited