

Dengue virus type 2 (DENV-2) capsid-Thai TaqMan assay (no longer in use; see Guidelines)

Ian Mackay, Judy Northill, Alyssa Pyke

Abstract

This protocol was designed and developed at this laboratory and incorporates a previously published oligoprobe (see below).

The protocol specifically aims to amplify DENV02 viruses and not other dengue viruses. The assay targets the capsid region and is designed as a qualitative test for investigating suspected human cases of DENV-2 infections.

This assay has been superseded by the Dengue virus type 2 (DENV-2) MGB TaqMan (DENV2-2016MGB) assay.

Citation: Ian Mackay, Judy Northill, Alyssa Pyke Dengue virus type 2 (DENV-2) capsid-Thai TaqMan assay (no longer in use; see Guidelines). **protocols.io**

dx.doi.org/10.17504/protocols.io.q4ydyxw

Published: 10 Jul 2018

Guidelines

This was a past assay that we no longer in use.

For the best DENV-2 TaqMan assay, please refer to our recommended protocol:

Dengue virus type 2 (DENV-2) MGB TaqMan (DENV2-2016MGB) assay

<https://www.protocols.io/view/dengue-virus-type-2-denv-2-mgb-taqman-denv2-2016mg-n7kdhkw>

Before start

- If using a different brand or model of real-time thermocycler, check the concentration of ROX is adequate.
- Method assumes the user is familiar with the thermocycler and software used to run the protocol and with PCR in general.

Materials

 SuperScript™ III Platinum™ One-Step qRT-PCR Kit 11732088 by Life Technologies

Protocol

Oligonucleotide sequences

Step 1.

Name	Sequence 5'-3'
------	----------------

D2TaqC(b)-f ¹	TTCATGGCCCT K GTGGC
D2Cor05r ¹	CCCCATCT Y TTY A RTATCCCTG
D2TaqCor-FAM ²	FAM - TCCTTCGTTTCCTAACAATCC- TAMRA


¹ Designed by Alyssa Pyke, Public Health Virology, 2005; ² Previously published, Callahan *et al.*
<http://jcm.asm.org/content/39/11/4119.abstract>

Reagents

Step 2.



REAGENTS

 SuperScript™ III Platinum™ One-Step qRT-PCR Kit 11732088 by Life Technologies

Reaction set-up

Step 3.

The assay has been used on both a Rotor-Gene 6000 and a Rotor-Gene Q real-time thermocycler

Prepare sufficient mix for the number of reactions.

Include a suitable 'dead volume' as necessary if using a robotic dispenser.

MIX PREPARATION

Reagent	Volume (µl) x1	Final reaction concentration
Nuclease-free water	4.44	N/A
D2TaqC(b)-f 150pmol/µl	0.04	300nM
D2Cor05r 150pmol/µl	0.04	300nM
D2TaqCor-FAM 75pmol/µl	0.04	150nM
2X Reaction Mix ¹	10	1X
SuperScript® III/Platinum® Taq Mix ¹	0.4	1X
ROX Reference Dye (25µM)	0.04	0.05µM
Template	5	N/A
TOTAL	20	

¹Superscript™III Platinum™ One-step qRT-PCR kit

- Dispense 15µL to each reaction well.
- Add 5µL of template (extracted RNA, controls or NTC [nuclease-free water]).
- Total reaction volume is 20µL

Amplification

Step 4.

CYCLING CONDITIONS

50°C	5min	1X
95°C	2min	1X
95°C	3sec	40X
60°C	30sec ¹	

¹Florescence acquisition step

Result Analysis

Step 5.

The definition used for a satisfactory positive result from a real-time fluorogenic PCR should include each of the following:

1. A **sigmoidal curve** – the trace travels horizontally, curves upward, continues in an exponential rise and followed by a curve towards a horizontal plateau phase
2. A **suitable level of fluorescence** intensity as measured in comparison to a positive control (y-axis)
3. A **defined threshold (C_T) value** which the fluorescent curve has clearly exceeded (Fig.1 arrow), which sits early in the log-linear phase and is <40 cycles
4. A flat or non-sigmoidal curve or a curve that crosses the threshold with a $C_T > 40$ cycles is considered a negative result.
5. NTCs should not produce a curve

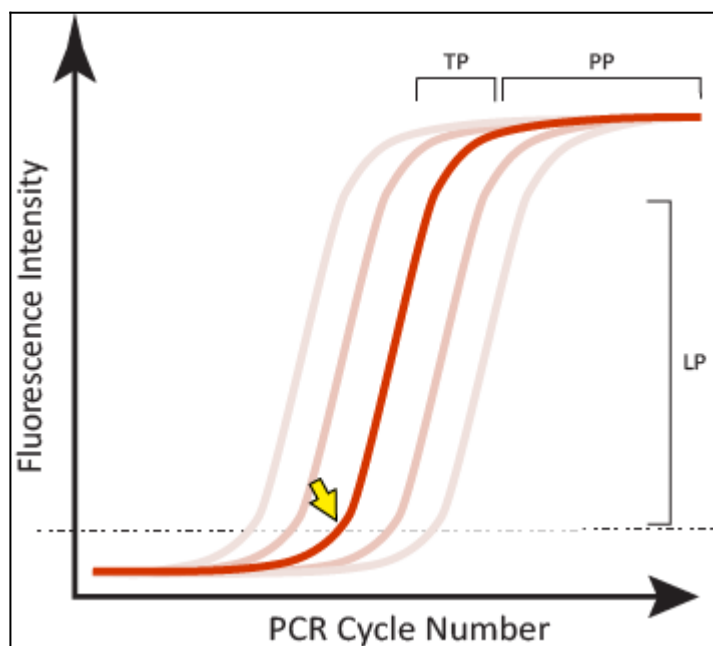


Figure 1. Examples of satisfactory sigmoidal amplification curve shape when considering an assay's fluorescent signal output. The crossing point or threshold cycle (C_T) is indicated (yellow arrow); it is the value at which fluorescence levels surpass a predefined (usually set during validation, or arbitrary) threshold level as shown in this normalized linear scale depiction. LP-log-linear phase of signal generated during the exponential part of the PCR amplification; TP-a slowing of the amplification and accompanying fluorescence signal marks the transition phase; PP-the plateau phase is reached when there is little or no increase in fluorescent signal despite continued cycling.

