

## Enterovirus (EV) A71 TaqMan 2018 (EV-A71-TM2018)

Version 1

Ian Mackay<sup>1</sup>, Judy Northill<sup>1</sup>

<sup>1</sup>Public Health Virology, Forensic and Scientific Services

[dx.doi.org/10.17504/protocols.io.w6ffhbn](https://doi.org/10.17504/protocols.io.w6ffhbn)

Public Health Virology, Forensic and Scientific Services



Ian Mackay

Public Health Virology, Forens...



### ABSTRACT

This protocol aims to amplify enterovirus (EV) A71 viruses but not other viruses.

This protocol was designed by us.

The oligonucleotides target the 5'UTR noncoding region. This is a qualitative test for investigating EV-A71 infection of humans.

The test has identified both historical EV-A71 strains and contemporary strains and has confirmed these using subgenomic sequencing of [partial 5'UTR-VP2](#) and partial VP1 sequencing. Culture of the virus is not required as this assay is capable of detecting EV-A71 direct from extracted clinical samples.

### PROTOCOL STATUS

#### Working

We use this protocol in our group and it is working

### STEPS MATERIALS

#### NAME

SuperScript™ III Platinum™ One-Step qRT-PCR Kit

#### CATALOG #

11732088

#### VENDOR

Life Technologies

### BEFORE STARTING

If using a different brand or model of real-time thermocycler, check the concentration of ROX is adequate.

Method assumes the user is familiar with the thermocycler and software used to run the protocol and with PCR in general.

### Oligonucleotide sequences

1	Name	Sequence 5'-3'
	EVA71-VP4-For1	TAYTAYAAAGAYTCBTATGCGY
	EVA71-VP4-Rev1	CCTTRACAGGRTTWCRAACTT
	EVA71-VP4-Rev2	CTTTRACAGGRTTWCAGAAATTT
	EVA71-VP4-Rev3	CCTTCACAGGGTTCGAAACTT
	EVA71-VP4-P1	FAM - ACAGCVGGCAAGCAGAGYCTCAA - BHQ1
	EVA71-VP4-P2	FAM - ACAGCRGGYAAACAGAGYCTCAA - BHQ1
	EVA71-VP4-P3	FAM - ACTGCTGGAAAGCAAAGTCTCAA - BHQ1

### The design philosophies.

- to ensure our assay detects old and new variants of EV-D68
- to reduce the total number of degenerate positions in any single primer

## Reagents

2



SuperScript™ III Platinum™ One-Step qRT-PCR Kit  
by Life Technologies  
Catalog #: 11732088

## Reaction set-up

3

The assay has been used on both a Rotor-Gene 6000 real-time thermocycler

- Prepare sufficient mix for the number of reactions.
- Include a suitable 'dead volume' as necessary if using a robotic dispenser.

Reagent	Volume (μl) x1	Final reaction concentration
Nuclease-free water	2.91	N/A
EVA71-VP4-For1 200pmol/μl	0.09	900nM
EVA71-VP4-Rev1 200pmol/μl	0.09	900nM
EVA71-VP4-Rev2 200pmol/μl	0.09	900nM
EVA71-VP4-Rev2 200pmol/μl	0.09	
EVA71-VP4-P1 100pmol/μl	0.03	150nM
EVA71-VP4-P2 100pmol/μl	0.03	150nM
EVA71-VP4-P3 100pmol/μl	0.03	150nM
MgSO4 (50mM)	1.2	6mM
2X Reaction Mix <sup>1</sup>	10	1X
SuperScript® III/Platinum® Taq Mix <sup>1</sup>	0.4	1X
ROX Reference Dye (25μM)	0.04	50nM
Template	5	N/A
<b>TOTAL</b>	20	

1-Superscript™ III Platinum™ One-step qRT-PCR kit

- Dispense 15μL to each reaction vessel.
- Add 5μL of template (extracted RNA, controls or no-template control [NTC; nuclease-free water]).
- Total reaction volume is 20μL

## Amplification

4

50°C	5min	1X
95°C	2min	1X
95°C	3sec	140X
60°C	30sec <sup>1</sup>	1


1-Fluorescence acquisition step

## Result Analysis

5

The definition used for a satisfactory positive result from a real-time fluorogenic PCR should include each of the following:

- A **sigmoidal curve** – the trace travels horizontally, curves upward, continues in an exponential rise and followed by a curve towards a horizontal plateau phase
- A **suitable level of fluorescence** intensity as measured in comparison to a positive control (y-axis)
- A **defined threshold ( $C_T$ ) value** which the fluorescent curve has clearly exceeded (Fig.1 arrow), which sits early in the log-linear phase and is <40 cycles
- A flat or non-sigmoidal curve or a curve that crosses the threshold with a  $C_T > 40$  cycles is considered a negative result.
- NTCs should not produce a curve

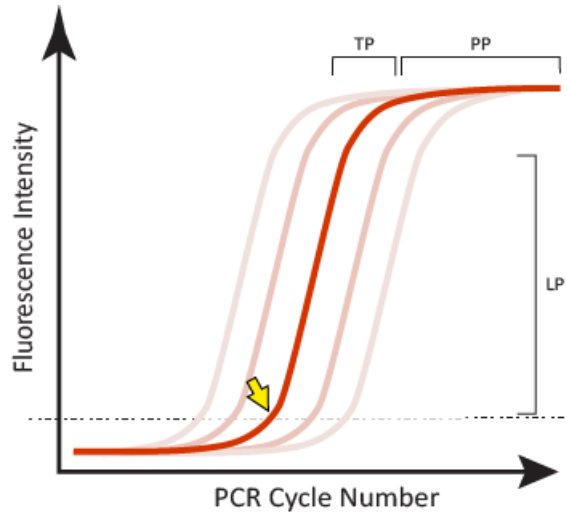

**Figure 1.** Examples of satisfactory sigmoidal amplification curve shape when considering an assay's fluorescent signal output. The crossing point or threshold cycle ( $C_T$ ) is indicated (yellow arrow); it is the value at which fluorescence levels surpass a predefined (usually set during validation, or arbitrary) threshold level as shown in this normalized linear scale depiction. LP-log-linear phase of signal generated during the exponential part of the PCR amplification; TP-a slowing of the amplification and accompanying fluorescence signal marks the transition phase; PP-the plateau phase is reached when there is little or no increase in fluorescent signal despite continued cycling.



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited