

Diatom Transformation via Bacterial Conjugation

Jeric Harper

Abstract

Citation: Jeric Harper Diatom Transformation via Bacterial Conjugation. **protocols.io**

dx.doi.org/10.17504/protocols.io.fp7bmrn

Published: 03 Oct 2016

Before start

- Make sure you are using fresh *E. coli* cells streaked for isolation on LB + antibiotics no more than 1 week from -80°C cryostock.

Protocol

Step 1.

Diatom cultivation (liquid)

Grow the diatom culture to mid-log phase ($\approx 8.0E^6$ cells/ml for *Phaeodactylum tricornutum* grown on F/2 media). We have investigated transformation efficiency throughout the growth curve and found this to be the sweetspot.

NOTES

Jeric Harper 30 Aug 2016

- [The original protocol](#) plated the culture (250uL of $1.0E8$ cells/ml) on 1/2 strength L1, 1% agar plates for 4 days prior to transformation.

- We have been investigating liquid cultivation because many diatom species do not survive on agar.

Jeric Harper 30 Aug 2016

Cultivation Conditions

- Cultivation in 0.2 um-filtered L1 media prepared using seawater collected from the Gulf of Mexico.

- 80 - 100 $\mu\text{mol photons m}^{-2} \text{ s}^{-1}$ provided by Philips Daylightt Delux 40 watt T12 flourscent tubes.
- Semi-continuous operation in 1L bulbble column bioreactors (500 - 800 ml working volume).
- 0.2 μm filtered air containing 1% CO_2 provides aeration/agitation at 0.2 vvm.

Step 2.

Prepare *E. coli* starter culture

Grow 1 mL of *E. coli* culture overnight (16-20 hrs) in LB+antibiotic for each planned transformation. We grow them at 37°C at 270 rpm in a shaking incubator.

Step 3.

Outgrow *E. coli*

Dilute cells 1:50 into 50 mL of fresh LB+antibiotic for each planned transformation.

On the day of transformation, use the overnight culture to inoculate 50 mL of fresh LB+antibiotic for each planned transformation.

- Grow to an OD_{600} of 0.8 - 1.0 (37°C with 270 rpm shaking).
- This takes about 3-4 hours.

Step 4.

Concentrate the diatom and *E. coli* cultures

Centrifuge both cultures at 3000 x g for 10 minutes.

- The diatom and *E. coli* cultures should be centrifuged at around the same time to minimize the amount of time they spend concentrated.

E. coli:

- Centrifuge 50 mL of culture for each transformation and resuspend in 500 μL of SOC medium.
-

Phaeodactylum tricornutum:

- Measure cell concentration with FlowCam or haemocytometer.
- Centrifuge 2.5×10^8 cells and resuspend in 500 μL of L1 medium.

🔗 NOTES

Jeric Harper 30 Aug 2016

[The original protocol](#) scraped the agar plates the culture was initiated on using 500 μL F/2, then adjusted the volume to attain 5.0×10^8 cells/ml. We have found no difference so far in transformation

efficiency between liquid- and plate-initiated cultures.

Step 5.

Conjugation

- **Mix** 200 µl of *E. coli* cells with 200 µl of *Phaeodactylum tricornutum* cells.
- **Spread** the mixture on transformation plates (0.5x BG-11 with 5% LB and 1% agar).
- **Incubate** plates for 90 minutes at 30°C in the dark.
- **Move** plates to light incubator (18°C and 80 µmol photons m⁻² s⁻¹) for 2 days.

Negative controls

- Mix 200 µl of SOC medium with 200 µl of *Phaeodactylum tricornutum* cells.
- Mix 200 µl of L1 medium with 200 µl of *E. coli* cells.
- Incubate and treat plates identically to conjugation plates.

■ ANNOTATIONS


Jeric Harper 30 Aug 2016

- The authors of the original protocol suggested spreading the mixture near, but not touching, the edges of the plate, as this mixture will be scraped and replated in 2 days.
- According to [Diner et. al](#), the conjugation occurs during the 2-day incubation.
- We have found greater success in using 0.5x BG-11 medium made with seawater instead of 0.5x L1 medium.

Step 6.

Selection

- **Collect** cells by adding 1 mL of L1 medium. Use a cell scraper to concentrate cells and medium to one side of the plate. Transfer resuspended cells to a 1.5 ml microcentrifuge tube with a P1000 pipette and filter tips.
- **Plate** 200 µl of the cell suspension on selection plates.
- Incubate at 18°C and 80 µmol photons m⁻² s⁻¹ until colonies appear.

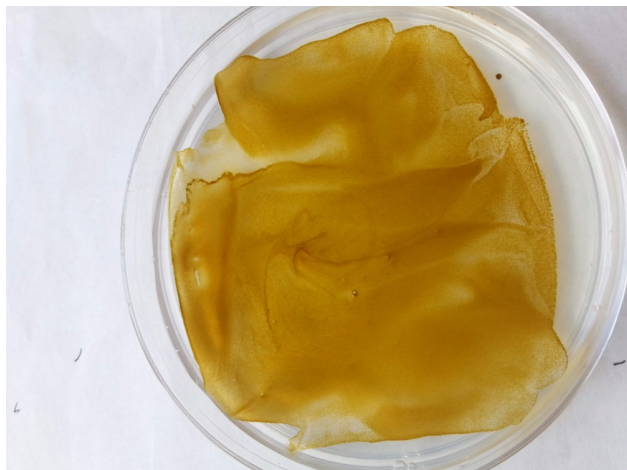

📌 NOTES


Jeric Harper 03 Oct 2016

It is important to use filter tips as small pieces of agar that are accidentally scraped up can cause the pipette to cavitate resulting in contamination.

Step 7.

Colony identification

After a minimum of 8-12 days, untransformed *Phaeodactylum tricornutum* cells die off, and colonies of transformed cells begin to appear – in some cases, this can take 3-4 weeks.

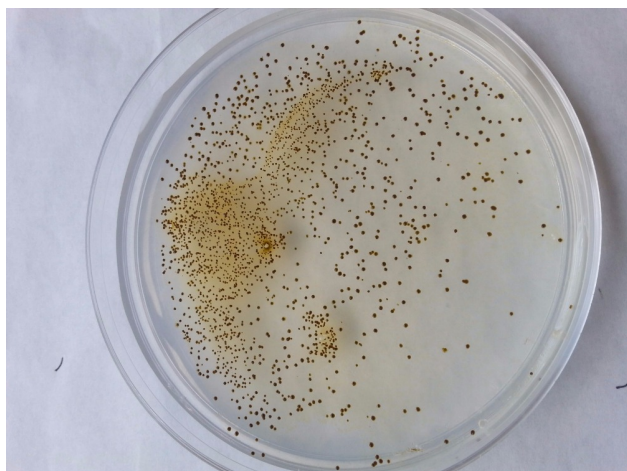

Step 8.

This protocol was modified from [the original procedure](#) and correspondence with the authors.