



Dec 03, 2019


## Cloning of standardized L0 parts into pUAP4 for Loop type IIS assembly

Forked from [Cloning of L0 parts into pUAP-ye for Loop type IIS](#)Eftychis Frangedakis<sup>1</sup>, Susana Sauret-Gueto<sup>1</sup><sup>1</sup>Plant Sciences, University of Cambridge, OpenPlant1 *Works for me* [dx.doi.org/10.17504/protocols.io.93dh8i6](https://doi.org/10.17504/protocols.io.93dh8i6)[OpenPlant Project](#)

Susana Sauret-Gueto

Plant Sciences, University of Cambridge, OpenPlant ⚡ 🌞 🌱

### Summary of primers and cloning into pUAP4

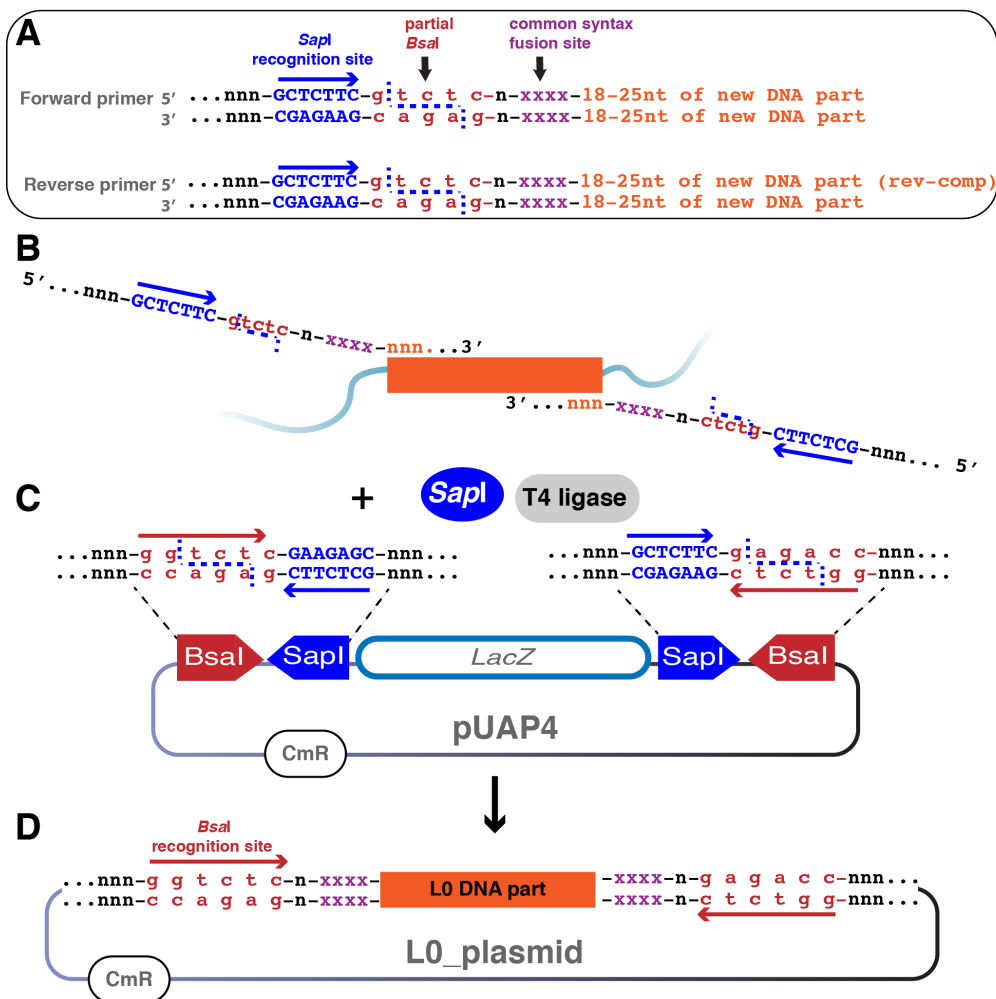

Design of primers for cloning of L0 parts into pUAP4 vectors. (A) To clone a standardized L0 part flanked by the 5' and 3' common syntax fusion sites specific to a type of part, the fusion sites are added to the DNA part by PCR, using specially designed primers. Primers include, in order from the 5' end: 3 random bp (in black), the 7 bp SapI recognition site (in blue), 5 bp that correspond to a partial BsaI recognition site (in red), one random bp for spacing (in black), and then the common syntax fusion site (in purple) followed by 18-25 bp complementary to the DNA part to be amplified (in orange). (B) The DNA part should be amplified with a high fidelity DNA polymerase and (C) it is cloned into pUAP4 using SapI type IIS assembly to produce (D) the desired L0 plasmid. Once cloned into pUAP4, the full BsaI recognition site is reconstituted. BsaI type IIS cloning will be used to assemble multiple L0 parts into a transcription unit (L1). Blue arrows: SapI recognition site. Blue dashed lines: SapI cleavage site. Red arrows: BsaI recognition site. CmR: chloramphenicol bacterial resistance cassette. LacZ: lacZα cassette for blue-white screening of colonies (negative blue colonies contain undigested pUAP4, positive white colonies contain a L0 part inserted into pUAP4).

#### PCR amplification of DNA part

- Design primers as described in step 1 to add the SapI recognition site, a partial BsaI recognition site, and the common syntax fusion site.
- L0 parts are PCR amplified from the source DNA part (e.g. plasmid DNA or genomic DNA), using a high fidelity DNA polymerase such as Phusion. Use 10 ng of DNA if the template is plasmid DNA and 100ng if the template is genomic DNA. The cycling conditions are: Denaturation at 98 °C for 30 s. 35 cycles of: denaturation at 98 °C for 10 s, annealing at primer annealing temperature for 30 s, and extension at 72 °C for 15 s/kb. Final extension at 72 °C for 10 min.

- 4 Run the PCR products on a 1.5% (w/v) agarose gel.
- 5 Gel extract the band that corresponds to the size of the amplified L0 part using a kit such as QIAquick Gel Extraction Kit.

#### SapI assembly

- 6 Determine L0 part and pUAP4 concentration with spectrophotometry (Nanodrop).
- 7 Prepare aliquots at a concentration of 15 nM for the L0 part and of 7.5 nM for the pUAP4 acceptor vector.
- 8 Set up a SapI Type IIS assembly reaction into a 0.2 mL tube according to Table:

Sterile water	5 µL
pUAP4	1 µL
Amplified L0 part	1 µL
10x Tango buffer (Thermo Fisher)	1 µL
1 mg/mL bovine serum albumin (NEB)	0.5 µL
T4 DNA ligase (5 U/µL) (Thermo Fisher)	0.25 µL
10mM ATP (SIGMA)	1 µL
SapI (LgI) (5 U/µL) (Thermo Fisher)	0.25 µL
Final volume	10 µL

- 9 Place samples on a thermocycler and use the following program:  
 Assembly: 26 cycles of 37 °C for 3 min and 16 °C for 4 min.  
 Termination and enzyme denaturation: 50 °C for 5 min and 80 °C for 10 min.
- 10 Transform 20 µL of chemically competent E. coli cells (transformation efficiency of  $1 \times 10^7$  transformants/µg plasmid DNA) using 2 µL of the assembly reaction and then spread on LB agar plates containing 25 µg/mL chloramphenicol and 40 µg/mL of X-gal for blue-white screening.
- 11 Incubate overnight at 37 °C.
- 12 Select white colonies for sequencing.
- 13 Confirm the presence of the correct insert with Sanger sequencing using the primers UAP\_F (CTCGAGTGGCCACCTGACGTCTAAGAAAC) and UAP\_R (CGAGGAAGCCTGCATAACGCGAAGTAATC) and any additional DNA part specific primers.



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited