

Measles virus TaqMan RT-PCR (no longer in regular use; see Guidelines)

Mitchell Finger, Michael Lyon, Judy Northill, Ian Mackay

Abstract

This real-time TaqMan-MGB RT-PCR protocol aimed to amplify measles virus (MeV) strains and not other viruses.

Michael Lyon and Mitchell Finger designed the assay in 2009 using Primer Express software.

The method was later published by Greg Smith in 2010 (see below).

The assay targets the fusion (F) gene region and is designed as a qualitative test for investigating MeV infection of humans.

This was a past assay that we no longer in use. For our favoured Measles virus TaqMan test, please refer to the *MeV N TaqMan* protocol.

Citation: Mitchell Finger, Michael Lyon, Judy Northill, Ian Mackay Measles virus TaqMan RT-PCR (no longer in regular use; see Guidelines). **protocols.io**

dx.doi.org/10.17504/protocols.io.rs6d6he

Published: 24 Aug 2018

Before start

- If using a different brand or model of real-time thermocycler, check the concentration of ROX is adequate.
- Method assumes the user is familiar with the thermocycler and software used to run the protocol and with PCR in general.

Protocol

Oligonucleotide sequences

Step 1.

Name	Sequence 5'-3'
Primer Measles MGB FP	GCTCAAATTGCTCAGATACTATACAGAAA
Primer Measles MGB RP	GCAGATATGGGGTCCCGTAA
Probe Measles MGB Probe	FAM - CCTGTCATTATTTGGCC - MGBNFQ

Reagents

Step 2.



REAGENTS

 SuperScript™ III Platinum™ One-Step qRT-PCR Kit [11732088](#) by [Life Technologies](#)

Reaction set-up

Step 3.

The assay has been used on both a Rotor-Gene 6000 and a Rotor-Gene Q real-time thermocycler

Prepare sufficient mix for the number of reactions.

Include a suitable 'dead volume' as necessary if using a robotic dispenser.

MIX PREPARATION

Reagent	Volume (µl) x1	Final reaction concentration
Nuclease-free water	4.45	N/A
Measles MGB FP 150pmol/µl	0.04	300nM
Measles MGB RP 150pmol/µl	0.04	300nM
Measles MGB Probe 100pmol/µl	0.03	155nM
2X Reaction Mix ¹	10	1X
SuperScript® III/Platinum® Taq Mix ¹	0.4	1X
ROX Reference Dye (25µM)	0.04	0.05µM
Template	5	N/A
TOTAL	20	

¹Superscript™III Platinum™ One-step qRT-PCR kit

- Dispense 15µL to each reaction well.
- Add 5µL of template (extracted RNA, controls or NTC [nuclease-free water]).
- Total reaction volume is 20µL

Amplification

Step 4.

CYCLING CONDITIONS

50°C	5min	1X
95°C	2min	1X
95°C	3sec	40X
60°C	30sec ¹	

¹Fluorescence acquisition step

Result Analysis

Step 5.

The definition used for a satisfactory positive result from a real-time fluorogenic PCR should include each of the following:

1. A **sigmoidal curve** – the trace travels horizontally, curves upward, continues in an exponential rise and followed by a curve towards a horizontal plateau phase
2. A **suitable level of fluorescence** intensity as measured in comparison to a positive control (y-axis)
3. A **defined threshold (C_T) value** which the fluorescent curve has clearly exceeded (Fig.1 arrow), which sits early in the log-linear phase and is <40 cycles
4. A flat or non-sigmoidal curve or a curve that crosses the threshold with a $C_T >40$ cycles is considered a negative result.
5. NTCs should not produce a curve

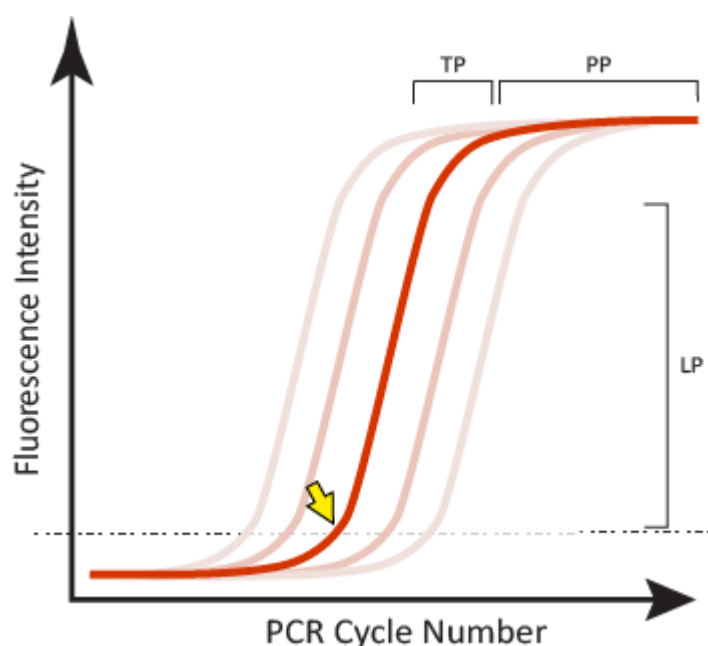


Figure 1. Examples of satisfactory sigmoidal amplification curve shape when considering an assay's fluorescent signal output. The crossing point or threshold cycle (C_T) is indicated (yellow arrow); it is the value at which fluorescence levels surpass a predefined (usually set during validation, or arbitrary) threshold level as shown in this normalized linear scale depiction. LP-log-linear phase of signal generated during the exponential part of the PCR amplification; TP-a slowing of the amplification and accompanying fluorescence signal marks the transition phase; PP-the plateau phase is reached when there is little or no increase in fluorescent signal despite continued cycling.