

# Human Parechovirus A real-time RT-PCR ["Nix assay"; 2008-2015] Version 3

Ian M Mackay

## Abstract

I and my team used this assay between 2008-2015; we dubbed it the "Nix assay". It targets the 5'UTR and employs quite degenerate oligos.

*In silico* sequence alignments indicated the oligonucleotides could theoretically detect at least HPeV 1-7, 17 and 18.

However during a period of assay comparison, another assay (see link below), the "Benschop assay" ([J.Clin.Virol. 2008. 41\(2\):69-74](#)), was found to produce more sigmoidal and higher curves and 1-5 cycle improvements to  $C_T$  values when compared among the sample sample extract set.

**Citation:** Ian M Mackay Human Parechovirus A real-time RT-PCR ["Nix assay"; 2008-2015]. **protocols.io**

[dx.doi.org/10.17504/protocols.io.krpcv5n](https://dx.doi.org/10.17504/protocols.io.krpcv5n)

**Published:** 11 Nov 2017

## Guidelines

- This protocol assumes the user is familiar with working in a laboratory, with PCR, the thermocycler and software used to run it
- This protocol should be re-evaluated if being used with different reagents, if the oligonucleotide sequences are changed or if the cycling conditions are altered

## Protocol

### Oligonucleotides...

#### Step 1.

| Name                               | 5'-3' oligonucleotide sequence       |
|------------------------------------|--------------------------------------|
| AN345_panHPeV/LV (sense primer)    | GTAACASWWGCCTCTGGGSCCAAAAG           |
| AN344_panHPeV/LV (antisense probe) | GGCCCCWGRTCAGATCCAYAGT               |
| AN257_HPeV/LV (probe)              | FAM-CCTRYGGGTACCTYCWGGGCATCCTTC-BHQ1 |

### Reagents

#### Step 2.



#### REAGENTS

SuperScript™ III Platinum™ One-Step qRT-PCR Kit [11732088](#) by [Life Technologies](#)

## Reaction setup...

### Step 3.

Below is the reaction setup for a single RT-PCR reaction.

Ideally, this work is conducted in a laboratory separate to any space used to perform PCR, molecular cloning or the analysis of high concentration DNA.

This volume has been used in 0.1-0.2ml tubes or various other connected tube configurations such as 100-place rings.

Multiply this according to the number of reactions you will need, remembering to include a positive control and at least two non-template controls (NTCs)

You may also need to allow some extra volume, depending on the method used to pipette mix into tubes for the run. For example, some robot-loaded tubes can require two reaction 'dead volumes'.

| <b>Reagent (stock concentration)</b>            | <b>Vol (μL) / reaction</b> | <b>Final concentration</b> |
|---|----------------------------|----------------------------|
| Nuclease free water                             | 4.47                       | N/A                        |
| AN345_panHPEV/LV (200pmol/ul)                   | 0.03                       | 300nM                      |
| AN344_panHPEV/LV (200pmol/ul)                   | 0.03                       | 300nM                      |
| AN257_HPEV/LV FAM-BHQ1 (100pmol/ul)             | 0.03                       | 150nM                      |
| 2X Reaction Mix <sup>1</sup>                    | 10                         | 1X                         |
| Rox Reference Dye 25mM <sup>1</sup>             | 0.04                       | 50nM                       |
| SuperScript® III/Platinum® Taq Mix <sup>1</sup> | 0.4                        | 1X                         |
| Template extract RNA                            | 5                          | N/A                        |
| Final volume                                    | 20μl                       |                            |

<sup>1</sup>SuperScript® III Platinum® One-Step qRT-PCR Kit, Cat No. 11732088

## Amplification...

### Step 4.

This assay has been optimized and validated for use with a Rotor-Gene 6000 or Rotor-Gene Q thermal cycler.

The cycling conditions are as follows:

| <b>RT-PCR</b> |      |     |
|---------------|------|-----|
| 50°C          | 5min |     |
| 95°C          | 2min |     |
| 95°C          | 3s   | 40X |
| 60°C          | 30s* |     |


\*Florescence acquisition step

Result calling...

### Step 5.

The definition used for a satisfactory positive result from a real-time fluorogenic PCR should include each of the following:

1. A **sigmoidal curve** – the trace travels horizontally, curves upward, continues in an exponential rise, curves and reaches a horizontal plateau phase
2. A **suitable level of fluorescence** intensity as measured in comparison to a positive control (y-axis)
3. A defined threshold ( $C_T$ ) value which the fluorescent curve has clearly exceeded (Fig.1 arrow) and which sits early in the log-linear phase
4. A flat or non-sigmoidal curve or a curve that crosses the threshold with a  $C_T$  value  $>40$  cycles is considered a negative result
5. No template controls (NTCs; water instead of specimen extract) should not produce a curve

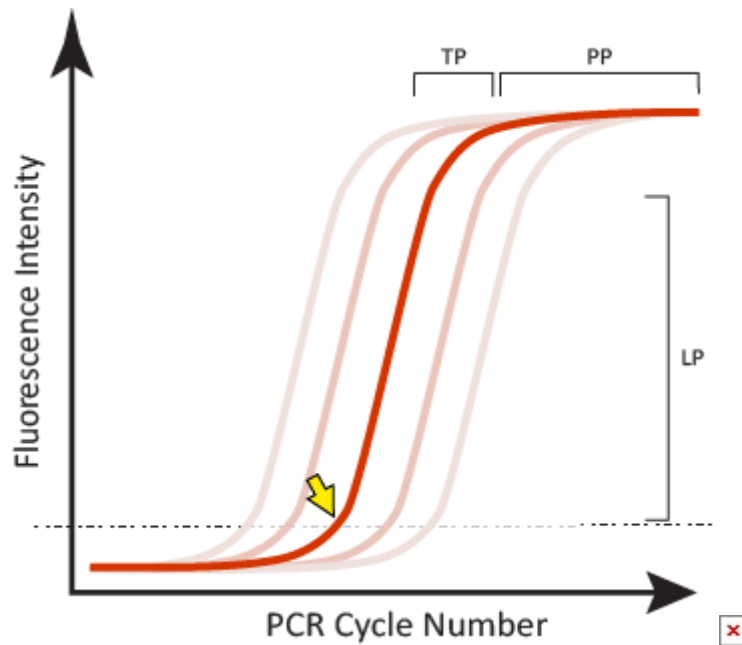

**Figure 1.** Examples of satisfactory sigmoidal amplification curve shape when considering an assay's fluorescent signal output. The crossing point or threshold cycle ( $C_T$ ) is indicated (yellow arrow); it is the value at which fluorescence levels surpass a predefined (usually set during validation, or arbitrary) threshold level as shown in this normalized linear scale depiction. LP-log-linear phase of signal generated during the exponential part of the PCR amplification; TP-a slowing of the amplification and accompanying fluorescence signal marks the transition phase; PP-the plateau phase is reached when there is little or no increase in fluorescent signal despite continued cycling.