

Nov 14,  
2019

## Enzymatic Assay of Trypsin Inhibition

Neilier Junior<sup>1</sup>

<sup>1</sup>Universidade Federal de Viçosa

1

Works for me

dx.doi.org/10.17504/protocols.io.9a4h2gw



Neilier Junior  
Universidade Federal de Viçosa



### ABSTRACT

This protocol has been standardized for analysis of protease inhibitors in soybean leaf extract, but can be easily adjusted for other biological samples.

### MATERIALS

NAME	CATALOG #	VENDOR
Tris		
DMSO	D1435	Sigma Aldrich
Calcium Chloride	C4904	Sigma Aldrich
Trypsin from bovine pancreas	T8003	Sigma-aldrich
N $\alpha$ -Benzoyl-L-arginine 4-nitroanilide hydrochloride (L-BAPNA)	B3279	Sigma-aldrich

### BEFORE STARTING

Check that all reagents and equipment are available. Plan the experiment!

### Prepare the solutions and the workspace

#### 1 Preparation of solutions

Trypsin solution:

Dilute 1.25 mg of bovine trypsin in 1 mL of water pH 3.0 (adjust with HCl).

L-BAPNA stock solution:

Dilute 130.47 mg L-BAPNA in 5 mL DMSO (concentration:  $6.0 \times 10^{-2}$  M). Store at -20 °C and protected from light.

L-BAPNA usage solution (freshly prepared):

Dilute 200  $\mu$ L of stock solution in 10 mL of 100 mM Tris-HCl buffer, pH 8.2 and 20 mM CaCl<sub>2</sub>. Protect from light.

- Separate three microtubes and name them "**blank**", "**control** (uninhibited test)" and "**test** (inhibited test)".  
Pipette the following reagents.

**Blank:** 500  $\mu$ L BAPNA usage solution and 500  $\mu$ L buffer.

**Control:** 100  $\mu$ L trypsin solution and 400  $\mu$ L buffer.

**Test:** 100  $\mu$ L enzyme, 100  $\mu$ L leaf extract (source of inhibitors) and 300  $\mu$ L buffer.

- Mix the three microtubes by inversion and equilibrate to 25°C for 5 min

Zero spectrophotometer with **blank** content at 410 nm

- 4 To the **control** microtube, add 500 µL of the BA<sub>p</sub>NA usage solution

Immediately mix by inversion and mark the time and pour the contents into a cuvette. After 30 s of reaction onset, monitor readings at 410 nm for 120 s

- 5 Add 500 µL of the BA<sub>p</sub>NA usage solution to the **test** tube

Immediately mix by inversion and mark the time and pour the contents into a cuvette. After 30 s of reaction onset, monitor readings at 410 nm for 120 s

## 6 Calculations

$$\% \text{ Inhibition} = (\Delta A_{410 \text{ nm}}^{\text{control}} - \Delta A_{410 \text{ nm}}^{\text{test}}) * 100 / (\Delta A_{410 \text{ nm}}^{\text{control}} - A_{410 \text{ nm}}^{\text{blank}})$$

or

$$\text{Trypsin Inhibitor Units / mL} = (\Delta A_{410 \text{ nm}}^{\text{control}} - \Delta A_{410 \text{ nm}}^{\text{test}}) / (8800 * \text{time} * \text{leaf extract volume})$$

$A_{410 \text{ nm}}^{\text{blank}}$  = Absorbance in the **blank** at 410 nm

$\Delta A_{410 \text{ nm}}^{\text{control}}$  = Absorbance variation in the **control** sample at 410 nm within 120 seconds


$\Delta A_{410 \text{ nm}}^{\text{test}}$  = Absorbance variation in the **test** sample at 410 nm within 120 seconds

8800 = extinction coefficient of *p*-nitroanilide at 410 nm

time = 120 seconds

leaf extract volume = Volume of inhibitor source used (in milliliters)

The presence of trypsin inhibitors in the leaf extract decreases the enzymatic activity and this inhibition can be represented in the Lineweaver-Burk graph, where, as the inhibitor concentration increases, the slope of the line also increases.

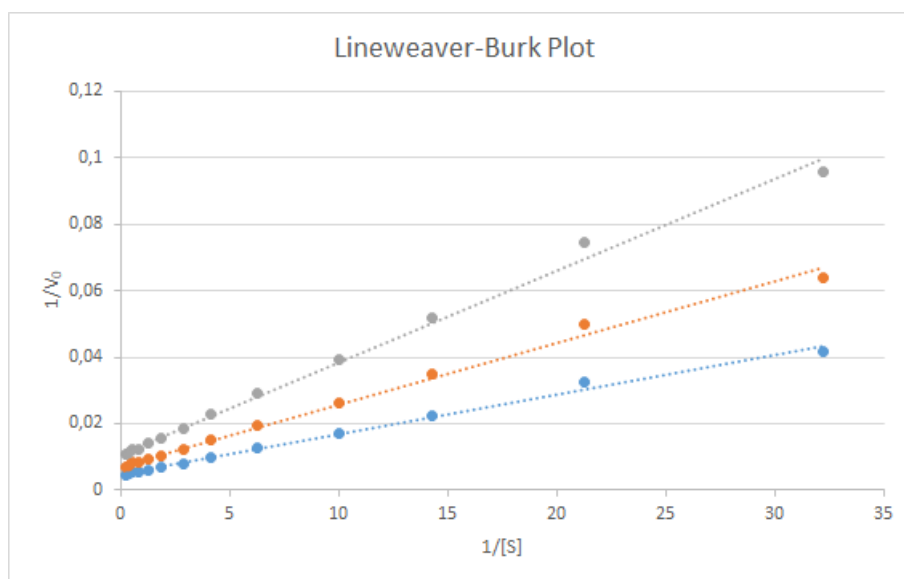

Lineweaver-Burk plot analysis of the inhibitory activity of soybean leaf extract toward trypsin. In blue, the kinetics in the absence of inhibitors (**control**). The kinetics in the presence of inhibitors (**test**) at  $0.5 K_i$  and  $K_i$  are represented in orange and gray, respectively.



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited