

Single nuclei isolation

Oct 28, 2019

Isolation of single nuclei from postmortem fresh frozen human brain and immunostaining for NeuN

Marcos Otero-Garcia¹, Inma Cobos¹¹Stanford University

1

Works for me

dx.doi.org/10.17504/protocols.io.6tuhenw

Neurodegeneration Method Development Community

Tech. support email: ndcn-help@chanzuckerberg.com

Inma Cobos



ABSTRACT

- Protocol based on Krishnaswami *et al.*, Nat Protoc. 2016, 3:499-524
- Used routinely in our lab for the isolation of single nuclei from fresh frozen postmortem brains from subjects with neurodegenerative disorders and healthy controls for single-nucleus RNA-seq
- Works also for isolation of nuclei from fresh frozen mouse brains
- Used for isolation of all nuclei or antibody-enriched populations (i.e., NeuN⁺ nuclei)
- We always perform and strongly recommend FACS if used for single-nucleus RNA-seq

THIS PROTOCOL ACCOMPANIES THE FOLLOWING PUBLICATION

Manuscript under review

Single nuclei isolation from postmortem frozen frozen brain_CobosLab_10282019.pdf

MATERIALS

NAME ▾	CATALOG # ▾	VENDOR ▾
DNase I (RNase-free) - 1,000 units	M0303S	New England Biolabs
Bovine Serum Albumin (BSA)	A7906	Sigma Aldrich
Magnesium chloride hexahydrate	M2670	Sigma Aldrich
UltraPure™ DNase/RNase-Free Distilled Water	10977023	Thermo Fisher Scientific
DAPI	D3571	Invitrogen - Thermo Fisher
Sucrose	S25590B	Fisher Scientific
Triton™ X-100	AC215680010	Fisher Scientific
Potassium chloride solution	60142	Millipore Sigma
UltraPure™ 1M Tris-HCl pH 8.0	15568025	Thermo Fisher Scientific
DL-Dithiothreitol	D9779-1G	Sigma Aldrich
cOmplete™ Protease Inhibitor Cocktail	11697498001	Sigma Aldrich
NxGen® RNase Inhibitor	30281-2	Lucigen
Anti-NeuN Antibody	MAB377	Merck Millipore


NAME ▾	CATALOG # ▾	VENDOR ▾
Goat anti-Mouse IgG (H L) Alexa Fluor 647	A-21235	Thermo Fisher Scientific
OptiPrep™ Density Gradient Medium	D1556)	Sigma Aldrich

STEPS MATERIALS

NAME ▾	CATALOG # ▾	VENDOR ▾
Anti-NeuN Antibody, clone A60	MAB377	Merck Millipore

MATERIALS TEXT

- Material and tools needed: Forceps, spatula, blades, dounce tissue grinder, cell strainers, petri dishes



Dounce All-Glass Tissue Grinders
Tissue Grinder
KIMBLE 885300-0007 [↗](#)

Buffer preparation (All solutions should be RNase-free for single-soma RNAseq experiments):

Isolation Medium #1 (IM1), 45 ml (optional)

Prepare in a 50 ml Falcon tube and store at **4 °C** up to 6 months.

Amount	Reagent	Final concentration
7500 µl	1.5 M Sucrose	250 mM
1125 µl	1 M KCl	25 mM
225 µl	1M MgCl ₂	5 mM
450 µl	1 M Tris (pH 8)	10 mM
35.7 ml	RNase-free H ₂ O	

Homogenization Buffer (3 ml per sample)

Prepare FRESH and keep ICED or at **4 °C** . Discard after use.

Amount	Reagent	Final concentration
2925 µl	IM1	
3 µl	DTT 1mM	1 µM
30 µl	50x Protease Inhibitor	0.5 x
15 µl	RNaseIN 40U/µl	0.2 U/µl
30 µl	TritonX100 10 %v/v	0.1 %

Iodixanol dilutions

Prepare in a 50 ml Falcon tubes and store at **4 °C** up to 6 months. Accuracy with Iodixanol and sucrose concentrations is critical.

■ Iodixanol medium (IDM), 45 ml

Amount	Reagent	Final concentration
7500 µl	1.5 M Sucrose	250 mM
6750 µl	1 M KCl	150 mM
1350 µl	1 M MgCl ₂	30 mM
2700 µl	1 M Tris (pH 8)	60 mM

26.7 ml	Nuclease-free H2O	
---------	-------------------	--

■ **Iodixanol 50%(v/v), 20 ml**

16.7 ml Iodixanol 60 %(v/v) + 3.3 ml IDM

■ **Iodixanol 29%(v/v), 30 ml**

14.5 ml Iodixanol 60 %(v/v) + 15.5 ml IDM

Freezing Storage Buffer (FSB), 15ml (optional)

Prepare in a 50 ml Falcon tube and store at  **4 °C** up to 6 months.

<i>Amount</i>	<i>Reagent</i>	<i>Final concentration</i>
1.665 ml	1.5M Sucrose	166.5 mM
75 µl	1M MgCl ₂	5 mM
150 µl	1M Tris (pH8)	10 mM
13.1 ml	Nuclease-free H2O	

Buffer for Immunostaining, 10ml

<i>Amount</i>	<i>Reagent</i>	<i>Final concentration</i>
9400 µl	RNase-free PBS (pH7.4)	
500 µl	RNase-free 10 % BSA	0.5 %
50 µl	1 M MgCl ₂	5 mM
10 µl	DNAse I (2000 U/ml)	2 U/ml

Buffer for Antibody Incubation

Add to 1 ml of Buffer for Immunostaining:

<i>Amount</i>	<i>Reagent</i>	<i>Final concentration</i>
5 µl	RNaseI 40 U/µl	0.2 U/µl

Collection medium for FACS (1 ml; 0.2 ml/vial)

<i>Amount</i>	<i>Reagent</i>	<i>Final concentration</i>
950 µl	RNase-free PBS (pH7.4)	
25 µl	RNase IN	
*	BSA 10% *After collecting!	1 %

SAFETY WARNINGS

Please see SDS (Safety Data Sheet) for hazards and safety warnings.
Uses fresh human brain tissue - Biosafety Level 2 lab work.

BEFORE STARTING

- All solutions and materials should be **RNase-free** and **kept iced or at 4 °C** at all times
- Glassware and metal tools are **sealed with aluminum foil** and **baked at 220 °C for 06:00:00**.

1 

Prepare **Homogenization Buffer** and cool it on ice.

Amount	Reagent	Final concentration
2925 µl	IM1	
3 µl	DTT 1mM	1 µM
30 µl	50x Protease Inhibitor	0.5 x
15 µl	RNaseIN 40U/µl	0.2 U/µl
30 µl	TritonX100 10 %v/v	0.1 %

Homogenization Buffer (3 ml per sample)

Prepare FRESH and keep ICED or at **4 °C**. Discard after use.



1.1 Add **2925 µl IM1**.

1.2 Add **3 µl DTT 1mM**.

1.3 Add **30 µl 50x Protease Inhibitor**.










1.4 Add **15 µl RNaseIN 40U/µl**.

1.5 Add **30 µl TritonX100 10 %v/v**.

- 2 Pre-cool the dounce tissue grinder on ice (**Kimble Kontes all glass tissue grinder**, 7 mL tubes, 0.02 - 0.056 mm clearance space between pestle and tubes)
- 3 Add  **2.4 ml** of **Homogenization buffer** to the dounce tissue grinder.
- 4 Collect the brain chunk (~  **100 mg**) and transfer to a Petri dish on ice.
- 5 Cut out into small pieces using a chilled scalpel or blade.
- 6 Transfer all pieces into the dounce tissue grinder.
- 7 Homogenize the tissue, on ice (~**30 firm strokes**).
Check on **hematocytometer** while homogenizing and adjust the number of strokes.



Too few strokes – you will see cell clumps; too many strokes – you will get damaged nuclei.

- 8 Filter homogenate using **40 µm** Corning cell strainer to remove clumps.
--> Take sample for hemocytometer.
- 9 Transfer the homogenate into two precooled 1.5 ml Eppendorf tubes.
- 10 
Centrifuge at  **1000 x g** for  **00:08:00** at  **4 °C** .
- 11 Slowly aspirate the supernatant from the side of the tube.
Avoid disturbing the pellets, up to  **50 µl** of supernatant can be left.
- 12 Gently **resuspend** each pellet in  **225 µl** (final volume) of cold **Homogenization buffer** and pool both tubes (final volume  **450 µl**).
- 13 
Add an **equal** volume of cold  **50 %v/v iodixanol** .



Critical! Be very exact with volumes.  **450 µl** suspension +  **450 µl** 50 % iod.
Final iodixanol concentration is **25 %**.

14 

Gently pipette mix.

15 Add an equal volume of **29 % iodixanol** ( **900 µl**) into a 2 ml Eppendorf tube.

16 Slowly layer off the 25 % iodixanol suspension mix over the 25 % iodixanol, without mixing them.

17 

Centrifuge at  **13500 x g** for  **00:20:00** at  **4 °C** .

18 

Prepare **buffer for immunostaining**.



<i>Amount</i>	<i>Reagent</i>	<i>Final concentration</i>
9400 µl	RNase-free PBS (pH7.4)	
500 µl	RNase-free 10 % BSA	0.5 %
50 µl	1 M MgCl ₂	5 mM
10 µl	DNAse I (2000 U/ml)	2 U/ml

18.1 Add  **9400 µl RNase-free PBS (pH7.4)** .

18.2 Add  **500 µl RNase-free 10 % BSA** .

18.3 Add  **50 µl 1 M MgCl₂** .

18.4 Add  **10 µl DNAse I (2000 U/ml)** .

19 Prepare **buffer for Antibody incubation** (1.5 ml per sample) by adding  **5 µl** of RNaseIN 40 U/µl to  **1 ml** of **buffer for Immunostaining** (final concentration 0.2U/µl)

20 After centrifugation, remove and discard the top myelin-rich debris layer (you can use a 1 ml pipette with the tip cut, or cotton swabs).



21 

Then, remove and discard the aqueous supernatant, without disrupting the **nuclei pellet**. Avoid contaminating with the top layer.

22 Use a small amount of **buffer for Antibody Incubation** to resuspend the pellet and transfer the solution to a new tube.

23 Gently resuspend in  **500 µl** of **buffer for Antibody Incubation**.

24 

Incubate for  **00:15:00** , at  **4 °C** or iced, for **blocking** nonspecific staining.

25 Take sample for hemocytometer.

26 Take a sample for **unstained** control. Take sample for **2AB-only** (nonspecific binding control).

27 Add primary antibody (**Ms-a-NeuN, 1:1000**) and incubate on a rotator in cold room ( **4 °C**) for  **00:40:00** .





Anti-NeuN Antibody, clone A60

by Merck Millipore

Catalog #: **MAB377**



Use Eppendorfs coated with BSA for collection (to coat the Eppendorf tubes, fill them with 10 % BSA solution in PBS for  **00:05:00** , rinse with PBS, and dry at  **4 °C** overnight).

Washing


28 

Add  **700 µl** of **buffer for Immunostaining** and invert several times.


29 


Centrifuge at  **500 x g** for  **00:05:00** at  **4 °C** .

30 

Carefully transfer all the supernatant (don't leave any supernatant! It's ok to drag some pellet) to a new Eppendorf tube and resuspend the pellet in  **400 µl** of **buffer for Antibody Incubation**.

31 

Centrifuge the supernatant at  **500 x g** for  **00:05:00** at  **4 °C** to recover non-pelleted nuclei.

32 

Resuspend them in  **200 µl** of **buffer for Antibody Incubation** and pool for a final volume of 600 µl.

33 Take samples for single color staining.

34 Add secondary antibody (**G-a-Ms Alexa Fluor 647, 1:1000**) and **Hoechst 33258** (5 µg/ml) to the same tube. (**Hoechst 1:1000 of 2,5 mg/ml stock**).

35 Check the quality of the sample and on hemocytometer.



For 100 mg of tissue, should have ~ 2 Million nuclei in ~0.6 ml.

36 Transfer to 7 ml culture tubes.

37 Place on ice and bring them to the FACS facility for sorting.



FACS sorting controls:

Non-stained

2AB only


Single stainings

38 Collect in BSA coated Eppendorfs containing **collection medium** (containing ~1/5 of the expected final volume after collection).

39 Add BSA after collection for final 1% BSA final concentration (tested for 10x Chromium v2 and v3).



Different Single Cell assays tolerate different BSA concentrations, but lowering it may increase nuclei aggregation.

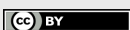

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited