



Oct 03, 2019

## MojoSort™ Isolation Kits Protocol - 1 V.2

Sam Li<sup>1</sup><sup>1</sup>BioLegend

1

Works for me

dx.doi.org/10.17504/protocols.io.7x6hpre

BioLegend

Sam Li  
BioLegend

## ABSTRACT

**Product description and procedure summary:**

Target cells are depleted by incubating the sample with the biotin antibody cocktail followed by incubation with magnetic Streptavidin Nanobeads (Cat. No. [480015/480016](#)). The magnetically labeled fraction is retained by the use of a magnetic separator. The untouched cells are collected. These are the cells of interest; do not discard the liquid. Some of the downstream applications include functional assays, gene expression, phenotypic characterization, etc.

**Note:** This protocol has been optimized to remove washing steps after antibody cocktail and nanobeads incubations, resulting in a shorter and more convenient protocol. This procedure is optimized for the isolation of  $10^7$  to  $2 \times 10^8$  cells per tube. If working with fewer than  $10^7$  cells, keep volumes as indicated for  $10^7$  cells. For best results, optimize the conditions to your specific cell number and tissue. Prepare fresh MojoSort™ Buffer solution by diluting the 5X concentrate with sterile distilled water. **Scale up volumes if using 14 mL tubes and Magnet, and place the tube in the magnet for 10 minutes.**

## EXTERNAL LINK

<https://www.biolegend.com/protocols/mojosort-isolation-kits-protocol-1/4599/>

## GUIDELINES

**Important Note**

MojoSort™ magnetic particles can be used with other commercially available magnetic separators, both free standing magnets and column-based systems. Because MojoSort™ protocols are optimized for the MojoSort™ separator, the protocols may need to be adjusted for other systems. Please contact BioLegend Technical Service ([tech@biolegend.com](mailto:tech@biolegend.com)) for more information and guidance. We do not recommend using MojoSort™ particles for BD's IMag™ or Life Technologies' DynaMag™.

**Application notes:** To use this product in magnetic separation columns, a titration of the Nanobeads should be performed. Optimal concentration for magnetic separation columns is lot-specific. Please contact BioLegend Technical Service ([tech@biolegend.com](mailto:tech@biolegend.com)) for further assistance on how to use MojoSort™ Nanobeads in magnetic separation columns.

## MATERIALS

NAME	CATALOG #	VENDOR
MojoSort™ Buffer	480017	BioLegend
MojoSort™ Magnet	480019	BioLegend
MojoSort™ Mouse CD3 T Cell Isolation Kit	480023, 480024, 480031	BioLegend
MojoSort™ Mouse CD4 T Cell Isolation Kit	480005, 480006, 480033	BioLegend
MojoSort™ Mouse CD8 T Cell Isolation Kit	480007, 480008, 480035	BioLegend
MojoSort™ Human CD4 T Cell Isolation Kit	480009, 480010	BioLegend
MojoSort™ Human CD3 T Cell Isolation Kit	480021, 480022	BioLegend
MojoSort™ Mouse CD4 Naïve T Cell Isolation Kit	480039, 480040	BioLegend


NAME ▾	CATALOG # ▾	VENDOR ▾
MojoSort™ Human CD4 Naïve T Cell Isolation Kit	480041, 480042	BioLegend
MojoSort™ Mouse CD8 Naïve T Cell Isolation Kit	480043, 480044	BioLegend
MojoSort™ Human CD8 Naïve T Cell Isolation Kit	480045, 480046	BioLegend
MojoSort™ Mouse Pan B Cell Isolation Kit	480051, 480052	BioLegend
MojoSort™ Human B Cell (CD43-) Isolation Kit	480061, 480062	BioLegend
MojoSort™ Human CD4 Memory T Cell Isolation Kit	480063, 480064	BioLegend
MojoSort™ Human Naïve B Cell Isolation Kit	480067, 480068	BioLegend
MojoSort™ Human Pan B Cell Isolation Kit	480081, 480082	BioLegend
MojoSort™ Mouse Pan B Cell Isolation Kit II	480087, 480088	BioLegend

#### MATERIALS TEXT

- Adjustable pipettes
- 70µm filters (one per sample)
- 5mL (12 x 75mm) or 14mL (17 x 100mm) polypropylene tubes
- Reagents for sample preparation
- Reagents and instruments (Flow cytometer) to determine yield and purity

- 1 Prepare cells from your tissue of interest or blood without lysing erythrocytes. Kits for human samples have been optimized for PBMCs, please prepare the cells using a suitable method.
- 2 In the final wash of your sample preparation, resuspend the cells in MojoSort™ Buffer by adding up to 4mL in a 5 mL (12 x 75 mm) polypropylene tube.  
**Note:** Keep MojoSort™ Buffer on ice throughout the procedure.
- 3 Filter the cells with a 70µm cell strainer, centrifuge at 300xg for 5 minutes, and resuspend in an appropriate volume of MojoSort™ Buffer. Count and adjust the cell concentration to  $1 \times 10^8$  cells/mL.
- 4 Aliquot 100µL of cell suspension ( $10^7$  cells) into a new tube. **Add 10µL of the Biotin-Antibody Cocktail.** Mix well and **incubate on ice for 15 minutes.** Scale up the volume accordingly if separating more cells. For example, add 100µL of biotinylated antibody cocktail for separating  $1 \times 10^8$  cells in 1 ml of MojoSort™ Buffer. When working with less than  $10^7$  cells, use indicated volumes for  $10^7$  cells.  
*Optional: Take an aliquot before adding the cocktail to monitor purity and yield.*
- 5 Resuspend the beads by vortexing, maximum speed, 5 touches. **Add 10µL of Streptavidin Nanobeads.** Mix well and **incubate on ice for 15 minutes.** Scale up the volume accordingly if separating more cells. For example, add 100µL of Nanobeads for separating  $1 \times 10^8$  cells in 1 ml of MojoSort™ Buffer. When working with less than  $10^7$  cells, use indicated volumes for  $10^7$  cells.
- 6 Add 2.5mL of MojoSort™ Buffer.  
**Note:** If you observe aggregates, filter the suspension. To maximize yield, you can disrupt the aggregates by pipetting the solution up and down.

- 7 Place the tube in the magnet for 5 minutes.  
*Optional: Take a small aliquot before placing the tube in the magnet to monitor purity and yield. Keep unused cells to be used as control or other applications if needed.*
- 8 Pour out and collect the liquid. These are your cells of interest; **DO NOT DISCARD**.
- 9 Repeat steps 6-8 with labeled cells once more for a total of **2 separations**. Pool the unlabeled fractions. The labeled cells may be useful as staining controls, to monitor purity/yield, or other purposes.
- Note:** Repeating the magnetic separation increases the yield, without a strong impact on the purity. The yield will typically increase about 8-10% with a second separation. The purity may decrease 1-2% with each separation.

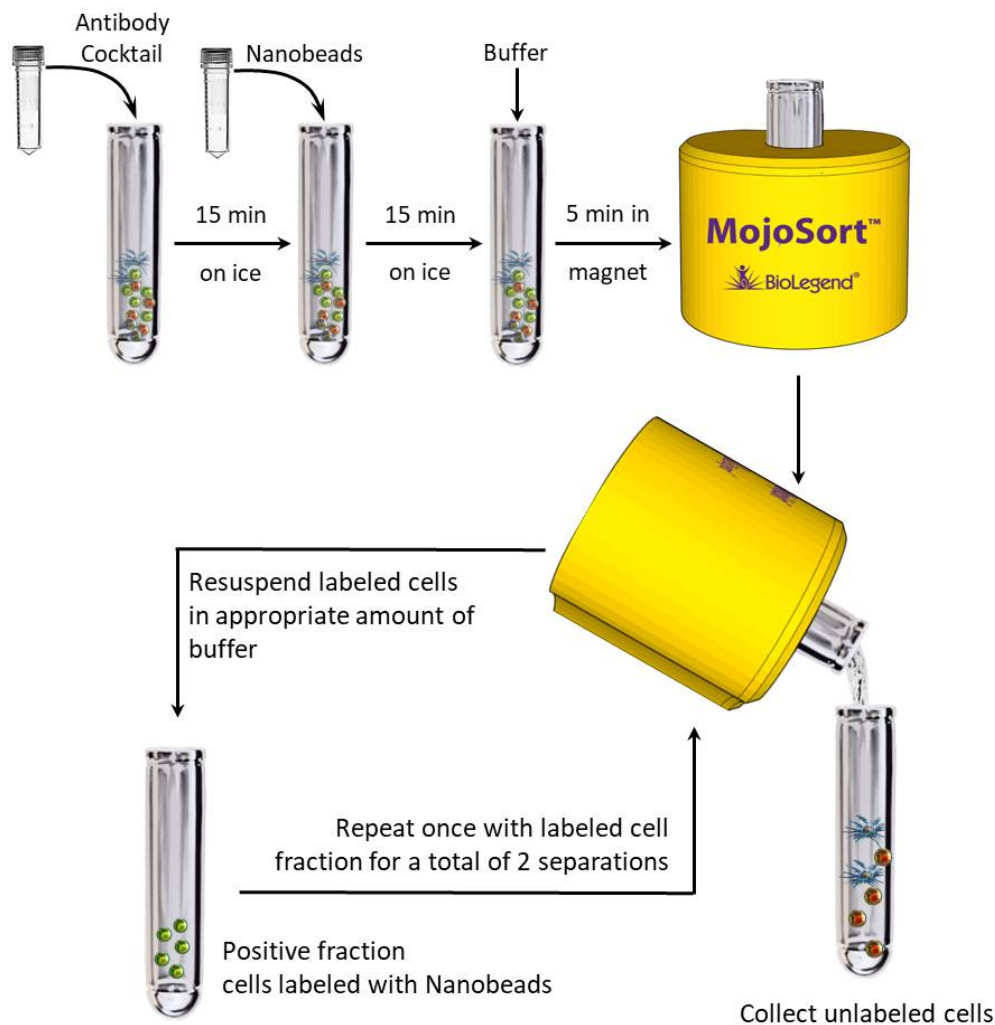

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited