

Human Parechovirus A real-time RT-PCR ["Nix assay"; 2008-2015]

Ian Mackay

Abstract

I and my team used this assay between 2008-2015; we dubbed it the "Nix assay". It targets the 5'UTR and employs quite degenerate oligos.

In silico sequence alignments indicated the oligonucleotides could theoretically detect at least HPeV 1-7, 17 and 18.

However during a period of assay comparison, another assay (see link below), the "Benschop assay" ([J.Clin.Virol. 2008. 41\(2\):69-74](#)), was found to produce more sigmoidal and higher curves and 1-5 cycle improvements to C_T values when compared among the sample sample extract set.

Citation: Ian Mackay Human Parechovirus A real-time RT-PCR ["Nix assay"; 2008-2015]. **protocols.io**

dx.doi.org/10.17504/protocols.io.k2jcycn

Published: 28 Nov 2017

Guidelines

- This protocol assumes the user is familiar with working in a laboratory, with PCR, the thermocycler and software used to run it
- This protocol should be re-evaluated if being used with different reagents, if the oligonucleotide sequences are changed or if the cycling conditions are altered

Protocol

Oligonucleotides...

Step 1.

Name	5'-3' oligonucleotide sequence
AN345_panHPeV/LV (sense primer)	GTAACASWWGCCTCTGGGSCCAAAAG
AN344_panHPeV/LV (antisense probe)	GGCCCCWGRTCAGATCCAYAGT
AN257_HPeV/LV (probe)	FAM-CCTRYGGGTACCTYCWGGGCATCCTTC-BHQ1

Reagents

Step 2.



REAGENTS

SensiFAST™ Probe Hi-ROX One-Step Kit [BIO-77001](#) by [Bioline](#)

Reaction setup...

Step 3.

Below is the reaction setup for a single RT-PCR reaction.

Ideally, this work is conducted in a laboratory separate to any space used to perform PCR, molecular cloning or the analysis of high concentration DNA.

This volume has been used in 0.1-0.2ml tubes or various other connected tube configurations such as 100-place rings.

Multiply this according to the number of reactions you will need, remembering to include a positive control and at least two non-template controls (NTCs)

You may also need to allow some extra volume, depending on the method used to pipette mix into tubes for the run. For example, some robot-loaded tubes can require two reaction 'dead volumes'.

Reagent (stock concentration)	Vol (μL) / reaction	Final concentration
Nuclease free water	0.8	N/A
Primers (combined at 2pmol/μl)	4	400nM
Probe (2pmol/μl)	1	100nM
MgCl ₂ (25mM)	1.6	2mM (+3mM in 2X mix)
2X Reaction Mix ¹	10	1X
RNase inhibitor	0.4	50nM
RT/Taq Mix ¹	0.2	1X
Template extract RNA	2	N/A
Final volume	20μl	

¹SensiFast OneStep Mix (2X), Cat No. 11732088

Amplification...

Step 4.

This assay has been optimized and validated for use with a Rotor-Gene 6000 or Rotor-Gene Q thermal cycler.

The cycling conditions are as follows:

RT-PCR		
45°C	20min	
94°C	2min	
94°C	15s	55X
58°C	30s*	
72°C	10s	

*Florescence acquisition step

Result calling...

Step 5.

The definition used for a satisfactory positive result from a real-time fluorogenic PCR should include each of the following:

1. A **sigmoidal curve** - the trace travels horizontally, curves upward, continues in an exponential rise, curves and reaches a horizontal plateau phase
2. A **suitable level of fluorescence** intensity as measured in comparison to a positive control (y-axis)
3. A defined threshold (C_T) value which the fluorescent curve has clearly exceeded (Fig.1 arrow) and which sits early in the log-linear phase
4. A flat or non-sigmoidal curve or a curve that crosses the threshold with a C_T value >40 cycles is considered a negative result
5. No template controls (NTCs; water instead of specimen extract) should not produce a curve

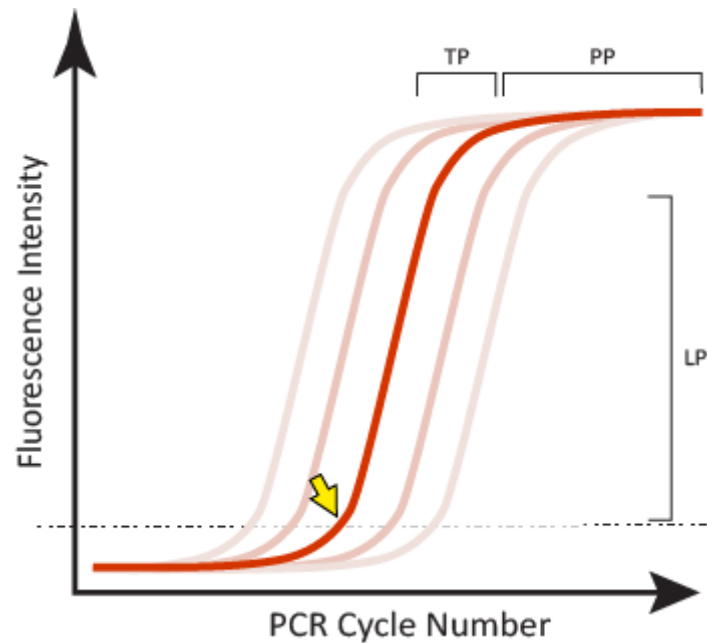


Figure 1. Examples of satisfactory sigmoidal amplification curve shape when considering an assay's fluorescent signal output. The crossing point or threshold cycle (C_T) is indicated (yellow arrow); it is the value at which fluorescence levels surpass a predefined (usually set during validation, or arbitrary) threshold level as shown in this normalized linear scale depiction. LP-log-linear phase of signal generated during the exponential part of the PCR amplification; TP-a slowing of the amplification and accompanying fluorescence signal marks the transition phase; PP-the plateau phase is reached when there is little or no increase in fluorescent signal despite continued cycling.