

Protocol for use with FFPE RNA, NEBNext rRNA Depletion Kit (Human/Mouse/Rat) (NEB #E6310) and NEBNext Ultra II Directional RNA Library Prep Kit for Illumina (E7760, E7765)

New England Biolabs¹

¹New England Biolabs

dx.doi.org/10.17504/protocols.io.tumenu6

New England Biolabs (NEB)

Tech. support phone: +1(800)632-7799 email: info@neb.com

ABSTRACT

The NEBNext Ultra II Directional RNA Library Prep Kit for Illumina contains the enzymes and buffers required to convert a broad range of input amounts of RNA into high quality directional (strand-specific) libraries for next-generation sequencing on the Illumina platform. The fast, user-friendly workflow has minimal hands-on time and is compatible with poly(A) mRNA enrichment and rRNA depletion methods.

PROTOCOL STATUS

Working

GUIDELINES

Section 3

RNA Sample Requirements

RNA Integrity:

Assess the quality of the input RNA by running the RNA sample on an Agilent Bioanalyzer RNA 6000 Nano/Pico Chip to determine the RNA Integrity Number (RIN). RNA with different RIN values require different fragmentation times or no fragmentation at all.

For intact (RIN > 7) or partially degraded RNA samples (RIN = 2 to 7) follow the library preparation protocol in Section 2. See Table 2.5.1. for the recommended fragmentation times, based on RIN.






For highly degraded samples (RIN = 1 to 2) (e.g. FFPE), which do not require fragmentation, follow the library preparation protocol in Section 3 (current Section).






RNA Sample Requirements:

The RNA sample should be free of salts (e.g., Mg^{2+} , or guanidinium salts) or organics (e.g., phenol and ethanol). RNA must be free of DNA. gDNA is a common contaminant from RNA preps. It may be carried over from the interphase of organic extractions or when the silica matrix of solid phase RNA purification methods is overloaded. If the total RNA sample may contain gDNA contamination, treat the sample with DNase I to remove all traces of DNA (not provided in this kit). After treatment with DNase I the enzyme should be removed from the sample. Any residual activity of DNase I will degrade the single stranded DNA probes necessary for the ribosomal depletion.


















DNase I can be removed from the extraction using phenol/ chloroform extraction and ethanol precipitation.


MATERIALS

NAME ▾	CATALOG # ▾	VENDOR ▾
NEBNext RNase H	E6318	 New England Biolabs
RNase H Reaction Buffer	E6312	 New England Biolabs
NEBNext rRNA Depletion Solution	E6313	 New England Biolabs
NEBNext Probe Hybridization Buffer	E6314	 New England Biolabs
DNase I (RNase-free)	E6316	 New England Biolabs

NAME ▾	CATALOG # ▾	VENDOR ▾
DNase I Reaction Buffer	E6315	 New England Biolabs
Nuclease-free Water	E6317	 New England Biolabs
NEBNext RNA Sample Purification Beads	E6315	 New England Biolabs
Magnetic Rack	View	
80% Ethanol (freshly prepared)	View	
Thermal cycler	View	
Agencourt RNAClean XP Beads	A63987	Beckman Coulter
DNase I (e.g., NEB #M0303) and DNase I Cleanup Reagents or Kit for Removal of DNA Prior to Depletion	View	 New England Biolabs
Random Primers	E7422	 New England Biolabs

STEPS MATERIALS

NAME ▾	CATALOG # ▾	VENDOR ▾
NEBNext RNase H	E6318	 New England Biolabs
RNase H Reaction Buffer	E6312	 New England Biolabs
Nuclease-free Water	E6317	 New England Biolabs
DNase I Reaction Buffer	E6315	 New England Biolabs
DNase I (RNase-free)	E6316	 New England Biolabs
Nuclease-free Water	E6317	 New England Biolabs
NEBNext Sample Purification Beads	E6315	 New England Biolabs
Nuclease-free Water	E6317	 New England Biolabs
Random Primers	E7422	 New England Biolabs
NEBNext rRNA Depletion Solution	E6313	 New England Biolabs
NEBNext Probe Hybridization Buffer	E6314	 New England Biolabs
80% Ethanol (freshly prepared)		
NEBNext Strand Specificity Reagent	E7766	 New England Biolabs
NEBNext First Strand Synthesis Enzyme Mix	E7761	 New England Biolabs
NEBNext Second Strand Synthesis Reaction Buffer with dUTP Mix	E7426	 New England Biolabs
NEBNext Second Strand Synthesis Enzyme Mix	E7425	 New England Biolabs
Nuclease-free Water	E7764	 New England Biolabs
NEBNext Sample Purification Beads	E7767	 New England Biolabs

NAME ▾	CATALOG # ▾	VENDOR ▾
Fresh 80% Ethanol		
(0.1X) TE Buffer	E7763	 New England Biolabs
NEBNext Ultra II End Prep Reaction Buffer	E7647	 New England Biolabs
NEBNext Ultra II End Prep Enzyme Mix	E7646	 New England Biolabs
NEBNext Ligation Enhancer	E7374	 New England Biolabs
NEBNext Ultra II Ligation Master Mix	E7648	 New England Biolabs
NEBNext USER Enzyme	E7458	 New England Biolabs
NEBNext Sample Purification Beads	E7767	 New England Biolabs
80% Ethanol (freshly prepared)		
NEBNext Sample Purification Beads	E6315	 New England Biolabs
80% Ethanol (freshly prepared)		

SAFETY WARNINGS

Please refer to the SDS (Safety Data Sheet) for safety warnings and hazard information.

BEFORE STARTING


Input Amount Requirement

10 ng – 100 ng FFPE RNA in up to 12 µl of Nuclease-free Water, quantified by Qubit Fluorometer and quality checked by Bioanalyzer.

The protocol is optimized for approximately 200 nt RNA inserts.

Keep all of the buffers on ice, unless otherwise indicated.

Probe Hybridization to RNA

1 Dilute the total RNA with Nuclease-free Water to a final volume of  12 µl in a PCR tube. Keep the RNA on ice.

2 Prepare a RNA/Probe master mix as follows:





RNA Probe Master Mix	Volume
NEBNext rRNA Depletion Solution	1 µl
Probe Hybridization Buffer	2 µl
Total Volume	3 µl



NEBNext rRNA Depletion Solution
by [New England Biolabs](#)
Catalog #: [E6313](#)



NEBNext Probe Hybridization Buffer
by [New England Biolabs](#)
Catalog #: [E6314](#)

- 3 Add  3 µl of the above mix to  12 µl total RNA (from Step 1), resulting in a total volume of  15 µl
- 4 Mix by pipetting up and down at least 10 times.
- 5 Spin down briefly in a microcentrifuge.
- 6 Place samples in a thermocycler with a heated lid set to approximately  105 °C, and run the following program, which will take approximately 15–20 minutes to complete:

Temperature	Time
95°C	2 min
Ramp down to 22°C	0.1°C/sec
Hold at 22°C	5 minutes
- 7 Spin down the samples in a microcentrifuge, and place on ice. Proceed immediately to RNase H Digestion Step.

RNase H Digestion

- 8 Assemble the RNase H master mix **on ice** as follows.

	Volume
RNase H Master Mix	
NEBNext RNase H	2 µl
RNase H Reaction Buffer	2 µl
Nuclease-free Water	1 µl
Total Volume	5 µl



NEBNext RNase H
by [New England Biolabs](#)
Catalog #: [E6318](#)








RNase H Reaction Buffer
by [New England Biolabs](#)
Catalog #: [E6312](#)



Nuclease-free Water
by [New England Biolabs](#)
Catalog #: [E6317](#)


- 9 Mix thoroughly by pipetting up and down 10 times.
- 10 Briefly spin down the samples in a microcentrifuge.

- 11 Add  3 µl of the RNase H master mix to the RNA sample from Step 7, resulting in a total volume of  20 µl
- 12 Mix thoroughly by pipetting up and down at least 10 times.
- 13 Incubate the sample in a thermocycler for  00:30:00 at  37 °C with the lid set to  40 °C (or on).
- 14 Briefly spin down the samples in a microcentrifuge, and place on ice. Proceed immediately to DNase I Digestion to prevent non-specific degradation of RNA.


DNase I Digestion

- 15 Assemble the DNase I master mix **on ice** in a nuclease-free tube.


DNase I Master Mix	Volume
DNase I Reaction Buffer	5 µl
DNase I (RNase-free)	2.5 µl
Nuclease-free Water	22.5 µl
Total Volume	30 µl









DNase I Reaction Buffer
by [New England Biolabs](#)
Catalog #: [E6315](#)



DNase I (RNase-free)
by [New England Biolabs](#)
Catalog #: [E6316](#)



Nuclease-free Water
by [New England Biolabs](#)
Catalog #: [E6317](#)


- 16 Mix thoroughly by pipetting up and down 10 times.
- 17 Briefly spin down the sample in a microcentrifuge.
- 18 Add  30 µl of DNase I master mix to  20 µl RNA sample from Step 14, resulting in a total volume of  50 µl
- 19 Mix thoroughly by pipetting up and down 10 times.
- 20 Incubate the sample in a thermocycler for  00:30:00 at  37 °C with the heated lid set to  40 °C (or on).

21 Briefly spin down the sample in a microcentrifuge, and place on ice. Proceed immediately to RNA Purification.

RNA Purification after rRNA Depletion Using Agencourt RNAClean XP Beads or NEBNext RNA Sample Purification Beads

22 Vortex Agencourt RNAClean XP Beads or **RNA** Sample Purification Beads to resuspend.



23 Add  **110 µl** (2.2X) beads to the RNA sample from Step 21 and mix thoroughly by pipetting up and down at least 10 times.




NEBNext Sample Purification Beads
by [New England Biolabs](#)
Catalog #: [E6315](#)

24 Incubate the sample for  **00:15:00** on ice to bind RNA to the beads.

25 Place the tube on a magnetic rack to separate beads from the supernatant. After the solution is clear, carefully remove and discard the supernatant. Be careful not to disturb the beads, which contain RNA.

26 Add  **200 µl** of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for  **00:00:30**, and then carefully remove and discard the supernatant. Be careful not to disturb the beads that contain the RNA.




80% Ethanol (freshly prepared)


27 Repeat Step 26 once for a total of 2 washing steps.

 [go to step #26 Repeat Step](#)


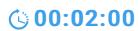
28 Completely remove residual ethanol, and air dry the beads for up to 5 minutes while the tube is on the magnetic rack with the lid open.

Caution: Do not over-dry the beads. This may result in lower recovery of RNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.

29 Remove the tube from the magnet. Elute the RNA from the beads by adding  **7 µl** Nuclease-free Water. Mix well by pipetting up and down at least 10 times and briefly spin the tube.



Nuclease-free Water
by [New England Biolabs](#)
Catalog #: [E6317](#)

30 Incubate for  00:02:00 at room temperature. Place the tube in the magnet until the solution is clear (~  00:02:00).


31 Remove  5 µl of the supernatant containing RNA and transfer to a nuclease-free tube.

32 Place the sample on ice and proceed to Priming of Highly Degraded RNA.

Priming of Highly Degraded RNA (FFPE) Which has a RIN ≤ 2 and Does not Require Fragmentation

33 Assemble the following fragmentation and priming reaction **on ice**:

Priming Reaction	Volume
rRNA Depleted Sample (Step 32)	5 µl
Random Primers	1 µl
Total Volume	6 µl



Random Primers
by [New England Biolabs](#)
Catalog #: [E7422](#)

34 Mix thoroughly by pipetting up and down 10 times.

35 Briefly spin down the samples in a microcentrifuge.

36 Incubate the sample in a preheated thermocycler as follows:

 00:05:00 at  65 °C , with heated lid set at  105 °C


Hold at  4 °C

37 Transfer the tube directly to ice and proceed to First Strand cDNA Synthesis.


First Strand cDNA Synthesis

38 Assemble the first strand synthesis reaction **on ice** by adding the following components to the fragmented and primed RNA from Step 37:

First Strand Synthesis Reaction	Volume
Primed RNA (Step 37)	6 µl
NEBNext First Strand Synthesis Reaction Buffer	4 µl
NEBNext Strand Specificity Reagent	8 µl
NEBNext First Strand Synthesis Enzyme Mix	2 µl
Total Volume	20 µl



NEBNext Strand Specificity Reagent
by [New England Biolabs](#)
Catalog #: [E7766](#)



NEBNext First Strand Synthesis Enzyme Mix
by [New England Biolabs](#)
Catalog #: E7761

39 Mix thoroughly by pipetting up and down 10 times.

40 [!] Incubate the sample in a preheated thermocycler with the heated lid set at $\geq 80^{\circ}\text{C}$ as follows:

Note: If you are following recommendations in Appendix A (Chapter 6), for libraries with longer inserts (> 200 bases), increase the incubation at 42°C from 15 minutes to 50 minutes at Step 2 below.

Step 1: 00:10:00 at 25°C

Step 2: 00:15:00 at 42°C

Step 3: 00:15:00 at 70°C


Step 4: Hold at 4°C

41 Proceed directly to Second Strand cDNA Synthesis.


Second Strand cDNA Synthesis

42 Assemble the second strand cDNA synthesis reaction **on ice** by adding the following components into the first strand synthesis product from Step 40.


Second Strand Synthesis Reaction	Volume
First Strand Synthesis Product (Step 40)	20 μl
NEBNext Second Strand Synthesis Reaction Buffer with dUTP (10X)	8 μl
NEBNext Second Strand Synthesis Enzyme Mix	4 μl
Nuclease-free Water	48 μl
Total Volume	80 μl



NEBNext Second Strand Synthesis Reaction Buffer with dUTP Mix
by [New England Biolabs](#)
Catalog #: [E7426](#)



NEBNext Second Strand Synthesis Enzyme Mix
by [New England Biolabs](#)
Catalog #: [E7425](#)



Nuclease-free Water
by [New England Biolabs](#)
Catalog #: [E7764](#)

Keeping the tube on ice, mix thoroughly by pipetting up and down at least 10 times.

43

44 Incubate in a thermocycler for  01:00:00 at  16 °C with the heated lid set at ≤  40 °C (or off).

Purification of Double-stranded cDNA Using SPRIselect Beads or NEBNext Sample Purification Beads



45 Vortex SPRIselect Beads or NEBNext Sample Purification Beads to resuspend.



NEBNext Sample Purification Beads

by [New England Biolabs](#)



Catalog #: [E7767](#)

46 Add  144 µl (1.8X) of resuspended beads to the second strand synthesis reaction (~  80 µl). Mix well on a vortex mixer or by pipetting up and down at least 10 times.

47 Incubate for  00:05:00 at room temperature.

48 Briefly spin the tube in a microcentrifuge to collect any sample on the sides of the tube. Place the tube on a magnet to separate beads from the supernatant. After the solution is clear, carefully remove and discard the supernatant. Be careful not to disturb the beads, which contain DNA.

Caution: Do not discard beads.

49 Add  200 µl of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for  00:00:30, and then carefully remove and discard the supernatant.





Fresh 80% Ethanol


50 Repeat Step 49 once for a total of 2 washing steps.

 [go to step #49](#)

51 Air dry the beads for up to 5 minutes while the tube is on the magnetic rack with lid open.

Caution: Do not over-dry the beads. This may result in lower recovery of DNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.

52 Remove the tube from the magnetic rack. Elute the DNA from the beads by adding  53 µl 0.1X TE Buffer (provided) to the beads. Mix well on a vortex mixer or by pipetting up and down at least 10 times. Quickly spin the tube and incubate for  00:02:00 at room temperature. Place the tube on the magnetic rack until the solution is clear.



(0.1X) TE Buffer
by [New England Biolabs](#)
Catalog #: [E7763](#)

53 Remove of the supernatant and transfer to a clean nuclease-free PCR tube.

NOTE


If you need to stop at this point in the protocol samples can be stored at -20°C .

End Prep of cDNA Library


54 Assemble the end prep reaction on ice by adding the following components to the second strand synthesis product from Step 53.

End Prep Reaction	Volume
Second Strand Synthesis Product (Step 52)	50 µl
NEBNext Ultra II End Prep Reaction Buffer	7 µl
NEBNext Ultra II End Prep Enzyme Mix	3 µl
Total Volume	60 µl

If a master mix is made, add of master mix to of cDNA for the End Prep reaction.



NEBNext Ultra II End Prep Reaction Buffer
by [New England Biolabs](#)
Catalog #: [E7647](#)



NEBNext Ultra II End Prep Enzyme Mix
by [New England Biolabs](#)
Catalog #: [E7646](#)

55 Set a 100 µl or 200 µl pipette to 50 µl and then pipette the entire volume up and down at least 10 times to mix thoroughly. Perform a quick spin to collect all liquid from the sides of the tube.

NOTE

It is important to mix well. The presence of a small amount of bubbles will not interfere with performance.

56 Incubate the sample in a thermocycler with the heated lid set at \geq as follows.

at

at

Hold at

57 Proceed immediately to Adaptor Ligation.

Adaptor Ligation

- 58 [!] Dilute the NEBNext Adaptor* prior to setting up the ligation reaction in ice-cold Adaptor Dilution Buffer and keep the adaptor on ice.

Total RNA Input	Dilution Required
100 ng– 10 ng	25–fold dilution in Adaptor Dilution Buffer

*The NEBNext adaptor is provided in NEBNext oligos kit. NEB has several oligo kit options, which are supplied separately from the library prep kit.

- 59 Assemble the ligation reaction **on ice** by adding the following components, in the order given, to the end prep reaction product from Step 5.

Ligation Reaction	Volume
End Prepped DNA (Step 57)	60 µl
Diluted Adaptor (Step 57)	2.5 µl
NEBNext Ligation Enhancer	1 µl
NEBNext Ultra II Ligation Master Mix	30 µl
Total Volume	93.5 µl

NOTE

The Ligation Master Mix and Ligation Enhancer can be mixed ahead of time and is stable for at least 8 hours @ 4°C. We do not recommend premixing the Ligation Master Mix, Ligation Enhancer and adaptor prior to use in the Adaptor Ligation Step.



NEBNext Ligation Enhancer

by New England Biolabs

Catalog #: [E7374](#)



NEBNext Ultra II Ligation Master Mix

by New England Biolabs

Catalog #: [E7648](#)

- 60 Set a 100 µl or 200 µl pipette to 80 µl and then pipette the entire volume up and down at least 10 times to mix thoroughly. Perform a quick spin to collect all liquid from the sides of the tube.

Caution: The NEBNext Ultra II Ligation Master Mix is very viscous. Care should be taken to ensure adequate mixing of the ligation reaction, as incomplete mixing will result in reduced ligation efficiency. The presence of a small amount of bubbles will not interfere with performance.

- 61 Incubate  00:15:00 at  20 °C in a thermocycler.

- 62 Add  3 µl (blue) USER Enzyme to the ligation mixture from Step 60, resulting in total volume of  96.5 µl



NEBNext USER Enzyme

by New England Biolabs


Catalog #: [E7458](#)

- 63 Mix well and incubate at  37 °C for  00:15:00 with the heated lid set to ≥  45 °C

Proceed immediately to Purification of the Ligation Reaction.

Purification of the Ligation Reaction Using SPRIselect Beads or NEBNext Sample Purification Beads


65 [!] Note: If you are selecting for libraries with larger insert size (> 200 nt) follow the size selection recommendations in Appendix A, Chapter 6.

Add  **87 µl** (0.9X) resuspended SPRIselect Beads or NEBNext Sample Purification Beads and mix well on a vortex mixer or by pipetting up and down at least 10 times.





NEBNext Sample Purification Beads
by New England Biolabs
Catalog #: [E7767](#)

66 Incubate for  **00:10:00** at room temperature.

67 Quickly spin the tube in a microcentrifuge and place the tube on an appropriate magnetic rack to separate beads from the supernatant. After the solution is clear (~  **00:05:00**), discard the supernatant that contains unwanted fragments.

Caution: Do not discard beads.

68 Add  **200 µl** of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for  **00:00:30**, and then carefully remove and discard the supernatant.



80% Ethanol (freshly prepared)



69 Repeat Step 68 once for a total of 2 washing steps.


 **go to step #68 Repeat Step**

70 Briefly spin the tube, and put the tube back in the magnetic rack.

71 Completely remove the residual ethanol, and air dry beads until the beads are dry for up to 5 minutes while the tube is on the magnetic rack with the lid open.

Caution: Do not over-dry the beads. This may result in lower recovery of DNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.

72 Remove the tube from the magnetic rack. Elute DNA target from the beads by adding  **17 µl** 0.1X TE (provided) to the beads. Mix well on a vortex mixer or by pipetting up and down. Quickly spin the tube and incubate for  **00:02:00** at room temperature. Put the tube in the magnet until the solution is clear.

73 Without disturbing the bead pellet, transfer  15 µl of the supernatant to a clean PCR tube and proceed to PCR enrichment.

NOTE

If you need to stop at this point in the protocol samples can be stored at -20°C.

PCR Enrichment of Adaptor Ligated DNA

74 [!] Check and verify that the concentration of your oligos is 10 µM on the label.

[!] Use **Option A** for any NEBNext oligos kit where index primers are supplied in tubes. These kits have the forward and reverse primers supplied in **separate** tubes.

Use **Option B** for any NEBNext oligos kit where index primers are supplied in a 96-well plate format. These kits have the forward and reverse (i7 and i5) primers **combined**.

75 Set up the PCR reaction as described below based on the type of oligos (PCR primers) used.

Option A: Forward and Reverse Primers Separate

Component	Volume Per One Library
Adaptor Ligated DNA (Step 73)	15 µl
NEBNext Ultra II Q5 Master Mix	25 µl
Universal PCR Primer/i5 Primer*,**	5 µl
Index (X) Primer/i7 Primer*,**	5 µl
Total Volume	50 µl

Option B: Forward and Reverse Primers Combined

Component	Volume Per One Library
Adaptor ligated DNA (Step 73)	15 µl
NEBNext Ultra II Q5 Master Mix	25 µl
Index (X)/i7 Primer Mix*	10 µl
Total Volume	50 µl

* NEBNext Oligos must be purchased separately from the library prep kit. Refer to the corresponding NEBNext Oligo kit manual for determining valid barcode combinations.

** Use only one i7 primer/ index primer per sample. Use only one i5 primer (or the universal primer for single index kits) per sample

76 Mix well by gently pipetting up and down 10 times. Quickly spin the tube in a microcentrifuge.

77 Place the tube on a thermocycler with the heated lid set to  105 °C and perform PCR amplification using the following PCR cycling conditions (refer to Table 77.A and Table 77.B):

Table 77.A:

Cycle Step	Temp	Time	Cycles
Initial Denaturation	98°C	30 seconds	1
Denaturation	98°C	10 seconds	12–16*,**
Annealing/Extension	65°C	75 seconds	
Final Extension	65°C	5 minutes	1
Hold	4°C	∞	

* The number of PCR cycles should be adjusted based on RNA input.



** It is important to limit the number of PCR cycles to avoid overamplification. If overamplification occurs, a second peak ~ 1,000 bp will appear on the Bioanalyzer trace (see Figure 7.2 in manual).

Table 77.B: Recommended PCR cycles based on input amount:

FFPE RNA Input	Recommended PCR Cycles
100 ng	12–13
10 ng	15–16

Purification of the PCR Reaction using SPRIselect Beads or NEBNext Sample Purification Beads

78 Vortex SPRIselect Beads or NEBNext Sample Purification Beads to resuspend.

79 Add  45 µl (0.9X) of resuspended beads to the PCR reaction (~  50 µl). Mix well on a vortex mixer or by pipetting up and down at least 10 times.




NEBNext Sample Purification Beads



by [New England Biolabs](#)

Catalog #: [E6315](#)

80 Incubate for  00:05:00 at room temperature.

81 Quickly spin the tube in a microcentrifuge and place the tube on an appropriate magnetic rack to separate beads from the supernatant. After the solution is clear (~  00:05:00), carefully remove and discard the supernatant. Be careful not to disturb the beads that contain DNA targets.

Caution: Do not discard beads.

82 Add  200 µl of freshly prepared 80% ethanol to the tube while in the magnetic rack. Incubate at room temperature for  00:00:30, and then carefully remove and discard the supernatant.



80% Ethanol (freshly prepared)

83 Repeat Step 82 once for a total of 2 washing steps.

 [go to step #82 Repeat Step](#)

84 Air dry the beads for up to 5 minutes while the tube is on the magnetic rack with the lid open.

Caution: Do not over-dry the beads. This may result in lower recovery of DNA target. Elute the samples when the beads are still dark brown and glossy looking, but when all visible liquid has evaporated. When the beads turn lighter brown and start to crack, they are too dry.

85 Remove the tube from the magnetic rack. Elute the DNA target from the beads by adding  23 µl 0.1X TE (provided) to the beads. Mix

well on a vortex mixer or by pipetting up and down ten times. Quickly spin the tube in a microcentrifuge and incubate for 00:02:00 at room temperature. Place the tube in the magnetic rack until the solution is clear.

- 86 Transfer 20 µl of the supernatant to a clean PCR tube, and store at -20°C .

Assess Library Quality on an Agilent Bioanalyzer DNA Chip


- 87 Run 1 µl library on a DNA High Sensitivity Chip. A dilution may be necessary for running on a Bioanalyzer High Sensitivity DNA Chip

- 88 Check that the electropherogram shows a narrow distribution with a peak size approximately 300 bp.

NOTE

Note: If a peak at ~ 80 bp (primers) or 128 bp (adaptor-dimer) is visible in the bioanalyzer traces, bring up the sample volume (from Step 86) to 50 µl with 0.1X TE buffer and repeat the SPRIselect Bead or NEBNext Sample Purification Bead Cleanup Step (Section "Purification of the PCR Reaction using SPRIselect Beads or NEBNext Sample Purification Beads").

- 89 **Figure 89: Example of RNA library size distribution on a Bioanalyzer.**

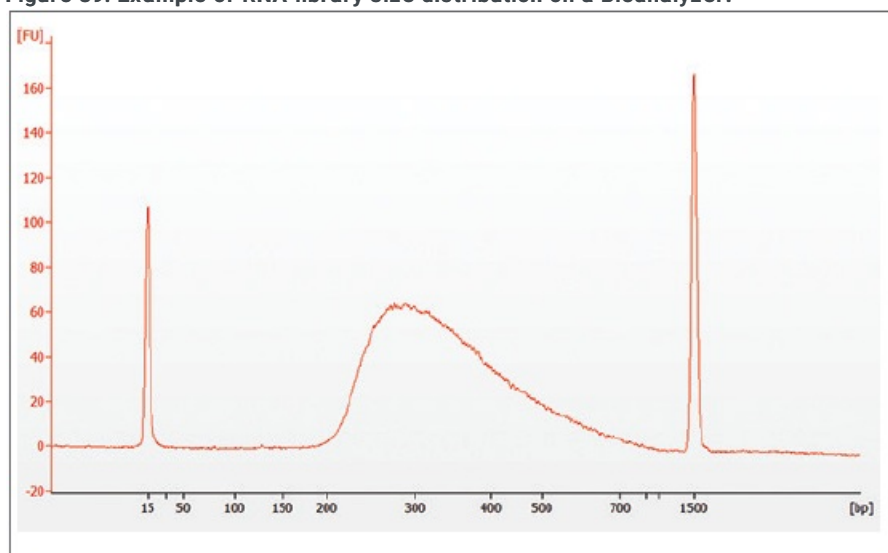

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited