

Modified ZN Staining Protocol Version 2

Asar Khan, Sumaira Shams, Saima Khan, Muhammad Iftikhar Khan, Sardar Khan, Abid Ali

Abstract










The Modified Ziehl-Neelsen stain (mZN stain) is a type of differential bacteriological stain used to identify acid-fast organisms, mainly *Mycobacteria*. Acid fast organisms are those which are capable of retaining the primary stain when treated with an acid (*fast=holding capacity*). Members of the Actinomycetes, genus *Nocardia* (*N. brasiliensis* and *N. asteroides* are opportunistic pathogens) are partially acid-fast. Oocysts of coccidian parasites, such as *Cryptosporidium* and *Isospora*, are also acid-fast. Hence they can also be detected and identified through mZN staining procedure.

Citation: Asar Khan, Sumaira Shams, Saima Khan, Muhammad Iftikhar Khan, Sardar Khan, Abid Ali Modified ZN Staining Protocol. **protocols.io**

dx.doi.org/10.17504/protocols.io.tb2eiqe

Published: 07 Sep 2018


Materials

-  Carbol-Fuchsin by Contributed by users
-  Distilled Water by Contributed by users
-  Methanol M3641 by [Sigma Aldrich](#)
-  Disposable Latex Gloves, Medium, 100/Box [GL002M.SIZE.1PK](#) by [Bio Basic Inc.](#)
-  Methylene Blue [M-680](#) by [Gold Biotechnology](#)
-  Microscope slides by Contributed by users
-  Compound Microscope by Contributed by users
-  ethanol by [BBI Biotech](#)
-  Acid Alcohol by Contributed by users

Protocol

Step 1.

The stool sample was Spread evenly on the middle of the slide with constant rotational movement.

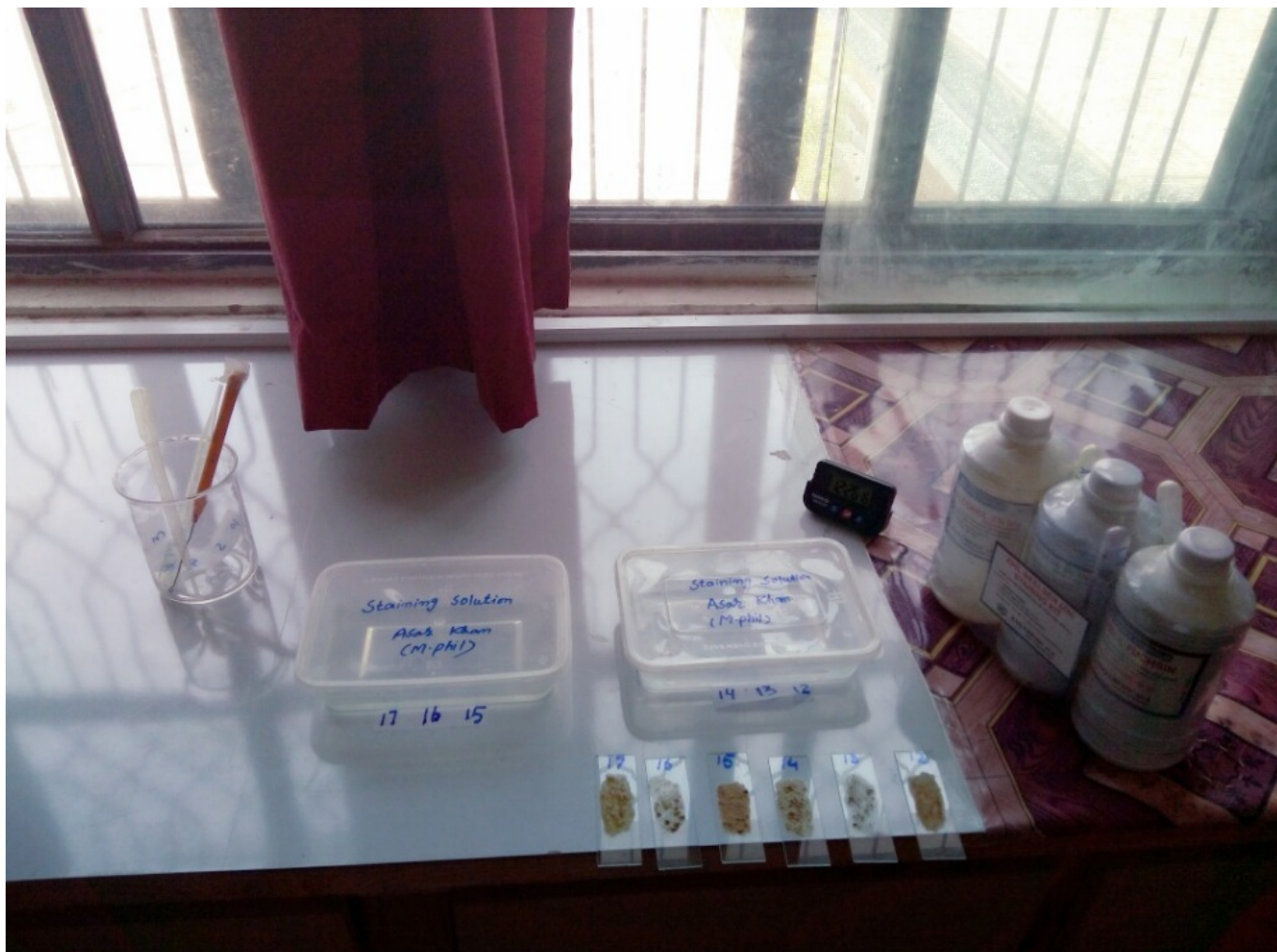

☐ AMOUNT

3 mg : (Amount of stool sample)

🕒 DURATION

00:10:00 : (5 to 10 minutes) for rotational movement

Step 2.

The slides were then placed on dryer with smeared surface upwards to air-dried them.

🌡️ TEMPERATURE

60 °C :

🕒 DURATION

00:10:00 : minutes

Step 3.

The dried smear was fixed with absolute methanol.

🕒 DURATION

00:05:00 : or (3-5 minutes)

Step 4.

Now, the Carbol-fuchsin solution was added to the slide to cover the whole smear.



REAGENTS

Carbol-Fuchsin by Contributed by users

DURATION

00:20:00 : minutes

Step 5.

The slides were washed gently with tap water with the help of a dropper.

SAFETY INFORMATION

Do not expose the slides to the high pressure of tap water directly, rather it will be better to use a dropper for washing the slides.

Step 6.

After washing the slide, decolorizer (Acid Alcohol) was added to the smear and the slide washed again with tap water.

AMOUNT

3 ml : or 4-6 drops


REAGENTS

Acid Alcohol by Contributed by users

Step 7.

Then the counter stain (Methylene Blue) was added and left for 5 minutes and then washed the slide with clean water.

REAGENTS

 Methylene Blue [M-680](#) by [Gold Biotechnology](#)

DURATION

00:05:00 : minutes wait for methylene blue

Step 8.


The back side of the slides were cleaned with a tissue paper and put in the draining rack to air-dry.

DURATION

00:05:00 : minutes, wait for slide to dry

Step 9.

The smear was examined with the help of a compound microscope with 40x and 100x (immersion oil lens) objective and scanned thoroughly for parasite identification.

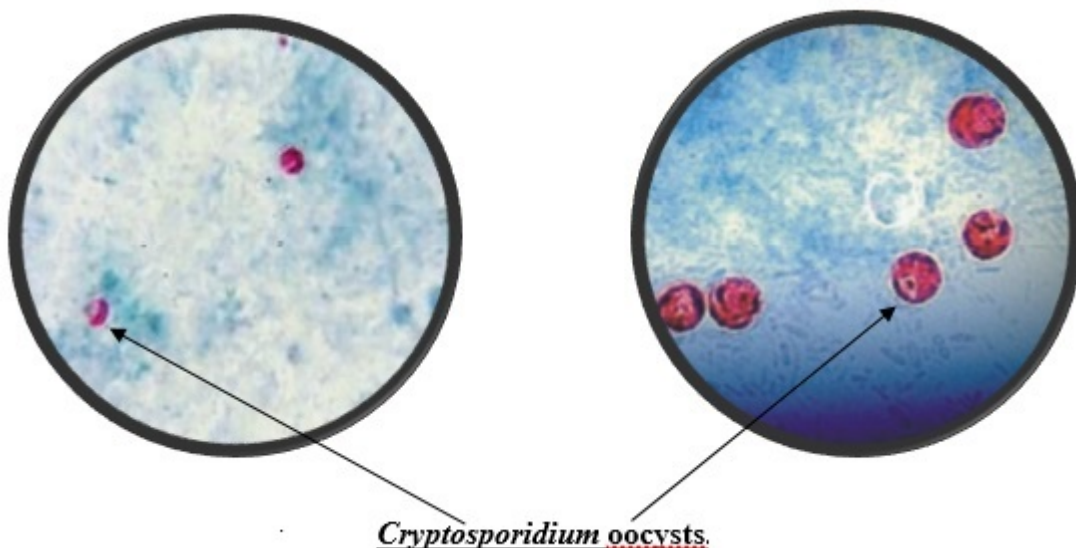

EQUIPMENT

Equipment brand:

Olympus

SKU:

CH20i

Specifications:

Biological microscope , Anti-fungus treated optics , Built to last- Superior build quality

 [31qe3qqpnzl.jpg](#)