



Oct 04, 2019

Human primary T cell culture media V.2

Bulent Arman Aksoy¹, Pinar Aksoy¹, Megan Wyatt¹, Chrystal M. Paulos¹, Jeff Hammerbacher¹¹Medical University of South Carolina

1 Works for me dx.doi.org/10.17504/protocols.io.7yihpue

Hammer Lab

Tech. support phone: +18437924527 email: arman@hammerlab.org

Bulent Arman Aksoy
Medical University of South Carolina

MATERIALS


NAME	CATALOG #	VENDOR
EMD Millipore™ Stericup™ Sterile Vacuum Filter Units	SCGPU05RE	Fisher Scientific
Corning™ RPMI 1640 Medium (Mod.) 1X with L-Glutamine	MT10041CV	Fisher Scientific
HyClone™ HEPES Solution	SH3023701	Fisher Scientific
Fetal Plus®	FP-0500-A	Atlas Biologicals
Penicillin-Streptomycin (10,000 U/mL)	15140122	Thermo Fisher Scientific
Gibco™ Sodium Pyruvate (100mM)	11-360-070	Fisher Scientific
HyClone™ Non Essential Amino Acids NEAA 100x solution	SH3023801	Fisher Scientific

Mix components

- Mix all in [a 500 ml 0.22 um filter bottle](#) in the hood:
- 500 ml RPMI-1640 ([Corning™ cellgro™ RPMI 1640 Medium \(Mod.\) 1X with L-Glutamine](#))
- 25 ul of [EMD Millipore™ Calbiochem™ β-Mercaptoethanol, Molecular Biology Grade](#) (715 uM)
- 12.5 ml of [HyClone™ HEPES Solution](#) (25 mM)
- 50 ml of FBS ([Fetal Plus®](#))
- 5 ml of [Penicillin-Streptomycin](#) (1%)
- 5 ml of [HyClone™ 100mM Sodium Pyruvate Solution](#) (1X) 5 ml of [HyClone™ Non Essential Amino Acids NEAA 100x solution](#) (1X)

Filter and aliquot

- 8 Sterile filter (0.22um) and aliquot in 50 ml falcon tubes. Keep the media in the fridge, protect from light.

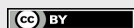

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited