

Oct 19, 2018

Working

# Dissociation of Single Cell Suspensions from Human Breast Tissues

Shanshan Bai<sup>1</sup>, Emi Sei<sup>1</sup>, Nicholas E. Navin<sup>1</sup><sup>1</sup>University of Texas MD Anderson Cancer Center

dx.doi.org/10.17504/protocols.io.t3geqjw

Nicholas Navin

## ABSTRACT

This protocol was developed for the human breast cell atlas (HBCA) project to obtain high-viability cell suspensions from freshly dissociated breast tissues from human patients. There are two options for performing this protocol: *rapid-dissociation* (15-30 min) or *exhaustive dissociation* (overnight). In our experience the rapid dissociation protocol will provide higher cell viability and better representation of breast cell types, however the total number of cells dissociated may be lower after enumeration. In contrast the exhaustive protocol provides higher cell numbers and maintains good cell viability but may result in skewing of cell type representations over the long time period. The rapid-dissociation protocol was used to generate most of the cell suspensions for the MD Anderson human breast cell atlas project and has been optimized and tested more extensively.

## TAGS

Single-cell

breast


 NL\_Breasttissue  
dissociation\_protocol\_  
v2.docx

## PROTOCOL STATUS

### Working

We use this protocol in our group and it is working

## MATERIALS

NAME	CATALOG #	VENDOR
Flowmi™ Cell Strainer 40 µm	H13680-0040	Bel-Art
Ultrapure BSA	AM2616	Ambion
10x MACS RBC lysis buffer	130-094-183	Miltenyi Biotec
Collagenase A	11088793001	Sigma
DMEM F12/HEPES	113300	Gibco - Thermo Fischer
DMEM medium	D5796	Sigma
BSA fraction V	15260037	Gibco - Thermo Fischer
Cell strainer 70um filter	352350	Falcon

## MATERIALS TEXT

### Consumables:

- Scalpels, 1.5ml tubes, 50ml conical, syringe, PBS.

## SAFETY WARNINGS

Please refer to the SDS (Safety Data Sheet) for safety warnings and hazard information.

### Making 10xcollagenase solution

- 1 Dissolve  in .

### Dissociation solution (DS)

- 2

### Mincing

- 3 In a 10 cm dish, place a 1-3cm<sup>3</sup> tissue section and add ~ . Mince the tissue with scalpels into 1 mm<sup>3</sup> pieces.
- 4 Transfer minced tissue to a 50 ml conical tube. Q.S. to 35 ml with DS, enough to cover all the tissue.

### Dissociation

- 5 Two options: **Rapid-dissociation** or **Exhaustive Dissociation**

step case

#### Rapid-dissociation

no description provided

- 6 Gently dissociate the suspension at  in a hybridization oven for  - .

#### NOTE

With rotator set to medium (optional: pipet up and down with a 5 ml pipette to aid digestion during the incubation period).

- 7 Place the 70 µm filter on top of a 50 ml conical tube and filter the tissue suspension. Use a syringe plunger flange to grind the leftover unfiltered tissue. Rinse the filter with DMEM to ensure any remaining single cells are filtered. If necessary, re-filter the solution one more time.

### Spin

- 8 Spin at 430 xg for .

Remove the top adipocyte layer and the rest of the supernatant.

9

**NOTE**

The adipocyte layer can be saved for future analysis by storage in 90% DMSO and 10% FBS. However, the adipocyte cells are usually larger than 50 microns and can therefore not be loaded into the 10X genomics microfluidic chips.

**Resuspending**

10 Resuspend the cell pellet in **1 ml DMEM**.

11 Pass cell suspension through a 40 µm flowmi filter to a 1.5 ml tube.

**Centrifugation**

12 Centrifuge cells at **4 °C** 300 xg for **00:05:00**. If no RBCs are visible, skip directly onto [step 19](#).

**RBC lysis**

13 Remove the supernatant and resuspend the cell pellet in: 1x MACS RBC lysis buffer.

**NOTE**

\*\*\*RBC lysis can be done up to 2 times.\*\*\* You can check the cells under a light microscope to see if RBCs are still present after each lysis step.

14 Incubate nutating at RT for **00:10:00** in **10 ml** - **20 ml 1X buffer**.

**Centrifugation**

15 Add 2X the volume of DMEM into the tube and centrifuge cells at **4 °C** 300 xg for **00:05:00**.

16 Remove the supernatant and resuspend the cell pellet in **1 ml cold PBS**.

17 Filter with a 40 µl flowmi into a 1.5 ml tube.

18 Spin at 300 xg for **00:05:00**.

## Wash

19 After centrifugation, remove the supernatant and wash cells with **1 ml cold PBS**.

20 Centrifuge cells at **4 °C** 300 xg for **00:05:00**.

21 Remove the supernatant.

## For 10X genomics RNA experiments

22 Resuspend cell pellet in PBS+ 0.04% BSA. Adjust cell concentration to 700-1200 cells/μl.

23 To check for **nucleated cells**: Make a 1:1 solution of Hoechst and PI, pipette **2 μl** of this mix to **10 μl filtered cells** and count in the Countess II FL (ThermoFisher) or a Hemocytometer.

24 To check for **viability**: Mix **10 μl Trypan blue** + **10 μl cells**. Load **10 μl of the cell suspension** for counting using a Hemocytometer or the Countess II (ThermoFisher)

## Dissociation

step case

### Exhaustive Dissociation


no description provided

6 Gently dissociate the suspension at **37 °C** in a hybridization oven for **12:00:00**.

#### NOTE

With rotator set to medium (optional: pipet up and down with a 5 ml pipette to aid digestion during the incubation period).

7 Place the 70 μm filter on top of a 50 ml conical tube and filter the tissue suspension. Use a syringe plunger flange to grind the leftover unfiltered tissue. Rinse the filter with DMEM to ensure any remaining single cells are filtered. If necessary, re-filter the solution one more time.

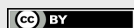

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited