



RNA in situ hybridization 👄

Nature Communications

Nicholas Leigh¹, Garrett Dunlap¹, Kimberly Johnson¹, Rachelle Mariano¹, Rachel Oshiro¹, Alan Y. Wong¹, Donald M. Bryant¹, Bess Miller¹, Jessica L. Whited¹

¹Brigham and Women's Hospital/Harvard Medical School

dx.doi.org/10.17504/protocols.io.p33dqqn



Nicholas Leigh



ABSTRACT

This protocol provides details on tissue harvest and fixation, RNA probe generation, and RNA in situ hybridization for use in sectioned axolotl tissue. This protocol was modified from Brent et al, A somitic compartment of tendon progenitors, Cell 2003

EXTERNAL LINK

https://doi.org/10.1038/s41467-018-07604-0

THIS PROTOCOL ACCOMPANIES THE FOLLOWING PUBLICATION

Leigh ND, Dunlap GS, Johnson K, Mariano R, Oshiro R, Wong AY, Bryant DM, Miller BM, Ratner A, Chen A, Ye WW, Haas BJ, Whited JL, Transcriptomic landscape of the blastema niche in regenerating adult axolotl limbs at single-cell resolution. Nature Communications doi: 10.1038/s41467-018-07604-0

PROTOCOL STATUS

Working

GUIDELINES

All tissue harvest, probe generation, and day 1 of RNA in situ hybridization should be performed using RNase-free reagents and tools.

| Ν. | ЛΛ | т | Е | DI | Λ | 10 |
|----|-----|----|---|-----|---|----|
| I۷ | 1/- | ١. | _ | L/I | М | LO |

| NAME V | CATALOG # | VENDOR ~ |
|--|---------------|--------------------------|
| OneTaq DNA Polymerase - 1,000 units | M0480L | New England Biolabs |
| RNeasy Mini Kit | 74104 | Qiagen |
| RNase Zap | R2020-250ML | Sigma Aldrich |
| Deoxynucleotide Solution Set - 25 umol of each | N0446S | New England Biolabs |
| Formamide, deionized | AB00600-00500 | American Bio |
| BCIP | B6149 | Sigma – Aldrich |
| Lamb Serum, New Zealand origin | 16070096 | Thermo Fisher Scientific |
| T7 RNA Polymerase | 10881767001 | Sigma – Aldrich |
| SP6 RNA Polymerase | 10810274001 | Sigma – Aldrich |
| DIG RNA Labeling Mix | 11277073910 | Sigma – Aldrich |
| Protector RNase Inhibitor | 3335402001 | Sigma – Aldrich |
| pGEM®-T Easy Vector Systems | A1360 | Promega |
| UltraPure™ DNase/RNase-Free Distilled Water | 10977023 | Thermo Fisher Scientific |
| | | |

STEPS MATERIALS



| NAME ~ | CATALOG # | VENDOR ~ |
|---|---------------|--------------------------|
| UltraPure™ DNase/RNase-Free Distilled Water | 10977023 | Thermo Fisher Scientific |
| Deoxynucleotide Solution Set - 25 umol of each | N0446S | New England Biolabs |
| OneTaq DNA Polymerase - 1,000 units | M0480L | New England Biolabs |
| QIAquick Gel Extraction Kit | 28704 | Qiagen |
| Protector RNase Inhibitor | 3335402001 | Sigma – Aldrich |
| T7 RNA Polymerase | 10881767001 | Sigma – Aldrich |
| SP6 RNA Polymerase | 10810274001 | Sigma – Aldrich |
| DIG RNA Labeling Mix | 11277073910 | Sigma – Aldrich |
| RNase Zap | R2020-250ML | Sigma Aldrich |
| RNeasy Mini Kit | 74104 | Qiagen |
| RNase-Free DNase Set | 79254 | Qiagen |
| UltraPure™ DNase/RNase-Free Distilled Water | 10977023 | Thermo Fisher Scientific |
| UltraPure TM DNase/RNase-Free Distilled Water | 10977023 | Thermo Fisher Scientific |
| UltraPure™ DNase/RNase-Free Distilled Water | 10977023 | Thermo Fisher Scientific |
| Tricaine methanesulfonate MS222 | | Sigma Aldrich |
| Tissue-Tek® O.C.T. Compound, Sakura® Finetek | 25608-930 | Vwr |
| RNase Zap | R2020-250ML | Sigma Aldrich |
| Tissue-Tek® O.C.T. Compound, Sakura® Finetek | 25608-930 | Vwr |
| Superfrost Plus Microscope Slies | 4951PLUS4 | Thermo Fisher Scientific |
| Proteinase K, recombinant, PCR grade | 3115836001 | Sigma Aldrich |
| Acetic anhydride | 320102-1L | Sigma – Aldrich |
| Triethanolamine hydrochloride | T1502 | Sigma – Aldrich |
| Formamide, deionized | AB00600-00500 | American Bio |
| Bel-Art™ Polyethylene Utility Bags | 22-260173 | Fisher Scientific |
| Texwipe™ KOGO™ TPR Cleanroom Adhesive Tapes | 18-366-432 | Fisher Scientific |
| Plastic staining dish | 25608-904 | VWR international Ltd |
| RNase A | 10109142001 | Sigma – Aldrich |
| Lamb Serum, New Zealand origin | 16070096 | Thermo Fisher Scientific |
| Blocking Reagent | 11096176001 | Sigma – Aldrich |
| Anti-Digoxigenin-AP, Fab fragments | 11093274910 | Sigma – Aldrich |
| 5-Bromo-4-chloro-3-indolyl phosphate disodium salt | B6149 | Sigma – Aldrich |
| Nitrotetrazolium Blue chloride | N6639 | Sigma – Aldrich |
| Polyvinyl alcohol mounting medium with DABCO®, antifading | 10981 | Sigma – Aldrich |
| Thermo Scientific™ Richard-Allan Scientific™ Cover Glass | 102460 | Thermo Fisher Scientific |

| NAME V | CATALOG # | VENDOR V |
|--------------------------------|---------------|--------------------------|
| Depc (Diethyl Pyrocarbonate) | 80730-888 | VWR international Ltd |
| Lamb Serum, New Zealand origin | 16070096 | Thermo Fisher Scientific |
| Dextran Sulfate, 50% Solution | AB00426-00200 | American Bio |
| Yeast tRNA | 15401029 | Thermo Fisher Scientific |
| Formamide, deionized | AB00600-00500 | American Bio |
| Ficoll® 400 | F2637 | Sigma – Aldrich |
| Polyvinylpyrrolidone | PVP40 | Sigma – Aldrich |
| Triethanolamine hydrochloride | T1502 | Sigma – Aldrich |

I. Tissue Preparation and Sectioning

Narcotize animals in 0.1% Tricaine prepared in 40% Holtfreter's. Remove tissue of interest. Deposit directly into 4% PFA prepared in 1X DEPC PBS.

NOTE

All solutions and tools/bottles should be RNase-free, being prepared with DEPC-treated water or being cleaned with RNase Zap, respectively.

After survival surgeries always keep animals overnight in 40% Holtfreter's supplemented with 0.5% sulfamerazine.



Tricaine methanesulfonate MS222

by Sigma Aldrich

- 7 Fix tissues in 4% PFA prepared in 1X DEPC PBS for 1-2 hr at room temperature with gentle rocking or overnight at 4°C with gentle rocking.
- 3 Remove 4% PFA/1X DEPC PBS and wash in 1X DEPC PBS by gently rocking tissue in 1X DEPC PBS for 30 min.
- Transfer tissue through a sucrose series, starting with 5% sucrose in 1X DEPC PBS, then 10%, 20%, and eventually 30%. Change to higher percent sucrose when tissue sinks.

■NOTE

Once tissue sinks in 30% sucrose in 1X DEPC PBS proceed to the next step. Alternatively, tissue can be rocked overnight in 30% at 4° C.

Add approximately equal volume of OCT compound to the 30% sucrose/1X DEPC PBS. Use a gentle inversion rocker to let it equilibrate for at least 30 minutes.



Tissue-Tek® O.C.T. Compound, Sakura®

Finetek

by Vwr

Catalog #: 25608-930

6 Prepare a dry-ice/ethanol "bath" and RNase Zap forceps.



7 Remove tissue from 30% sucrose/OCT solution and move into tissue mold. Fill mold with OCT and orient tissue as desired.



- 8 Place mold with tissue into dry ice/ethanol bath to freeze OCT.
- Q After freezing, move tissue into plastic bags and then store at -80°C until sectioning
- 10 Section tissue on a cryostat. Sections can be 10-20 microns (usually 16 microns). Collect sections on Superfrost Plus slides.



11 Air-dry the sections 30 min-2 hrs. Store at -80°C or continue on to hybridization.

II. Probe Preparation

Perform PCR using OneTaq and with T7 and Sp6 primers using pGEM-T easy with probe inserted as template. Do 5+ reactions to ensure ample PCR product. Annealing temperature 40.5°C, extension time time appropriate for length of probe (typically around 30 seconds).

| Reagent | Volume (μL) |
|---------------------------------------|-------------|
| 5x OneTaq Standard Buffer | 10 |
| 10mM dNTP | 1 |
| 10μm T7 primer (TAATACGACTCACTATAGGG) | 1 |
| 10μm Sp6 primer (ATTTAGGTGACACTATAG) | 1 |
| OneTaq | 0.25 |
| 2ng/μL pGEM-T Easy | 1 |
| w/ probe inserted | |
| Ultrapure H ₂ 0 | 35.75 |
| Reaction Volume | 50 |

| Ė | UltraPure™ DNase/RNase-Free Distilled |
|---|---------------------------------------|
| | Water |
| | by Thermo Fisher Scientific |
| | Catalog #: 10977023 |
| | |

Deoxynucleotide Solution Set - 25 umol of each by New England Biolabs Catalog #: N0446S

OneTaq DNA Polymerase - 1,000 units by New England Biolabs Catalog #: M0480L

13 Gel extract PCR product and calculate concentration of gel extracted product.




Perform in vitro transcription (3 reactions) with either T7 or Sp6 polymerase depending on anti-sense orientation of probe in pGEM-T Easy. Prepare below reaction and incubate at 37°C for 2 hours.

| Reagent | Volume (μL) |
|----------------------------|-------------|
| 10x Transcription Buffer | 2 |
| 10x DIG | 2 |
| Protector RNase Inhibitor | 0.5 |
| Sp6 or T7 RNA polymerase | 1 |
| template (500ng) | variable |
| Ultrapure H ₂ 0 | to 20µL |

Protector RNase Inhibitor
by Sigma – Aldrich
Catalog #: 3335402001

T7 RNA Polymerase
by Sigma – Aldrich
Catalog #: 10881767001


SP6 RNA Polymerase
by Sigma – Aldrich
Catalog #: 10810274001

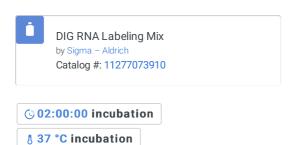


NOTE

Use RNase-free reagents for all steps and thoroughly RNase Zap all surfaces and pipettes for this step and all subsequent steps of probe preparation.

Prepare a gel mold and electrophoresis box by spraying thoroughly with RNase Zap followed by washing out with MilliQ water. In addition, RNase Zap and rinse with MilliQ the interior of a bottle and prepare a 2% agarose gel in this RNase-free bottle.

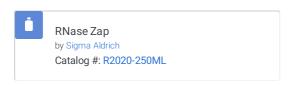

16 After 2 hour incubation, proceed to RNeasy RNA cleanup with on column DNase I digestion.

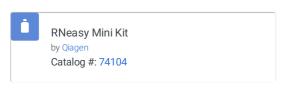



17 $\,\,$ Combine all three reactions and volume up to $100\mu L$ with Ultrapure $H_20.$



- 18 Add $350\mu L$ RLT buffer and mix well.
- 19 Add $250\mu L100\%$ ethanol and mix well by pipetting. Proceeded immediately to the next step.
- Transfer 700μ L to a RNeasy Mini spin column that is in a 2mL collection tube, close lid and centrifuge for 15s at $\underline{10,000 \text{ rpm}}$ (8,000g). Discard flow through.
- 21 Add $350\mu L$ Buffer RW1, close lid, and centrifuge for 15s at 10,000 rpm.

Combine 10µL DNase I stock solution and 70µL Buffer RDD and mix gently by inverting the tube. 22 Add 80µL of DNase I/Buffer RDD solution directly to the RNeasy column membrane, close lid, and let it sit on the bench for 15min at room 23 temperature. **© 00:15:00** incubation Add 350µL Buffer RW1 to column, close lid, centrifuge for 15s at 10,000 rpm. Discard flow-through. 24 Add 500µL Buffer RPE to column, close lid, and centrique for 15s at 10,000 rpm. Discard flow-through. 25 Add 500µL Buffer RPE to column, close lid, and centrifuge for 2 min at 10,000 rpm. Discard flow-through. 26 **© 00:02:00** centrifugation Place RNeasy column into a new 2mL collection tube, close lid, and centrifuge at max speed for 1 min. 27 Place RNeasy column into a 1.5mL collection tube (supplied with kit) and add 50μL of Ultrapure H₂0 directly to the column membrane. Close 28 the lid and centrifuge for 1 min at 10,000 rpm. UltraPure™ DNase/RNase-Free Distilled by Thermo Fisher Scientific Catalog #: 10977023 29 Remove eluate from collection and re-apply to the column for a second wash. 30 Remove 2ul of eluate and put into a new 1.5mL tube. Store the remainder of the eluate at -80°C until use. 31 Add $8\mu L$ of Ultrapure H_20 and incubate at $60^{\circ}C$ for 10 minutes. **© 00:10:00** incubation

32 Add a clear loading dye (e.g. 40% sucrose) and run on the 2% RNase-free agarose gel to confirm size and quality of product.

UltraPure™ DNase/RNase-Free Distilled

by Thermo Fisher Scientific Catalog #: 10977023

III. Hybridization



Depc (Diethyl Pyrocarbonate)

by VWR international Ltd Catalog #: 80730-888



Lamb Serum, New Zealand origin

by Thermo Fisher Scientific Catalog #: 16070096

■NOTE

SOLUTIONS

DEPC water

In a 4L Erlenmeyer flask, put 3L deionized water in, add 3 mL DEPC, stir overnight (in the hood, with a little foil cap loosely on top). Autoclave the next day. Two 3L flasks should be enough for an in situ.

4% PFA in DEPC PBS

Make a 5X stock = 20% PFA in 5X DEPC PBS. Freeze in aliquots. When defrosting, some PFA/salts might come out of solution. Defrost in a warm bath and shake once in a while. Dilute 1:5 in a Falcon tube to give the correct concentration. Can be repeatedly freeze-thawed.

10X PBS (1 L)

80g NaCl

2g KCl

14.4 g Na₂PO₄

2.4 g KH₂PO₄

adjust to pH 7.4 with HCl and bring volume to 1 L with water

<u>PBT</u>

1X PBS +0.1% Tween-20

Hybridization Solution

10 mM Tris pH 7.5

600 mM NaCl

1 mM EDTA

0.25% SDS

10% Dextran Sulfate

1x Denhardt's (recipe follows)

200 mg/mL yeast tRNA

50% formamide

**Dissolve Tris and EDTA into DEPC-treated water after DEPC has been autoclaved. Adding Tris and EDTA prior to autoclaving will interfere with the ability of DEPC to eliminate RNases.

store at -20°C (in aliquots)

100X Denhardt's solution (100 mL)

2g BSA

2g Ficoll (MW 400,000)

2g Polyvinylpyrolidone (MW 40,000)

15mL 20X SSC (3x final conc.)

volume to 100 mL

20X SSC (1L)

8

175.4g NaCl 88.2g Sodium citrate dihydrate

pH to 7.0, volume, autoclave

10X TEA (triethanolamine HCL)

1M stock pH 8.0

<u>TNE</u>

10mM Tris pH 7.5 500mM NaCl 1mM EDTA

5X MABT (1L)

58g Maleic acid 43.8 NaCl volume up to 1L

pH to 7.5 with NaOH

filter sterilize

when ready to use, dilute 1:5 to make 1X MAB, then add Tween to equal 0.1% Tween

NTM pH 9.5

100mM NaCl 100mM Tris pH 9.5 50mM MgCl₂

NBT stock

50mg/mL in 70% dimethylformamide (store at -20°C, dark)

BCIP stock

25mg/mL in water (store at -20°C, dark)

<u>HISS</u>


Sheep serum, heat inactivated at 56°C for 30 min, store at -20°









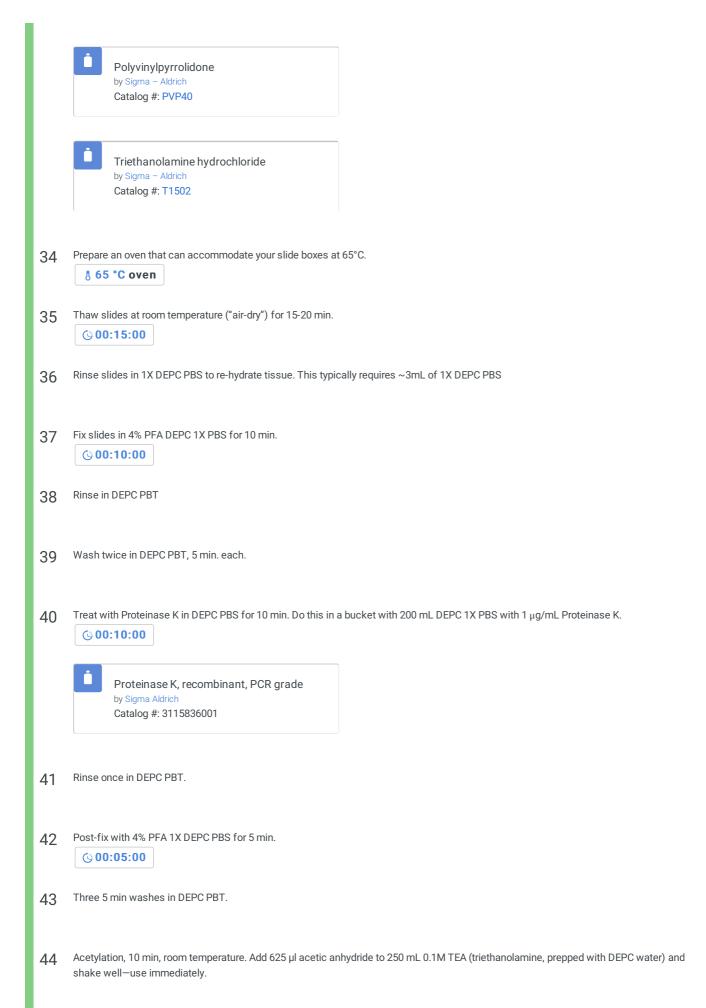



Triethanolamine hydrochloride
by Sigma – Aldrich
Catalog #: T1502

- 45 Rinse in DEPC PBT.
- 46 Wash twice in DEPC PBT, 5 min each.
- 47 Rinse in DEPC H₂O.
- Prepare a humidified slide box. To do this, stand a small slide box on its end and roll a cylinder of paper towels that can fit into the bottom of the box (this will take up about a fifth of the slots). Soak the paper towels in ~25mL of 5X SSC/50% formamide. Set this on your bench with the remaining slots ready to accept the slides.






Take each slide individually and gently let the water wick off onto a Kimwipe. Wipe the edges of the slide with the Kimwipe so that a border of dry slide surrounds the tissue. Continue on to the next step before the tissue sections dry. This step should only take ~1 min/ slide.



Apply the probe. Use 150μl of probe/hybridization buffer for each slide. Gently apply a coverslip-sized piece of plastic bag material (polyethylene) on each slide. Use 2 μl of probe for every 100 μl of hybridization buffer, so here use 3 μl, which should make about 150μl probe. Insert slide into humidified slide box.

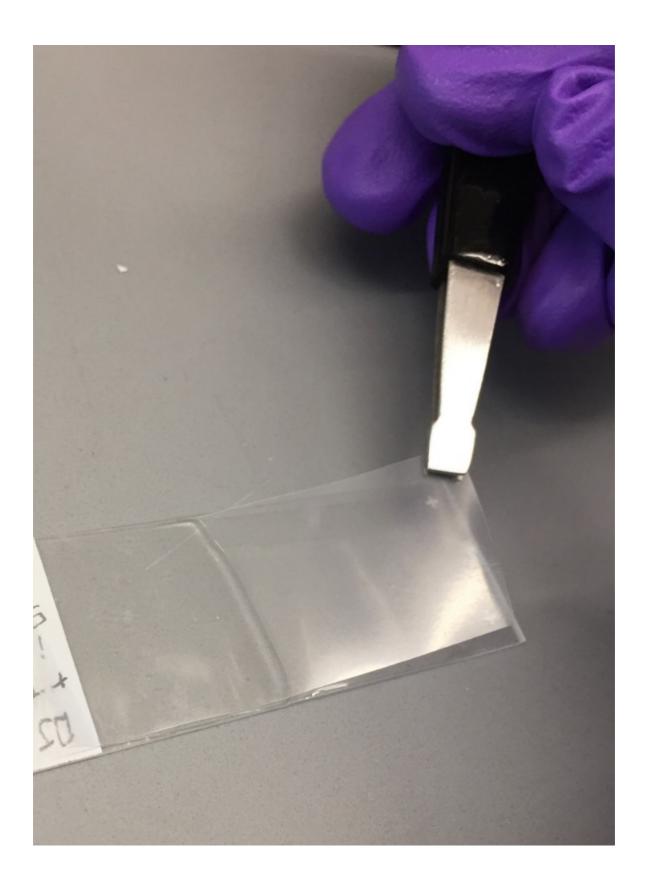



Bel-Art™ Polyethylene Utility Bags by Fisher Scientific Catalog #: 22-260173



Steps 49 and 50 should be repeated for each slide individually. Placing the slides into the slide holder after applying the plastic coverslip.



 $51 \qquad \text{Seal the box with plastic tape along all the edges, being careful not to let the slides tip forward at all.}$





TexwipeTM KOGOTM TPR Cleanroom

Adhesive Tapes by Fisher Scientific Catalog #: 18-366-432


52 Gently transfer the slide box to the prepared 65°C oven. Stand it upright inside. Incubate overnight.

IV. Post-hybridization washes

The next day, prewarm all the solutions to the appropriate temperature. **65°C:** 5X SSC, 2X SSC, 0.2X SSC, 1X SSC/50% formamide.

37°C: TNE

Prepare block: 20% sheep serum, volumed with MABT, with 2% (w/v) blocking reagent dissolved. Start dissolving block in MABT first, then add the sheep serum. Incubate in a water bath at 55°C, shake intermittently to get it into solution. To determine how much block to make multiply your total number of slides by 2 to determine mL of block. Block solution will be used for the block step as well as the incubation.

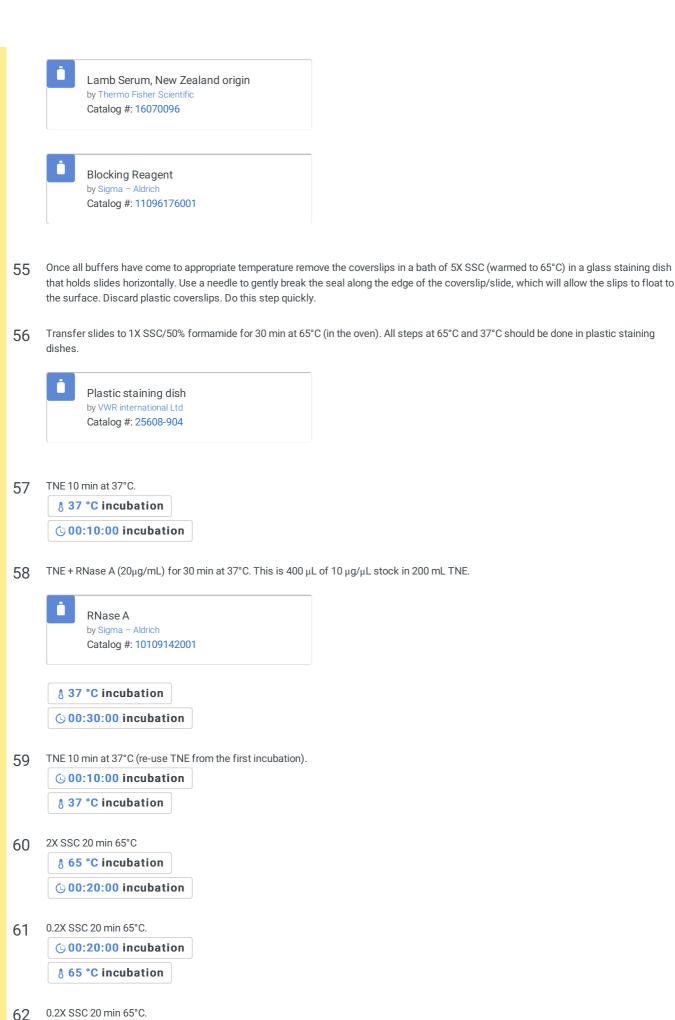

§ 65 °C incubation

V. Antibody Incubation and Detection

- Wash twice in MABT for 5 min. each at room temperature.
- Move slides to a staining tray and add 1mL block/ slide. Block must be left on for a minimum of 1 hour.
- Incubate in a-DIG-AP (1:2000) in 20% HISS/2% blocking reagent(w/v) in MABT (i.e. blocking solution made in step 53), 1mL/slide, overnight, 4°C. Do this in a humidified chamber.




- The next day, rinse once in MABT, then wash three times for 5 min. each in MABT (all room temperature).
- Wash in NTM pH 9.5 for 10 min., room temperature.

© 00:10:00

Develop in NTM pH 9.5 + NBT + BCIP. To each mL NTM, add 4.5 µL NBT and 7 µL BCIP. Prepare 1mL/ slide. Once developing keep slides out of direct light.




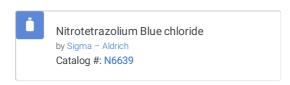

Develop until you see the signal without background coming up (different times for different probes). If necessary to develop overnight, keep covered from light in a humidified chamber in the at 4°C. Developing in the cold can help reduce background even when the probe is strong enough to show good signal after a few hours. Prior to moving to 4°C and slides should be washed with NTM. New developing solution should be prepared when incubating at 4°C. This same wash and re-applying new developing solution should be performed if developing continues at room temperature the following day. This can go out for a 2-3 days as long as background doesn't appear.

VI. Post-detection

70 Wash in NTM pH 9.5 for 5 min.

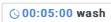

71 Wash in PBS for 5 min.

() 00:05:00 wash

72 Fix in 4%PFA/ PBS for 30 min.

© 00:30:00 fix

- 73 Wash: 1X PBS, 1X PBS, water (5 min. each).
- 74 Dry each slide slightly as described earlier in the protocol (step 48). While sections are still slightly wet, apply mounting media (warmed to 65°C) and coverslip. Avoid bubbles.
 - Polyvinyl alcohol mounting medium with DABCO®, antifading by Sigma Aldrich Catalog #: 10981
 - Thermo Scientific™ Richard-Allan
 Scientific™ Cover Glass
 by Thermo Fisher Scientific
 Catalog #: 102460
- 75 Cure the slides in a dry place overnight.

This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited