

Sep 20, 2019

Collection and Post-Surgical Excision of Human Kidney Tissue through the Cooperative Human Tissue Network

Maya Brewer¹, Jamie Allen¹, Elizabeth Neumann¹, Agnes Fogo², Raymond Harris², Danielle Gutierrez¹, Mark de Caestecker³, Jeff Spraggins¹

¹Vanderbilt University, ²Vanderbilt University Medical Center, ³Division of Nephrology, Vanderbilt University Medical Center



VU Biomolecular Multimodal Imaging Center

Human BioMolecular Atlas Program (HuBMAP) Method Development Community



ABSTRACT

Scope:

Obtain kidney tissue and metadata about tissue location within the whole kidney for storage and analysis.

Expected Outcome:

A portion of kidney tissue and a series of images that provide information about the original location of the smaller tissue.


MATERIALS TEXT

Digital Camera


STERILAB Inc Grossing Boards 12X18IN, Fisher 50-131-7556

- 1 Kidney tissues are selected based on patient and clinical metadata (e.g. age<65 years old). All patients have consented to allow their tissue to be used for research purposes as part of this study.
- 2 Receive excised kidney from full nephrectomy.
- 3 Cut kidney in half along the longest plane and open like a book.
- 4 Cut out a portion of the kidney (~80 x 60 x 4 mm) that is farthest away from the tumor location.
- 5 Take a picture of the smaller piece of kidney in reference to its original location within the entire kidney if possible.

6 Take only a picture of the smaller piece of kidney.

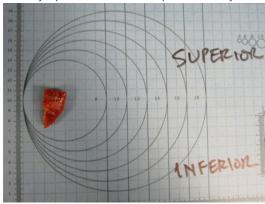

- 7 Cut the piece of kidney into ~4 smaller, rectangular pieces (~2 x 15 x 4 mm).
- 8 Image each of these smaller pieces in reference to their original position.

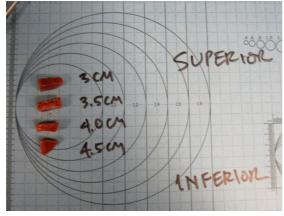

9 Wrap kidneys in gauze and transfer tissue to those preparing tissue for long-term storage.

This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited