



Oct 21, 2019

DIFFPREP for DTI Processing and Corrections

[Courtney Comrie](#)¹¹University of Arizona*In Development*[dx.doi.org/10.17504/protocols.io.7fahjie](https://doi.org/10.17504/protocols.io.7fahjie)

Courtney Comrie



ABSTRACT

This protocol will provide a basic guide to utilizing the DIFFPREP tool in TORTIOSE.

Note: Steps may vary based upon image.

Introduction

- 1 DIFFPREP is a distortion and motion correction module in TORTIOSE. When processing MR images you will commonly start with DIFFPREP for corrections.

Beginning Steps

- 2 Build your directory in the terminal, and make sure to use copied raw data for processing

- 2.1 In the terminal go to your new directory with the data

Import DTI Data

- 3 Type the following command in your terminal to read information and options available:

**ImportBrukerAnatomical**

Type the following command to import your blip_up and blip_down data:

**ImportBrukerAnatomical -i directory/file_name_blip_up -c 1**
ImportBrukerAnatomical -i directory/file_name_blip_down -c 1

- 3.1 Blip up is the forward phase encoding of a scan while blip down is the same exact scan but with reversed phase encoding.

You may not always will have both blip up and blip down, but for the purpose of this protocol both will be shown.

- 4 Change directories to the containers.



```
cd /rsgrps/hutchinsone/Singularity_Containers
```



- 5 Load the singularity containers



```
module load singularity  
singularity run nklab-neurotools-v0.4.sif
```



- 5.1 You should be able to see the contents of the container.

- 6 Start by loading MIPAV and looking at the orientation of your image

Note: You can also use ITK-SNAP and MRtrix to view your imgs, but MIPAV has tools more helpful for this stage of processing

- 6.1 Type:



```
singularity run nklab-neurotools-v0.4.sif mipav
```



- 6.2 On the MIPAV GUI, scroll over to file and import your DTI NIFTY file.
 - 6.3 View the raw data image, and make note of the orientation of the blip down image.

This will become important in step 8.

(Note: The image should naturally be in RPS, and for our purposes of DTI images you want LPS)

Here are a couple sources for understanding orientation.


https://www.slicer.org/wiki/Coordinate_systems

<http://www.grahamwideman.com/gw/brain/orientation/orientterms.htm>

ImportNIFTI

- 7 This step is crucial for obtaining the *.list* file which will be used in DIFFPREP and DRBUDDI

Type the following to get information and options of the command.

```
 ImportNIFTI
```

- 7.1 The general full command is the following:


```
 ImportNIFTI -i directory/file_name.nii -p phase_direction -o output_folder -b  
file_name.bval -v file_name.bvec
```

Note:

- directory: where your NIFTI file is stored
- phase_direction: can either be vertical or horizontal based upon what your phase encoding direction is
- output_folder: where you want your data to be stored when it's done importing


Reorient DTI Images

- 8 In the terminal type:

```
 ReorientImage
```

This will give you information about the command, and the different options available.

- 8.1 Type:


```
 ReorientImage -i directory_name_DTI/file_name_blip_down.list -r RPS
```



- 8.2 Note: This should be done with the list file in your proc folder.

- 8.3 Note: You only need to reorient the blip down data. The blip down data is inverse when it is acquired from the scanner, changing the orientation of the image.

Extract DTI Images


- 9 Extract your reoriented image in the terminal with



 **ExtractImage**

This will give you information about the command, and the different options available.

- 9.1 Type:


 **ExtractImage -i directory_name_DTI/file_name_blip_down_reoriented.nii**



- 10 You can view and edit your extracted image on MIPAV.

Edit Settings File

- 11 Before you run DIFFPREP, you must make sure your settings file is
- Compatible with size and resolution of your data
 - In your user's directory
 - Will give you the corrections needed in DIFFPREP
- 12 In your user's directory make a folder titled "DIFF_PREP_WORK". By doing it in the file browser GUI or by navigating in your terminal and typing:

 **cd /home/u26/cjoy1895**

mkdir DIFF_PREP_WORK


12.1 Copy the setting templates to YOUR OWN DIFF_PREP_WORK



```
cp /rsgrps/hutchingsone/Programming/DIFF_PREP_WORK/IVmouse_2D.dmc  
/home/u26/cjoy1895/DIFF_PREP_WORK
```



- 12.2
- Only need to make the directory and copy the settings file if you haven't done it already
 - There are multiple settings.dmc templates on the HPC, it would be a good idea to copy all of them to your DIFF_PREP_WORK directory
- 13 In your home DIFF_PREP_WORK directory, right click on the IVmouse_2D.dmc and select Open With > gedit

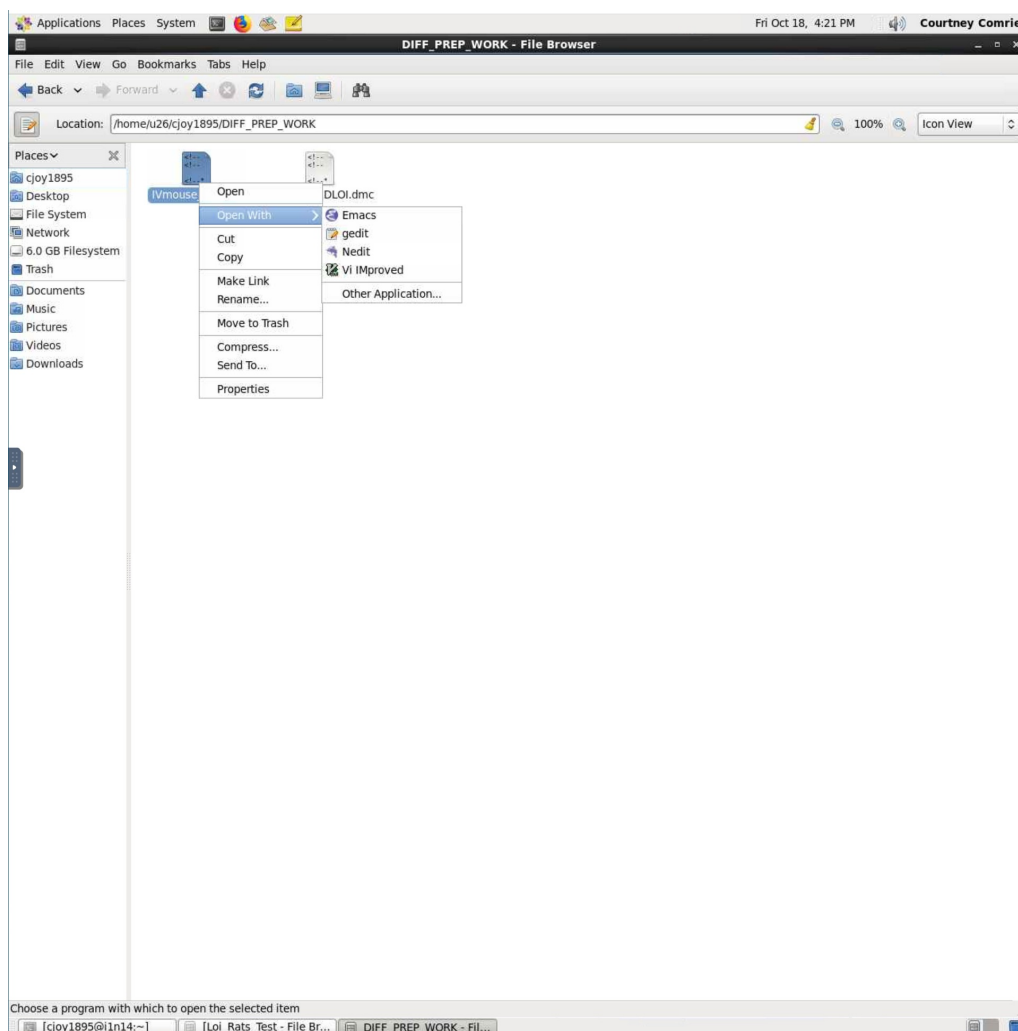

13.1 You should see a bunch of text like the following:

```

@ IVmouse_2D.dmc
<!-- DTI initialization file -->
<!-- do not remove the above comment line -->

<!--*****-->
<!-- CORRECTION MODE -->
<!-- Specifies which undesirable effects will be corrected. Predefined -->
<!-- motion & eddy distortion correction optimization settings (each -->
<!-- setting points to a file in the software settings folder: -->
<!-- off :No DWI registration. -->
<!-- rigid :Corrects only motion. -->
<!-- quadratic :Eddy currents are modeled with quadratic functions. -->
<!-- Quadratic -->
<!-- model is sufficient 99% of the time. -->
<!-- cubic :Eddy currents are modeled with cubic functions. -->
<!-- multistart :Eddy currents are modeled with quadratic functions. -->
<!-- The rigid search space is VERY LARGE. This can be used -->
<!-- to register sagittal images to axial for example. -->
<!-- You can also create a custom settings file in ~/DIFF_PREP_WORK and -->
<!-- put its name here -->
<correction_mode>quadratic</correction_mode>

<!--*****-->
<!-- EPI CORRECTION -->
<!-- Predefined EPI distortion correction optimization settings. Each -->
<!-- setting points to a file in the software settings folder. -->
<!-- off :Do not perform EPI distortion correction. -->
<!-- ITKBSP :Perform BSplines transformation based correction with -->
<!-- ITK. Similar to previous version of DIFFPREP. -->
<!-- ITKBSPv4 :Perform BSplines transformation based correction with -->
<!-- itk4. Multiresolution. More powerful than ITKBSP but -->
<!-- more sensitive. -->
<!-- ANTSSyN :Perform SyN transformation based correction with ANTS. -->
<!-- Most powerful but most sensitive too. -->
<!-- The parameters for each method is described in their corresponding -->
<!-- files. -->
<EPI_optimization_settings>off</EPI_optimization_settings>

<!--*****-->
<!-- REFERENCE VOL INDEX -->
<!-- Among possibly many b=0 s/mm2 images, the index of the b=0 image -->
<!-- (within only b=0 images) to be used as template. Starting index is -->
<!-- zero, so that b0 id=2 will select the 3rd b=0 image. Defaults to the -->
<!-- first b=0 image with id of 0. -->
<b0_id>0</b0_id>

<!--*****-->
<!-- GIBBS RINGING CORRECTION -->
<!-- For the original unprocessed DWIs. See Kellner, Dhital, Kiselev and -->
<!-- Resiart, MRM 2016, 76:1574-1581. -->
<!-- off :Do not perform Gibbs ringing correction. -->
<!-- on :Perform Gibbs ringing correction. -->
<gibbs_ringing_correction>on</gibbs_ringing_correction>

<!--*****-->
<!-- DENOISING -->

```

The DARK BLUE text tells you information about the sections and options you can use, the BLACK text is what you can edit.

Correction Mode

- Change to the registration to **rigid**

EPI Correction

- Make sure it is turned off.
- Only turn it on if you don't plan on using DRBUDDI later.

Gibbs Ringing Correction

- on

Denoising

- for_final

Humanity

- 0
- Basically is asking if your brain is human or not

Upsampling

- all

High B-value Registration

- You need to know the 3 values for the scan you are processing. If you don't know it, you can find it in the methods file.
- Enter your three values

Pre-Smoothing

- Turn off!!! Very important!

FOV

- You need to know the 3 values for the scan you are processing. If you don't know it, you can find it in the methods file.

Resolution

- You need to know the 3 values for the scan you are processing. If you don't know it, you can find it in the methods file.

These are the main settings you will need to change on a scan to scan basis, but keep in mind you may change other settings not mentioned here or will need different inputs than what are mentioned above.

14 Type in terminal:



```
DIFFPREP -i directory/file_name_blip_up.list --reg_settings IVmouse_2D.dmc
```

```
DIFFPREP -i directory/file_name_blip_down_reoriented.list --reg_settings  
IVmouse_2D.dmc
```



The specific settings.dmc file will change based upon the animal and type of image you are processing. This specific .dmc file is for ferrets, but be aware mice will use a different settings file.



This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited